



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
270 S. Stone Avenue Tucson, AZ 85701-1917

TS4  
3-17-09

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-202265 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 3-6-09 Cal 1 0.101 AC Cal 2 0.101 AC

Between the dates of 2-6-09 and 3-6-09 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A  
The instrument was returned to service on N/A

Clare Conley

Clare L. Conley


3-17-09

ADPS QAS Permit # 33314

SUBSCRIBED AND SWORN to before me this 17, day of March, 2009.

A.M. Hernandez

Notary Public

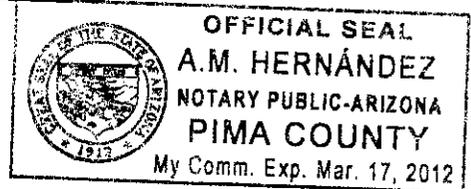

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

1067  
cu

ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000  
STANDARD QUALITY ASSURANCE PROCEDURE

TS4  
3-17-09

QA SPECIALIST Clare Conley AGENCY Tucson Police Department  
DATE 3-6-09 TIME 1443  
INTOXILYZER SERIAL # 80-502265 LOCATION ODD

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 0.101 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS None.  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE Clare Conley

287  
ew  
TSH

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-002265  
03/06/2009

14:45:06

---

Standard Lot#: 801002I  
Last Changed By: S. CIONE #50616

---

QAS: C. CONLEY #34517  
TUCSON PD

---

Operator: C. CONLEY #34517  
TUCSON PD

---

Subject: TEST, SQAP  
DOB: 11/11/1111  
Sex: M Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	14:45:50
Diagnostic Test	Pass	14:46:18
Air Blank	0.000	14:46:45
0.100 Cal Check	0.102	14:47:04
Air Blank	0.000	14:47:33
Subject Test	0.000	14:47:50
Air Blank	ABT*	14:47:57
Air Blank	0.000	14:48:25

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

*Pursuant to section 28-1527, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department*

3 of 7 cu

154

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-002265  
03/06/2009 14:48:30

Standard Lot#: 801002I  
Last Changed By: S. CIONE #50616

QAS: C. CONLEY #34517  
TUCSON PD

Operator: C. CONLEY #34517  
TUCSON PD

Subject: TEST, SQAP  
DOB: 11/11/1111  
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	14:49:01
Air Blank	IPS**	14:49:12

\*RFI Detect  
\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

*Pursuant to section 26-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department*

467  
en

154

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-002265  
03/06/2009

14:49:17

---

Standard Lot#: 801002I  
Last Changed By: S. CIONE #50616

---

QAS: C. CONLEY #34517  
TUCSON PD

---

Operator: C. CONLEY #34517  
TUCSON PD

---

Subject: TEST, SQAP  
DOB: 11/11/1111  
Sex: M Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	14:50:02
Diagnostic Test	Pass	14:50:30
Air Blank	0.000	14:50:57
0.100 Cal Check	0.102	14:51:16
Air Blank	0.000	14:51:45
Subject Test	0.422*	14:54:54
Air Blank	0.000	14:55:26
Air Blank	0.000	14:55:53
Subject Test	XXX**	14:56:16
Air Blank	0.000	14:56:49

\*Deficient Sample  
 \*\*Invalid Sample  
 Begin new deprivation period

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-002265  
03/06/2009

14:57:21

-----  
Standard Lot#: 801002I  
Last Changed By: S. CIONE #50616  
-----

QAS: C. CONLEY #34517  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	14:57:48
Cal Check	0.101	14:58:07
Air Blank	0.000	14:58:37

*Clare Conley*

\_\_\_\_\_  
Operator's Signature

5067  
cu

156

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

677  
cu

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

TS4  
3-17-09

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Clare Conley AGENCY Tucson Police Department  
DATE 3-6-09 TIME 1443  
Intoxilyzer Serial # 80-002265 Location ODD

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 101 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Clare Conley

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-002265  
03/06/2009 14:43:41

767  
cu

154

-----  
Standard Lot#: 801002I  
Last Changed By: S. CIONE #50616  
-----

QAS: C. CONLEY #34517  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	14:44:09
Cal Check	0.101	14:44:27
Air Blank	0.000	14:44:57

*Clare Conley*

\_\_\_\_\_  
Operator's Signature