



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

TS4
7-1-09

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002265 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 6-26-09 Cal 1 0.101 AC Cal 2 0.099 AC

Between the dates of 5-28-09 and 6-26-09 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.
The instrument was returned to service on N/A.

Clare Conley

Clare L. Conley


6-30-09

ADPS QAS Permit # 33314

SUBSCRIBED AND SWORN to before me this 30th, day of June, 2009.

A.M. Hernandez

Notary Public

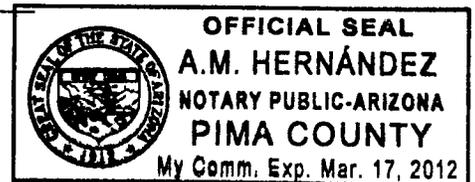

*Includes documents specified in A.A.C. R13-10-104.A.5

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

1 of 7
cu

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

By
7-1-09

QA SPECIALIST Clare Conley AGENCY Tucson Police Department
DATE 6-26-09 TIME 1157
INTOXILYZER SERIAL # 80-002265 LOCATION ODD

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 099 AC

Instrument is operating properly and accurately. YES NO

COMMENTS None.

SIGNATURE Clare Conley

2567 Ca

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002265
06/26/2009

11:58:32

154

Standard Lot#: 801002I
Last Changed By: S. CIONE #50616

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:59:52
Diagnostic Test	Pass	12:00:19
Air Blank	0.000	12:00:47
0.100 Cal Check	0.101	12:01:05
Air Blank	0.000	12:01:35
Subject Test	0.000	12:01:51
Air Blank	ABT*	12:01:57
Air Blank	0.000	12:02:25

*Sequence Aborted

Not a Successfully Completed Test Sequence

3587CA

BY

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002265
06/26/2009

12:02:30

Standard Lot#: 801002I
Last Changed By: S. CIONE #50616

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	12:02:57
Air Blank	IPS**	12:03:05

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

Pursuant to section 28-1527, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

4587
en

PH

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002265
06/26/2009 12:03:09

Standard Lot#: 801002I
Last Changed By: S. CIONE #50616

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	12:03:56
Diagnostic Test	Pass	12:04:23
Air Blank	0.000	12:04:50
0.100 Cal Check	0.101	12:05:09
Air Blank	0.000	12:05:39
Subject Test	0.074*	12:08:47
Air Blank	0.000	12:09:21
Air Blank	0.000	12:09:48
Subject Test	XXX**	12:10:13
Air Blank	0.000	12:10:47

*Deficient Sample
**Invalid Sample
Begin new deprivation period

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

587
en

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-002265
06/26/2009 12:11:03

TSY

Standard Lot#: 801002I
Last Changed By: S. CIONE #50616

QAS: C. CONLEY #34517
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	12:11:30
Cal Check	0.099	12:11:49
Air Blank	0.000	12:12:18

Clare Conley

Operator's Signature

6 of 7
cu

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

BSG
7-1-09

QA Specialist Clare Conley AGENCY Tucson Police Department
DATE 6-26-09 TIME 1157
Intoxilyzer Serial # 80-002265 Location ODD

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 101 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Clare Conley

DPS Form Exh G-2 (Rev. 05-01)

7 of 7
cm

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-002265
06/26/2009 11:57:08

TS4

Standard Lot#: 801002I
Last Changed By: S. CIONE #50616

QAS: C. CONLEY #34517
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	11:57:35
Cal Check	0.101	11:57:54
Air Blank	0.000	11:58:24

Clare Conley
Operator's Signature