



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002265 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

SC
6/7/10

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 6.3.2010 Cal 1 0.100 AC Cal 2 0.099 AC

Between the dates of 5.6.2010 and 6.3.2010 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on .
The instrument was returned to service on .

Terry S. Gallegos

Terry S. Gallegos


6-7-2010

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 7, day of June, 2010.

A.M. Hernandez

Notary Public

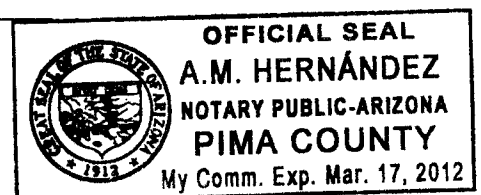

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1 of 2
TS4
SC
6/7/10

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department
DATE 6.3.2010 TIME 1605
INTOXILYZER SERIAL # 80-002265 LOCATION ODD

(☒) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (☒) 1. Clock time check.
(☒) 2. Date check.

OPERATIONAL TESTS

- (☒) 1. Alcohol-free subject test result 0. 000 AC.
(☒) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed ✓✓
(☒) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed ✓
Deficient sample printed. ✓
(☒) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 099 AC

Instrument is operating properly and accurately. YES ☒ NO ☐

COMMENTS

SIGNATURE

Terry Gallegos

DPS Form Exh G-4 (Rev. 05-01)

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002265
06/03/2010

16:07:29

Standard Lot#: 642845
Last Changed By: C. STRAUB #50772

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Subject: TEST, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	16:08:22
Diagnostic Test	Pass	16:08:50
Air Blank	0.000	16:09:17
0.100 Cal Check	0.099	16:09:36
Air Blank	0.000	16:10:05
Subject Test	0.000	16:10:29
Air Blank	ABT*	16:10:42
Air Blank	0.000	16:11:11

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002265
06/03/2010

16:11:27

Standard Lot#: 642845
Last Changed By: C. STRAUB #50772

QAS: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Subject: TEST, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	16:12:14
Air Blank	IPS**	16:12:22

*RFI Detect

**Improper Sample

3087
174

SC

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-002265

06/03/2010

16:12:30

Standard Lot#: 642845

Last Changed By: C. STRAUB

#50772

QAS: T. GALLEGOS

#28847

TUCSON PD CRIME LAB

Operator: T. GALLEGOS

#28847

TUCSON PD CRIME LAB

Subject: TEST, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	16:13:30
Diagnostic Test	Pass	16:13:57
Air Blank	0.000	16:14:24
0.100 Cal Check	0.100	16:14:43
Air Blank	0.000	16:15:13
Subject Test	0.094*	16:18:21
Air Blank	0.000	16:18:53
Air Blank	0.000	16:19:21
Subject Test	XXX**	16:19:49
Air Blank	0.000	16:20:23

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-002265

06/03/2010

16:20:41

Standard Lot#: 642845

Last Changed By: C. STRAUB #50772

QAS: T. GALLEGOS #28847

TUCSON PD CRIME LAB

Test	g/210L	Time
Air Blank	0.000	16:21:09
Cal Check	0.099	16:21:27
Air Blank	0.000	16:21:57



Operator's Signature

5087
TSH
SC

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6087
154
SC
6/7/10

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department
DATE 6.3.2010 TIME 1605
Intoxilyzer Serial # 80-002265 Location ODD

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C \pm 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Terry Gallegos

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-002265

06/03/2010

16:05:44


Standard Lot#: 642845

Last Changed By: C. STRAUB #50772

QAS: T. GALLEGOS #28847

TUCSON PD CRIME LAB

Test	g/210L	Time
Air Blank	0.000	16:06:12
Cal Check	0.100	16:06:30
Air Blank	0.000	16:07:00

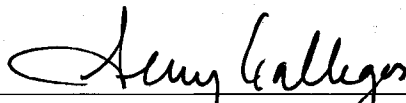

Operator's Signature