



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
270 S. Stone Avenue Tucson, AZ 85701-1917

By  
9-2-2010

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002265 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 08/26/10 Cal 1 0.099 AC Cal 2 0.098 AC

Between the dates of 07/29/10 and 08/26/10 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.  
The instrument was returned to service on N/A.

Cheryl K. Straub

Cheryl K. Straub


08/31/10

ADPS QAS Permit # 33408

SUBSCRIBED AND SWORN to before me this 31, day of August, 2010.

A.M. Hernandez

Notary Public

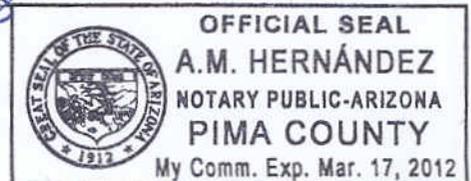

\*Includes documents specified in A.A.C. R13-10-104.A.5

177  
CS

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

TS4  
9-2-2010

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Cheryl STRAUB AGENCY Tucson Police Department  
DATE 08/26/10 TIME 09:16  
INTOXILYZER SERIAL # 80-002265 LOCATION TEAM 5

( ✓ ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- ( ✓ ) 1. Clock time check.
- ( ✓ ) 2. Date check.

OPERATIONAL TESTS

- ( ✓ ) 1. Alcohol-free subject test result 0. 000 AC.
- ( ✓ ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- ( ✓ ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- ( ✓ ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 098 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE

Cheryl Straub

297  
CS

B4

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-002265  
08/26/2010 09:18:16

---

Standard Lot#: 642845  
Last Changed By: C. STRAUB #50772

---

QAS: C. STRAUB #50772  
TUCSON PD

---

Operator: C. STRAUB #50772  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	09:19:11
Diagnostic Test	Pass	09:19:38
Air Blank	0.000	09:20:06
0.100 Cal Check	0.100	09:20:25
Air Blank	0.000	09:20:54
Subject Test	0.000	09:21:11
Air Blank	ABT*	09:21:18
Air Blank	0.000	09:21:46

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

377  
CS

TS4

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-002265  
08/26/2010 09:21:51

Standard Lot#: 642845  
Last Changed By: C. STRAUB #50772

QAS: C. STRAUB #50772  
TUCSON PD

Operator: C. STRAUB #50772  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	09:22:28
Air Blank	IPS**	09:22:36

\*RFI Detect  
\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

497  
CS

B4

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-002265  
08/26/2010 09:22:41

Standard Lot#: 642845  
Last Changed By: C. STRAUB #50772

QAS: C. STRAUB #50772  
TUCSON PD

Operator: C. STRAUB #50772  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:23:33
Diagnostic Test	Pass	09:24:00
Air Blank	0.000	09:24:28
0.100 Cal Check	0.099	09:24:46
Air Blank	0.000	09:25:16
Subject Test	0.159*	09:28:24
Air Blank	0.000	09:28:57
Air Blank	0.000	09:29:25
Subject Test	XXX**	09:29:56
Air Blank	0.000	09:30:32

\*Deficient Sample  
\*\*Invalid Sample  
Begin new deprivation period

### Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

587  
CS


TS4

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-002265  
08/26/2010 09:33:08

-----  
Standard Lot#: 642845  
Last Changed By: C. STRAUB #50772  
-----

QAS: C. STRAUB #50772  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	09:33:36
Cal Check	0.098	09:33:55
Air Blank	0.000	09:34:24

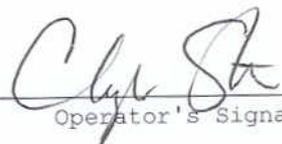  
-----  
Operator's Signature

697  
CSTHIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

B4

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000


9-2-2010

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Cheryl Straub AGENCY Tucson Police Department  
 DATE 08/26/10 TIME 09:16  
 Intoxilyzer Serial # 80-002265 Location TEAM 5

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
 OR  
 Pour a standard alcohol concentration solution N/A AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 099 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

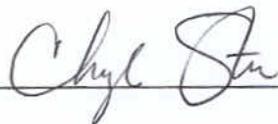

DPS Form Exh G-2 (Rev. 05-01)

797  
CS

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number: 80-002265

08/26/2010

09:16:40

BS4


-----  
Standard Lot#: 642845

Last Changed By: C. STRAUB #50772

-----  
QAS: C. STRAUB #50772

TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	09:17:07
Cal Check	0.099	09:17:26
Air Blank	0.000	09:17:56

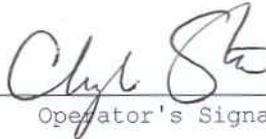

Operator's Signature