



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
270 S. Stone Avenue Tucson, AZ 85701-1917

*TS4  
9.27.2010*

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002265 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 09/23/10 Cal 1 0.099 AC Cal 2 0.098 AC

Between the dates of 08/26/10 and 09/23/10 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.  
The instrument was returned to service on N/A.

*Cheryl K. Straub*

Cheryl K. Straub


09/27/10

ADPS QAS Permit # 33408

SUBSCRIBED AND SWORN to before me this 27, day of September, 20 10.

*A.M. Hernandez*

Notary Public

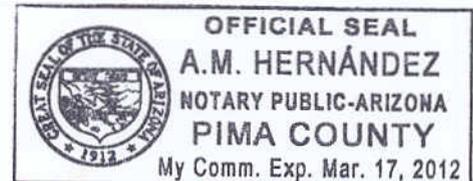

\*Includes documents specified in A.A.C. R13-10-104.A.5

197  
CS

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000  
STANDARD QUALITY ASSURANCE PROCEDURE

TSy  
9-27-2010

QA SPECIALIST Cheryl STRAUB AGENCY Tucson Police Department

DATE 09/23/10 TIME 08:37

INTOXILYZER SERIAL # 80-002265 LOCATION TEAM 5

( ✓ ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- ( ✓ ) 1. Clock time check.
- ( ✓ ) 2. Date check.

OPERATIONAL TESTS

- ( ✓ ) 1. Alcohol-free subject test result 0. 000 AC.
- ( ✓ ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- ( ✓ ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- ( ✓ ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 098 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE Cheryl Straub

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

277  
CS

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-002265  
09/23/2010

08:39:12

TS4

Standard Lot#: 642845  
Last Changed By: C. STRAUB #50772

QAS: C. STRAUB #50772  
TUCSON PD

Operator: C. STRAUB #50772  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:40:16
Diagnostic Test	Pass	08:40:44
Air Blank	0.000	08:41:11
0.100 Cal Check	0.100	08:41:30
Air Blank	0.000	08:41:59
Subject Test	0.000	08:42:18
Air Blank	ABT*	08:42:24
Air Blank	0.000	08:42:53

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

377  
CS

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-002265  
09/23/2010

08:43:02

BY

Standard Lot#: 642845  
Last Changed By: C. STRAUB . #50772

QAS: C. STRAUB #50772  
TUCSON PD

Operator: C. STRAUB #50772  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	08:43:38
Air Blank	IPS**	08:43:49

\*RFI Detect  
\*\*Improper Sample

### Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

437  
CS

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-002265  
09/23/2010

08:43:58

BSY

Standard Lot#: 642845  
Last Changed By: C. STRAUB. #50772

QAS: C. STRAUB #50772  
TUCSON PD

Operator: C. STRAUB #50772  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:44:50
Diagnostic Test	Pass	08:45:17
Air Blank	0.000	08:45:44
0.100 Cal Check	0.100	08:46:03
Air Blank	0.000	08:46:32
Subject Test	0.228*	08:49:41
Air Blank	0.000	08:50:15
Air Blank	0.000	08:50:43
Subject Test	XXX**	08:51:09
Air Blank	0.000	08:51:44

\*Deficient Sample  
\*\*Invalid Sample  
Begin new deprivation period

**Not a Successfully Completed Test Sequence**

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

507  
CS


INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-002265  
09/23/2010 08:51:59

184

-----  
Standard Lot#: 642845  
Last Changed By: C. STRAUB #50772  
-----

QAS: C. STRAUB #50772  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	08:52:26
Cal Check	0.098	08:52:45
Air Blank	0.000	08:53:14

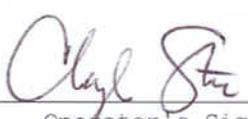

Operator's Signature

EXHIBIT G-2

677  
CS

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000  
STANDARD CALIBRATION CHECK PROCEDURE

By  
9-27-2010

QA Specialist Cheryl STRAUB AGENCY Tucson Police Department  
DATE 09/23/10 TIME 08:37  
Intoxilyzer Serial # 90-002265 Location TEAMS

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution N/A AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 099 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Cheryl Straub

DPS Form Exh G-2 (Rev. 05-01)


797  
CS

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-002265  
09/23/2010 08:37:27

104

-----  
Standard Lot#: 642845  
Last Changed By: C. STRAUB #50772  
-----  
QAS: C. STRAUB #50772  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	08:37:55
Cal Check	0.099	08:38:13
Air Blank	0.000	08:38:43

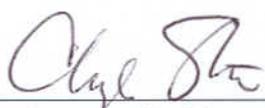

Operator's Signature