



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002265 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

TSY
4.13.2011

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 04/07/11 Cal 1 0.100 AC Cal 2 0.099 AC

Between the dates of 03/10/2011 and 04/07/2011 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.
The instrument was returned to service on N/A.

Cheryl K. Straub

Cheryl K. Straub


04/07/2011

ADPS QAS Permit # 33408

SUBSCRIBED AND SWORN to before me this 7th day of April, 2011.

A.M. Hernandez

Notary Public

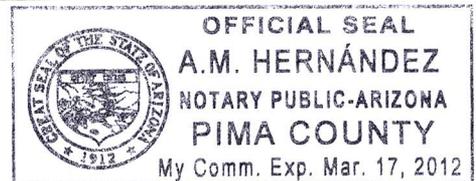

*Includes documents specified in A.A.C. R13-10-104.A.5

187
CS

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

TSY
4.13.2011

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Cheryl STRAUB AGENCY Tucson Police Department
DATE 04/07/2011 TIME 09:19
INTOXILYZER SERIAL # 80-002265 LOCATION TEAM 5

() 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- () 1. Clock time check.
- () 2. Date check.

OPERATIONAL TESTS

- () 1. Alcohol-free subject test result 0. 000 AC.
- () 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- () 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- () 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 099 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Cheryl Straub

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

287
CS

134

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-002265
04/07/2011 09:24:10

Standard Lot#: 20910100A1
Last Changed By: C. STRAUB #50772

OAS: C. STRAUB #50772
TUCSON PD

Operator: C. STRAUB #50772
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:25:09
Diagnostic Test	Pass	09:25:38
Air Blank	0.000	09:26:05
0.100 Cal Check	0.100	09:26:24
Air Blank	0.000	09:26:53
Subject Test	0.000	09:27:22
Air Blank	ABT*	09:27:29
Air Blank	0.000	09:27:57

*Sequence Aborted

Not a Successfully
Completed
Test Sequence

377
CS

TS4

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-002265
04/07/2011 09:28:05

Standard Lot#: 20910100A1
Last Charged By: C. STRAUB #50772

OAS: C. STRAUB #50772
TUCSON PD

Operator: C. STRAUB #50772
TUCSON PD

Subject: SOAP, SOAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	09:28:38
Air Blank	IPS**	09:28:47

*RFI Detect
**Improper Sample

Not a Successfully
Completed
Test Sequence

497
CS

104

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-002265
04/07/2011 09:29:00

Standard Lot#: 20910100A1
Last Changed By: C. STRAUB #50772

QAS: C. STRAUB #50772
TUCSON PD

Operator: C. STRAUB #50772
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:29:49
Diagnostic Test	Pass	09:30:18
Air Blank	0.000	09:30:45
0.100 Cal Check	0.100	09:31:04
Air Blank	0.000	09:31:33
Subject Test	0.238*	09:34:42
Air Blank	0.000	09:35:15
Air Blank	0.000	09:35:43
Subject Test	XXX**	09:36:09
Air Blank	0.000	09:36:43

*Deficient Sample
**Invalid Sample
Begin new deprivation period

Not a Successfully
Completed
Test Sequence

577
CS

754

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-002265
04/07/2011 09:37:52

Standard Lot#: 20910100A1
Last Changed By: C. STRAUB #50772

QAS: C. STRAUB #50772
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:38:20
Cal Check	0.099	09:38:38
Air Blank	0.000	09:39:08

C. Straub
Operator's Signature

687
CS

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

T54
4.13.2011

QA Specialist Cheryl Straub AGENCY Tucson Police Department
DATE 04/07/11 TIME 09:19
Intoxilyzer Serial # 80-002265 Location TEAMS

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution N/A AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Cheryl Straub

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

797
08

154

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-002265
04/07/2011 09:19:51

Standard Lot#: 20910100A1
Last Changed By: C. STRAUB #50772

QAS: C. STRAUB #50772
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:20:19
Cal Check	0.100	09:20:38
Air Blank	0.000	09:21:07

C. Straub

Operator's Signature