



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

TS4
2.23.2012

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002265 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 2/2/12 Cal 1 0.100 AC Cal 2 0.100 AC

Between the dates of 1/5/12 and 2/2/12 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.
The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione


2/10/12

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 10th, day of February, 2012.

Elma Escarcega

Notary Public

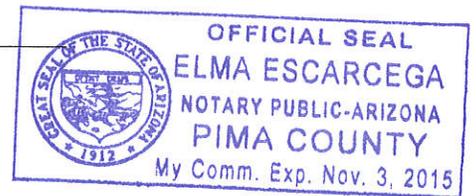

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

TS4
2-23-2012

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 2/2/12 TIME 1326

INTOXILYZER SERIAL # 80-002265 LOCATION DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 100 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Sarah Cione

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002265
02/02/2012

13:28:59

2 of 7sc

TS4

Standard Lot#: 20910100A1
Last Changed By: C. STRAUB #50772

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	13:30:01
Diagnostic Test	Pass	13:30:29
Air Blank	0.000	13:30:56
0.100 Cal Check	0.100	13:31:15
Air Blank	0.000	13:31:45
Subject Test	0.000	13:32:49
Air Blank	0.000	13:33:16
Five Minute Wait	ABT*	
Air Blank	0.000	13:33:51

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002265
02/02/2012 13:34:56

3 of 7 sc

TJ

Standard Lot#: 20910100A1
Last Changed By: C. STRAUB #50772

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	13:35:38
Diagnostic Test	Pass	13:36:05
Air Blank	0.000	13:36:33
0.100 Cal Check	RFI*	13:36:48
Air Blank	IPS**	13:36:56

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

4 of 7 sc

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002265
02/02/2012 13:38:24

TS4

Standard Lot#: 20910100A1
Last Changed By: C. STRAUB #50772

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	13:39:20
Diagnostic Test	Pass	13:39:47
Air Blank	0.000	13:40:14
0.100 Cal Check	0.100	13:40:33
Air Blank	0.000	13:41:03
Subject Test	0.342*	13:44:11
Air Blank	0.000	13:44:44
Air Blank	0.000	13:45:11
Subject Test	XXX**	13:45:42
Air Blank	0.000	13:46:16

*Deficient Sample
**Invalid Sample
Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-002265
02/02/2012 13:51:46

5 of 7sc

Standard Lot#: 20910100A1
Last Changed By: C. STRAUB #50772

TS4

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	13:52:14
Cal Check	0.100	13:52:33
Air Blank	0.000	13:53:03



Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

6 of 7 sc

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

TS4
2.23.2012

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 2/2/12 TIME 1326
Intoxilyzer Serial # 80-002265 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0.100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0.100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-002265
02/02/2012 13:26:52

7 of 7 sc

TS4

Standard Lot#: 20910100A1
Last Changed By: C. STRAUB #50772

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	13:27:20
Cal Check	0.100	13:27:38
Air Blank	0.000	13:28:08



Operator's Signature