



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

TS4
9-20-2012

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002265 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 9/13/12 Cal 1 0.098 AC Cal 2 0.098 AC

Between the dates of 8/16/12 and 9/13/12 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.
The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione

9/17/12


ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 17th, day of September, 2012.

Adma Escarcega

Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

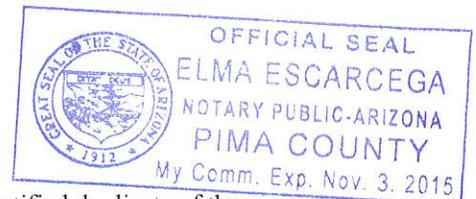

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

TS4
9.20.2012

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 9/13/12 TIME 1010

INTOXILYZER SERIAL # 80-002265 LOCATION DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 098 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Sarah Cione

2 of 9sc

TS4

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002265
09/13/2012

10:12:27

Standard Lot#: 04911100A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|------------------|--------|----------|
| Air Blank | 0.000 | 10:13:14 |
| Diagnostic Test | Pass | 10:13:41 |
| Air Blank | 0.000 | 10:14:09 |
| 0.100 Cal Check | 0.099 | 10:14:27 |
| Air Blank | 0.000 | 10:14:57 |
| Subject Test | 0.000 | 10:15:24 |
| Air Blank | 0.000 | 10:15:51 |
| Five Minute Wait | ABT* | |
| Air Blank | 0.000 | 10:16:26 |

*Sequence Aborted

Not a Successfully Completed Test Sequence

3 of 9 sc

TS4

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002265
09/13/2012

10:18:25

Standard Lot#: 04911100A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 10:19:12 |
| Diagnostic Test | Pass | 10:19:40 |
| Air Blank | RFI* | 10:19:47 |
| Air Blank | IPS** | 10:19:54 |

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

4 of 9 sc

134

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002265
09/13/2012 10:25:02

Standard Lot#: 04911100A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 10:25:45 |
| Diagnostic Test | Pass | 10:26:12 |
| Air Blank | 0.000 | 10:26:40 |
| 0.100 Cal Check | 0.099 | 10:26:59 |
| Air Blank | 0.000 | 10:27:29 |
| Subject Test | 0.119* | 10:30:37 |
| Air Blank | 0.000 | 10:31:10 |
| Air Blank | 0.000 | 10:31:37 |
| Subject Test | XXX** | 10:32:07 |
| Air Blank | 0.000 | 10:32:43 |

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

5 of 9 sc


INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-002265
09/13/2012 10:34:09

TS4

Standard Lot#: 04911100A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

| Test | g/210L | Time |
|-----------|--------|----------|
| Air Blank | 0.000 | 10:34:36 |
| Cal Check | 0.098 | 10:34:55 |
| Air Blank | 0.000 | 10:35:26 |

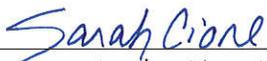

Operator's Signature

6 of 9sc

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000
STANDARD CALIBRATION CHECK PROCEDURE

TS4
9.20.2012

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 9/13/12 TIME 1010
Intoxilyzer Serial # 80-002265 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 098 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione


INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-002265
09/13/2012 10:10:42

7 of 9 sc
134

Standard Lot#: 04911100A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

| Test | g/210L | Time |
|-----------|--------|----------|
| Air Blank | 0.000 | 10:11:09 |
| Cal Check | 0.098 | 10:11:28 |
| Air Blank | 0.000 | 10:11:58 |

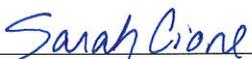

Operator's Signature



TUCSON POLICE DEPARTMENT
CRIME LABORATORY

8 of 9_{sc}

T54
9.20.2012

INTOXILYZER SERVICE RECORD
FOR
CALIBRATION GAS TANK CHANGE

Date: 9/13/12
QAS: Sarah Cione
Intoxilyzer Serial Number: 80-002265
New Lot Number: 09312100A1
New Expiration Date: 6/1/2014

- (✓) 1. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 2. Remove the old gas tank and install the new one.
- (✓) 3. Hit the "ESC" key twice and swipe the QAS card.
- (✓) 4. Go to the "M-Maintenance" menu and select "S-Standard Change."
Select "D" for dry gas and enter the following information:
 - Enter Target Value: 0.100
 - Enter Control Lot #: From the tank label
 - Enter Exp Date: From the tank label (mm/dd/yyyy)
 - Edit info Y/N? Y or N
 - Update Std Info? Y ENTER
 - "Please wait saving..."
- (✓) 5. Hit the "ESC" key twice.
- (✓) 6. Go to the "M-Maintenance" menu and select "T-Control Testing."
Select "D" for dry control test. After selection is made press ENTER.
- (✓) 7. Air blank completed.
- (✓) 8. Calibration check completed. Test results 0. 099 AC.
- (✓) 9. Air blank completed.
- (✓) 10. Remove and sign printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

9 of 9 sc


154

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-002265
09/13/2012 10:37:57

Standard Lot#: 09312100A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

| Test | g/210L | Time |
|-----------|--------|----------|
| Air Blank | 0.000 | 10:38:25 |
| Cal Check | 0.099 | 10:38:44 |
| Air Blank | 0.000 | 10:39:14 |

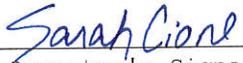

Operator's Signature