



CERTIFICATION by CUSTODIAN of RECORDS  
 Tucson Police Department Crime Laboratory  
 1306 W. Miracle Mile, Tucson, AZ 85705

*TSK*  
*11-19-2012*

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002265 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 11/8/12 Cal 1 0.099 AC Cal 2 0.099 AC

Between the dates of 10/11/12 and 11/8/12 the instrument was operating properly and accurately.


The instrument was taken out of service for repair on -.  
 The instrument was returned to service on -.

*Sarah Cione*  
 Sarah L. Cione

11/14/12  
 ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 14<sup>th</sup>, day of November, 2012.

*Elma Escarcega*  
 Notary Public

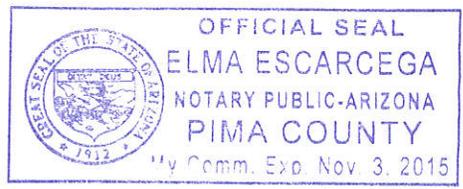

\*Includes documents specified in A.A.C. R13-10-104.A.5

1 of 7 sc

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

By  
11.19.2012

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 11/8/12 TIME 1008

INTOXILYZER SERIAL # 80-002265 LOCATION DUI HQ

( ✓ ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- ( ✓ ) 1. Clock time check.
- ( ✓ ) 2. Date check.

OPERATIONAL TESTS

- ( ✓ ) 1. Alcohol-free subject test result 0. 000 AC.
- ( ✓ ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- ( ✓ ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- ( ✓ ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 099 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

---



---



---



---

SIGNATURE Sarah Cione

2 of 7<sub>5</sub>

15h

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-002265  
11/08/2012 10:11:33

---

Standard Lot#: 09312100A1  
Last Changed By: S. CIONE #50616

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1980  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	10:12:20
Diagnostic Test	Pass	10:12:47
Air Blank	0.000	10:13:14
0.100 Cal Check	0.100	10:13:33
Air Blank	0.000	10:14:03
Subject Test	0.000	10:14:43
Air Blank	0.000	10:15:11
Five Minute Wait	ABT*	
Air Blank	0.000	10:16:10

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

3 of 7sc  
BH

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-002265  
11/08/2012 10:17:20

---

Standard Lot#: 09312100A1  
Last Changed By: S. CIONE #50616

---

QAS: S. CIONE #50616  
TUCSON PD

---

Operator: S. CIONE #50616  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1980  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	10:18:21
Diagnostic Test	Pass	10:18:48
Air Blank	0.000	10:19:15
0.100 Cal Check	0.100	10:19:34
Air Blank	RFI*	10:19:59
Air Blank	IPS**	10:20:06

\*RFI Detect  
\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

4 of 7<sub>sc</sub>

134

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number : 80-002265  
11/08/2012 10:24:36

Standard Lot#: 09312100A1  
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1980  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:25:36
Diagnostic Test	Pass	10:26:03
Air Blank	0.000	10:26:30
0.100 Cal Check	0.100	10:26:49
Air Blank	0.000	10:27:19
Subject Test	0.155*	10:30:28
Air Blank	0.000	10:31:00
Air Blank	0.000	10:31:27
Subject Test	XXX**	10:31:51
Air Blank	0.000	10:32:26

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-002265  
11/08/2012 10:33:17

5 of 7 sc  
Bl

-----  
Standard Lot#: 09312100A1  
Last Changed By: S. CIONE #50616  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	10:33:44
Cal Check	0.099	10:34:03
Air Blank	0.000	10:34:33

  
-----  
Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6 of 7 sc  
TS4  
11.19.2012

QA Specialist Sarah Cione AGENCY Tucson Police Department  
DATE 11/8/12 TIME 1008  
Intoxilyzer Serial # 80-002265 Location DUI HQ

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 099 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-002265  
11/08/2012 10:08:57

7 of 7 sec  
124

-----  
Standard Lot#: 09312100A1  
Last Changed By: S. CIONE #50616  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	10:09:24
Cal Check	0.099	10:09:43
Air Blank	0.000	10:10:13

  
-----  
Operator's Signature