



CERTIFICATION by CUSTODIAN of RECORDS
 Tucson Police Department Crime Laboratory
 1306 W. Miracle Mile, Tucson, AZ 85705

TS4
 1-7-2013

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002265 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 12/7/12 Cal 1 0.100 AC Cal 2 0.100 AC

Between the dates of 11/8/12 and 12/7/12 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.
 The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione


12/17/12

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 17th, day of December, 2012.

Elma Escarcega

Notary Public

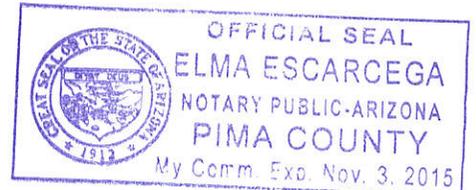

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

1 of 7 sc

ARIZONA DEPARTMENT OF PUBLIC SAFETY

TS4

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

1.7.2013

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 12/7/12 TIME 947

INTOXILYZER SERIAL # 80-002265 LOCATION Crime Lab

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

(✓) 1. Clock time check.

(✓) 2. Date check.

OPERATIONAL TESTS

(✓) 1. Alcohol-free subject test result 0. 000 AC.

(✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed

(✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.

(✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 100 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

2 of 7_{sc}

TS4

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002265
12/07/2012 09:48:53

Standard Lot#: 09312100A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:49:35
Diagnostic Test	Pass	09:50:03
Air Blank	0.000	09:50:30
0.100 Cal Check	0.101	09:50:49
Air Blank	0.000	09:51:19
Subject Test	0.000	09:51:45
Air Blank	0.000	09:52:13
Five Minute Wait	ABT*	
Air Blank	0.000	09:52:57

*Sequence Aborted

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

3 of 7 sc
154

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002265
12/07/2012 09:53:46

Standard Lot#: 09312100A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	09:54:09
Air Blank	IPS**	09:54:17

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

4 of 7 sc

D4

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002265
12/07/2012

09:54:21

Standard Lot#: 09312100A1 #50616
Last Changed By: S. CIONE

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:55:05
Diagnostic Test	Pass	09:55:33
Air Blank	0.000	09:56:00
0.100 Cal Check	0.101	09:56:19
Air Blank	0.000	09:56:49
Subject Test	0.082*	09:59:57
Air Blank	0.000	10:00:29
Air Blank	0.000	10:00:56
Subject Test	XXX**	10:01:40
Air Blank	0.000	10:02:14

*Deficient Sample
 **Invalid Sample
 Begin new deprivation period

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

5 of 7_{sc}

B4

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-002265
12/07/2012 10:02:41

Standard Lot#: 09312100A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	10:03:09
Cal Check	0.100	10:03:28
Air Blank	0.000	10:03:58

Sarah Cione
Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

6 of 7_{sc}

ARIZONA DEPARTMENT OF PUBLIC SAFETY

TS4
1.7.2013

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department

DATE 12/7/12 TIME 947

Intoxilyzer Serial # 80-002265 Location Crime Lab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-002265
12/07/2012 09:47:23

7 of 7 sc
TS4

Standard Lot#: 09312100A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:47:50
Cal Check	0.100	09:48:09
Air Blank	0.000	09:48:39

Sarah Cione

Operator's Signature