



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
1306 W. Miracle Mile, Tucson, AZ 85705

TS4  
2.19.2013

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002265 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 2/5/13 Cal 1 0.099 AC Cal 2 0.099 (2/11/13) AC

Between the dates of 1/4/13 and 2/5/13 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on 2/11/13.  
The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione


2/14/13

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 14th, day of February, 20 13.

Elma Escarcega

Notary Public

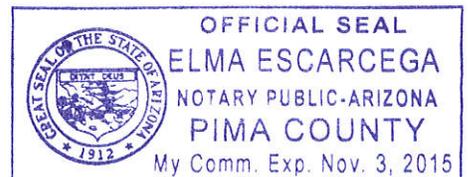

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

1 of 9 sc  
TBH  
2-19-2013

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 2/11/13 TIME 1050

INTOXILYZER SERIAL # 80-002265 LOCATION Crime Lab

( ✓ ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- ( ✓ ) 1. Clock time check.
- ( ✓ ) 2. Date check.

OPERATIONAL TESTS

- ( ✓ ) 1. Alcohol-free subject test result 0. 000 AC.
- ( ✓ ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- ( ✓ ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- ( ✓ ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 099 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS G-2 (Cal Check) performed on 2/5/13 (1 day past 31 days)  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

2 of 9 sc  
TS4

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-002265  
02/11/2013 10:50:17

-----  
Standard Lot#: 09312100A1  
Last Changed By: S. CIONE #50616  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Operator: S. CIONE #50616  
TUCSON PD  
-----

Subject: SQAP, SQAP  
DOB: 11/11/1980  
Sex: F Weight: 111  
-----

15 Minute Deprivation Period? Yes  
-----

Test	g/210L	Time
Air Blank	0.000	10:51:01
Diagnostic Test	Pass	10:51:30
Air Blank	0.000	10:51:58
0.100 Cal Check	0.100	10:52:16
Air Blank	0.000	10:52:46
Subject Test	0.000	10:53:17
Air Blank	0.000	10:53:45
Five Minute Wait	ABT*	
Air Blank	0.000	10:54:20

\*Sequence Aborted

Not a Successfully  
Completed  
Test Sequence

3 of 9 sc  
B4

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-002265  
02/11/2013 10:54:36

-----  
Standard Lot#: 09312100A1  
Last Changed By: S. CIONE #50616  
-----

OAS: S. CIONE #50616  
TUCSON PD  
-----

Operator: S. CIONE #50616  
TUCSON PD  
-----

Subject: SOAP, SOAP  
DOB: 11/11/1111  
Sex: F Weight: 111  
-----

15 Minute Deprivation Period? Yes  
-----

Test	g/210L	Time
Air Blank	0.000	10:55:20
Diagnostic Test	Pass	10:55:49
Air Blank	RFI*	10:55:59
Air Blank	IPS**	10:56:07

\*RFI Detect  
\*\*Improper Sample

Not a Successfully  
Completed  
Test Sequence

4 of 9 sc  
456  
2.19.2013

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-002265  
02/11/2013 10:56:13

-----  
Standard Lot#: 09312100A1  
Last Changed By: S. CIONE #50616  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Operator: S. CIONE #50616  
TUCSON PD  
-----

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111  
-----

15 Minute Deprivation Period? Yes  
-----

Test	g/210L	Time
Air Blank	0.000	10:56:57
Diagnostic Test	Pass	10:57:26
Air Blank	0.000	10:57:53
0.100 Cal Check	0.099	10:58:12
Air Blank	0.000	10:58:42
Subject Test	0.200*	11:01:50
Air Blank	0.000	11:02:23
Air Blank	0.000	11:02:50
Subject Test	XXX**	11:03:15
Air Blank	0.000	11:03:50

\*Deficient Sample  
\*\*Invalid Sample  
Begin new deprivation period

Not a Successfully  
Completed  
Test Sequence

5 of 9 sc  
TS4



INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-002265  
02/11/2013 11:04:19

-----  
Standard Lot#: 09312100A1  
Last Changed By: S. CIONE #50616  
-----

QAS: S. CIONE #50616  
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	11:04:46
Cal Check	0.099	11:05:05
Air Blank	0.000	11:05:36

*Sarah Cione*  
-----  
Operator's Signature



EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

6 of 9 sc  
TBH  
2.19.2013

ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000  
STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department  
DATE 2/5/13 TIME 1055  
Intoxilyzer Serial # 80-002265 Location Crime Lab

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 099 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

7 of 9 SC  
TS4

INTOXILYZER 8000  
Location: TUCSON POLICE DEPT  
Serial Number: 80-002265  
02/05/2013 10:55:56

-----  
Standard Lot#: 09312100A1  
Last Changed By: S. CIONE #50616  
-----

OAS: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	10:56:23
Cal Check	0.099	10:56:42
Air Blank	0.000	10:57:12

-----  
*Sarah Cione*  
Operator's Signature  
-----



# Service Evaluation Form

8 of 9 SC  
TS4  
2-19-2013

This form MUST be completed and enclosed with instrument to be serviced.  
Failure to complete and return this form may cause delays in service.

(Note: please ship items in their original shipping container(s) or a similar protective box.)

1. Contact information: Customer Number 105075 (contact Customer Service)  
Name Terry Gallegos Phone: (520) 837-7677  
Fax: (520) 791-4642 Email: Terry.Gallegos@tucsonaz.gov

2. Bill to Address: Ship to Address:  
Finance / Accounts Payable Tucson Police Department Crime Lab  
Po Box 27450 1306 W. Miracle Mile  
Tucson, AZ 85726-7450 Tucson, AZ 85705

3. Serial Number: 80-002265 Instrument Model: I 8000

4. Detailed Description of Problem:  
Takes a long time to initially go to Ready Mode (at room temp). Also needs software update to 8105.48.

\*\*\*Hazardous Material Warning! - DO NOT return gas cylinders with instrument!\*\*\*

I Authorize Repairs Up To:  All  \$250  \$500  \$750  Other \$ \_\_\_\_\_  
Purchase Order Number (attach a copy of P.O. if applicable) 34770

Authorized By:

Terry Gallegos Toxicology Coordinator  
Name (Please Print) Title  
[Signature] 2-11-2013  
Signature Date

Ship to:  
CMI, Inc.  
Attn: Service Dept.  
316 East Ninth Street  
Owensboro, KY 42303

No, please send estimate before repairs are made.  
Note: An estimate will be faxed before performing any repairs and may cause delays in service.  
An evaluation fee (\$79.00 or actual costs) will apply to estimates that are not repaired.



TUCSON POLICE DEPARTMENT  
CRIME LABORATORY

9 of 9 sc

INTOXILYZER SERVICE RECORD *TS4*

2.19.2013

Date: 2/11/13

QAS: Sarah Cione

Intoxilyzer Serial Number: 80-002265

Reason For Service: At room temp, takes a long time to get to Ready Mode

Date Out of Service: 2/11/13

Calibration Check completed prior to removal from service

Date: 2/11/13

SQAP completed prior to removal from service

Date: 2/11/13

Corrective Action:

Sent to CMI, Inc.

downloaded before sending to CMI