



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

TSy
7-24-2013

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002265 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 7/18/13 Cal 1 0.101 AC Cal 2 0.100 AC

Between the dates of 6/20/13 and 7/18/13 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.
The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione


7/22/13

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 22nd, day of July, 2013.

Elma Escarcega

Notary Public

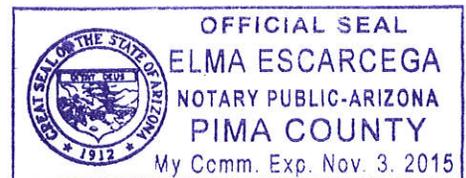

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1077
JC
TS4
7.24.2013

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 7/18/13 TIME 1009

INTOXILYZER SERIAL # 80-002265 LOCATION DUI HQ

() 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

() 1. Clock time check.

() 2. Date check.

OPERATIONAL TESTS

() 1. Alcohol-free subject test result 0. 000 AC.

() 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed

() 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.

() 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 100 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE

Sarah Cione

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002265
Core Version: 8105.48
07/18/2013

10:11:42

2 of 7_{sc}
B4

Standard Lot#: 09312100A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 10:12:28 |
| Diagnostic Test | Pass | 10:13:03 |
| Air Blank | 0.000 | 10:13:31 |
| 0.100 Cal Check | 0.101 | 10:13:51 |
| Air Blank | 0.000 | 10:14:21 |
| Subject Test | 0.000 | 10:14:40 |
| Air Blank | 0.000 | 10:15:07 |
| Wait | ABT* | |
| Air Blank | 0.000 | 10:15:43 |

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002265
Core Version: 8105.48
07/18/2013

10:19:35

3 of 7
TS4 SC

Standard Lot#: 09312100A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 10:20:17 |
| Diagnostic Test | Pass | 10:20:52 |
| Air Blank | 0.000 | 10:21:19 |
| 0.100 Cal Check | 0.100 | 10:21:40 |
| Air Blank | 0.000 | 10:22:10 |
| Subject Test | RFI* | 10:22:20 |
| Air Blank | IPS** | 10:22:28 |

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002265
Core Version: 8105.48
07/18/2013

10:29:46

4 of 7
sc
TS4

Standard Lot#: 09312100A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1980
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 10:30:41 |
| Diagnostic Test | Pass | 10:31:16 |
| Air Blank | 0.000 | 10:31:44 |
| 0.100 Cal Check | 0.100 | 10:32:04 |
| Air Blank | 0.000 | 10:32:34 |
| Subject Test | 0.133* | 10:35:42 |
| Air Blank | 0.000 | 10:36:13 |
| Air Blank | 0.000 | 10:36:41 |
| Subject Test | XXX** | 10:36:53 |
| Air Blank | 0.000 | 10:37:26 |

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002265
Core Version: 8105.48
07/18/2013 10:37:46

5 of 7
sc
TSy

Standard Lot#: 09312100A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 10:38:14 |
| 0.100 Cal Check | 0.100 | 10:38:35 |
| Air Blank | 0.000 | 10:39:05 |

Preventative Maintenance
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000
STANDARD CALIBRATION CHECK PROCEDURE

6 of 7
5c
Tdy
7.24.2013

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 7/18/13 TIME 1009
Intoxilyzer Serial # 80-002265 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 101 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002265
Core Version: 8105.48
07/18/2013 10:09:11

7 of 7 sc
TS4

Standard Lot#: 09312100A1
Last Changed By: S. CIONE #50616

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 10:09:40 |
| 0.100 Cal Check | 0.101 | 10:10:01 |
| Air Blank | 0.000 | 10:10:31 |

Preventative Maintenance
Performed -- Timer Reset