



CERTIFICATION by CUSTODIAN of RECORDS
 Tucson Police Department Crime Laboratory
 1306 W. Miracle Mile, Tucson, AZ 85705

SC
 5/19/14

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002265 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 2-7-2014 Cal 1 0.101 AC Cal 2 0.100 AC

Between the dates of 1-22-2014 and 2-7-2014 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on 2-7-2014.
 The instrument was returned to service on —.

Terry Gallegos

Terry S. Gallegos


2-13-2014

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 13th, day of February, 20 14.

Edna Escarcega

Notary Public

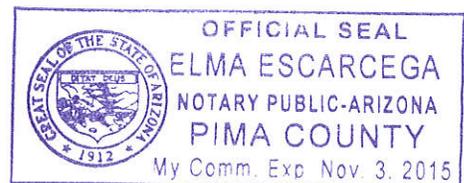

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

10f9
TS4
SC
5/19/14

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department

DATE 2-7-2014 TIME 1515

INTOXILYZER SERIAL # 80-002265 LOCATION Crime Lab

() 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- () 1. Clock time check.
- () 2. Date check.

OPERATIONAL TESTS

- () 1. Alcohol-free subject test result 0. 000 AC.
- () 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- () 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- () 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 100 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Terry Gallegos

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TPD UAPD
Serial Number: 80-002265
Core Version: 8105.48
02/07/2014

15:16:33

20f9
TS4
SL

Standard Lot#: 30412100A1
Last Changed By: S. CIONE #50616

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:17:23
Diagnostic Test	Pass	15:17:58
Air Blank	0.000	15:18:26
0.100 Cal Check	0.101	15:18:46
Air Blank	0.000	15:19:16
Subject Test	0.000	15:19:34
Air Blank	0.000	15:20:01
Wait	ABT*	
Air Blank	0.000	15:20:36

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TPD UAPD
Serial Number: 80-002265
Core Version: 8105.48
02/07/2014

15:20:44

3 of 9
TS4

SC

Standard Lot#: 30412100A1
Last Changed By: S. CIONE #50616

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:21:42
Diagnostic Test	Pass	15:22:17
Air Blank	0.000	15:22:45
0.100 Cal Check	RFI*	15:23:00
Air Blank	IPS**	15:23:09

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TPD UAPD
Serial Number: 80-002265
Core Version: 8105.48
02/07/2014

15:23:15

4 of 9
TSh

Standard Lot#: 30412100A1
Last Changed By: S. CIONE #50616

SL

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:24:04
Diagnostic Test	Pass	15:24:39
Air Blank	0.000	15:25:06
0.100 Cal Check	0.101	15:25:27
Air Blank	0.000	15:25:57
Subject Test	0.065*	15:29:04
Air Blank	0.000	15:29:36
Air Blank	0.000	15:30:04
Subject Test	XXX**	15:30:30
Air Blank	0.000	15:31:04

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TPD UAPD
Serial Number: 80-002265
Core Version: 8105.48
02/07/2014 15:31:31

Sof9
TS4
SC

Standard Lot#: 30412100A1
Last Changed By: S. CIONE #50616

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	15:32:02
0.100 Cal Check	0.100	15:32:23
Air Blank	0.000	15:32:53

Preventative Maintenance
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

60f9
754

sc
5/19/14

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department
DATE 2-7-2014 TIME 1515
Intoxilyzer Serial # 80-002265 Location Crime Lab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 101 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Terry Gallegos

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: TPD UAPD
Serial Number: 80-002265
Core Version: 8105.48
02/07/2014 15:15:04

70f9
TS4
SC

Standard Lot#: 30412100A1
Last Changed By: S. CIONE #50616

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD CRIME LAB

Test	g/210L	Time
Air Blank	0.000	15:15:32
0.100 Cal Check	0.101	15:15:53
Air Blank	0.000	15:16:23



TUCSON POLICE DEPARTMENT
CRIME LABORATORY

80f9
T56
SC

INTOXILYZER SERVICE RECORD

Date: 1-22-2014

QAS: Terry Gallegos

Intoxilyzer Serial Number: 80-002265

Reason For Service: Failed to detect RFI

Date Out of Service: 2-7-2014

Calibration Check completed prior to removal from service

Date: 2-7-2014

SQAP completed prior to removal from service

Date: 2-7-2014

Corrective Action:

Sent to CMI, Inc.



Service Evaluation Form

9 of 9
TS4

This form MUST be completed and enclosed with instrument to be serviced.
Failure to complete and return this form may cause delays in service.

SC

(Note: please ship items in their original shipping container(s) or a similar protective box.)

1. Contact information: Customer Number 105075 (contact Customer Service)
 Name Terry Gallegos Phone: (520) 837-7677
 Fax: (520) 791-4642 Email: Terry.Gallegos@tucsonaz.gov

2. Bill to Address: Ship to Address:
Finance / Accounts Payable Tucson Police Department Crime Lab
Po Box 27450 1306 W. Miracle Mile
Tucson, AZ 85726-7450 Tucson, AZ 85705

3. Serial Number: 8U-002265 Instrument Model: I8000

4. Detailed Description of Problem:
Failed to detect RFI

Hazardous Material Warning! - DO NOT return gas cylinders with instrument!

I Authorize Repairs Up To: All \$250 \$500 \$750 Other \$ _____
 Purchase Order Number (attach a copy of P.O. if applicable) 34770

Authorized By:

Terry Gallegos Toxicology Coordinator
 Name (Please Print) Title
Terry Gallegos 2.7.2014
 Signature Date

Ship to:

CMI, Inc.
 Attn: Service Dept.
 316 East Ninth Street
 Owensboro, KY 42303

No, please send estimate before repairs are made.
 Note: An estimate will be faxed before performing any repairs and may cause delays in service.
 An evaluation fee (\$79.00 or actual costs) will apply to estimates that are not repaired.

Pursuant to Section 28.1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.