



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
1306 W. Miracle Mile, Tucson, AZ 85705

SC  
5/19/14

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002265 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 4-2-2014 Cal 1 0.096 AC Cal 2 0.095 AC

Between the dates of N/A and N/A the instrument was operating properly and accurately.

The instrument was taken out of service for repair on —.  
The instrument was returned to service on 4-2-2014.

Terry S. Gallegos

Terry S. Gallegos


5-15-2014

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 15<sup>th</sup>, day of May, 2014.

Elma Escarcega

Notary Public

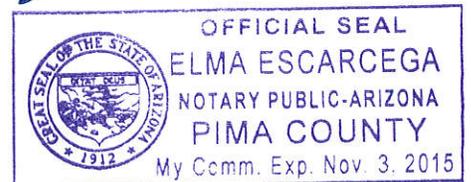

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

10f12  
TS4  
SC  
5/19/14

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department

DATE 4.2.2014 TIME 1805

INTOXILYZER SERIAL # 80-002265 LOCATION CrimeLab

( ✓ ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- ( ✓ ) 1. Clock time check.
- ( ✓ ) 2. Date check.

OPERATIONAL TESTS

- ( ✓ ) 1. Alcohol-free subject test result 0. 000 AC.
- ( ✓ ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- ( ✓ ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- ( ✓ ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 095 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

---



---



---



---

SIGNATURE

Terry Gallegos

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002265  
Core Version: 8105.48  
04/02/2014 18:06:42

20f12  
TS4

-----  
Standard Lot#: 30412100A1  
Last Changed By: T. GALLEGOS #28847  
-----

SC

-----  
QAS: T. GALLEGOS #28847  
TUCSON PD  
-----

-----  
Operator:  
-----

Test	g/210L	Time
Air Blank	0.000	18:07:10
Subject Test	ABT*	18:07:22
Air Blank	0.000	18:07:51

\*Sequence Aborted

\_\_\_\_\_  
Operator's Signature

Configure Start test = Card Swipe required  
& Custom test sequence  
TS4 28847

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002265  
Core Version: 8105.48  
04/02/2014

18:08:41

30f12  
TS4

5L

---

Standard Lot#: 30412100A1  
Last Changed By: T. GALLEGOS #28847

---

QAS: T. GALLEGOS #28847  
TUCSON PD

---

Operator: T. GALLEGOS #28847  
TUCSON PD

---

Subject: TEST, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	18:09:32
Diagnostic Test	Pass	18:10:07
Air Blank	0.000	18:10:35
0.100 Cal Check	0.095	18:10:56
Air Blank	0.000	18:11:25
Subject Test	0.000	18:11:50
Air Blank	0.000	18:12:18
Wait	ABT*	
Air Blank	0.000	18:12:52

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002265  
Core Version: 8105.48  
04/02/2014

18:12:55

40F12  
TS6

SC

---

Standard Lot#: 30412100A1  
Last Changed By: T. GALLEGOS #28847

---

QAS: T. GALLEGOS #28847  
TUCSON PD

---

Operator: T. GALLEGOS #28847  
TUCSON PD

---

Subject: TEST, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	18:13:47
Diagnostic Test	Pass	18:14:23
Air Blank	0.000	18:14:50
0.100 Cal Check	RFI*	18:15:08
Air Blank	IPS**	18:15:17

\*RFI Detect  
\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002265  
Core Version: 8105.48  
04/02/2014

18:15:21

50F12  
TS4

5c

---

Standard Lot#: 30412100A1  
Last Changed By: T. GALLEGOS #28847

---

QAS: T. GALLEGOS #28847  
TUCSON PD

---

Operator: T. GALLEGOS #28847  
TUCSON PD

---

Subject: TEST, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	18:16:11
Diagnostic Test	Pass	18:16:46
Air Blank	0.000	18:17:13
0.100 Cal Check	0.096	18:17:34
Air Blank	0.000	18:18:04
Subject Test	0.230*	18:21:11
Air Blank	0.000	18:21:45
Air Blank	0.000	18:22:12
Subject Test	XXX**	18:22:28
Air Blank	0.000	18:23:01

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002265  
Core Version: 8105.48  
04/02/2014 18:23:21

60F12  
TS4

-----  
Standard Lot#: 30412100A1  
Last Changed By: T. GALLEGOS #28847  
-----

SC

QAS: T. GALLEGOS #28847  
TUCSON PD  
-----

Operator: T. GALLEGOS #28847  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	18:23:51
0.100 Cal Check	0.095	18:24:12
Air Blank	0.000	18:24:42

Preventative Maintenance  
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

70F12  
TS4

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

SC  
5/19/14

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department  
DATE 4.2.2014 TIME 1805  
Intoxilyzer Serial # 80-002265 Location CrimeLab

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 096 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE 

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002265  
Core Version: 8105.48  
04/02/2014 18:05:17

8 of 12  
Dh

-----  
Standard Lot#: 30412100A1  
Last Changed By: T. GALLEGOS #28847  
-----

sc

QAS: T. GALLEGOS #28847  
TUCSON PD  
-----

Operator: T. GALLEGOS #28847  
TUCSON PD CRIME LAB  
-----

Test	g/210L	Time
Air Blank	0.000	18:05:45
0.100 Cal Check	0.096	18:06:05
Air Blank	0.000	18:06:35



TUCSON POLICE DEPARTMENT  
CRIME LABORATORY

9 of 12  
TS4

INTOXILYZER SERVICE RECORD  
FOR  
CALIBRATION GAS TANK CHANGE

SC  
5/19/14

Date: 4.2.2014  
QAS: Terry Gallegos 28847  
Intoxilyzer Serial Number: 80-002265  
New Lot Number: 30412100A1  
New Expiration Date: 12.1.2014

- (✓) 1. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 2. Remove the old gas tank and install the new one.
- (✓) 3. Hit the "ESC" key twice and swipe the QAS card.
- (✓) 4. Go to the "M-Maintenance" menu and select "S-Standard Change."  
Select "D" for dry gas and enter the following information:  
Enter Target Value: 0.100  
Enter Control Lot #: From the tank label  
Enter Exp Date: From the tank label (mm/dd/yyyy)  
Edit info Y/N? Y or N  
Update Std Info? Y ENTER  
"Please wait saving..."
- (✓) 5. Hit the "ESC" key twice.
- (✓) 6. Go to the "M-Maintenance" menu and select "T-Control Testing."  
Select "D" for dry control test. After selection is made press ENTER.
- (✓) 7. Air blank completed.
- (✓) 8. Calibration check completed. Test results 0. 096 AC.
- (✓) 9. Air blank completed.
- (✓) 10. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Terry Gallegos

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002265  
Core Version: 8105.48  
04/02/2014 18:03:48

10 of 12  
TS4

-----  
Standard Lot#: 30412100A1  
Last Changed By: T. GALLEGOS #28847  
-----

sc

-----  
QAS: T. GALLEGOS #28847  
TUCSON PD  
-----

-----  
Operator: T. GALLEGOS #28847  
TUCSON PD CRIME LAB  
-----

Test	g/210L	Time
Air Blank	0.000	18:04:16
0.100 Cal Check	0.096	18:04:37
Air Blank	0.000	18:05:06



12 of 12  
TS4 SC

# Certificate of Calibration

This is to certify the calibration of **Intoxilyzer**® serial number 80-000065, manufactured by CMI Inc., a subsidiary of MPD, Inc. of Owensboro, Kentucky, was tested and found to conform to the National Highway Traffic Safety Administration (NHTSA) Standard for Devices to Measure Breath Alcohol (Federal Register, Vol.58 No.179, pp 48705-48710, Sept. 17, 1993) for accuracy and precision. Reference materials are traceable through the National Institute of Standards and Technology (NIST) to the International System of Units (SI).

Date March 12, 2014 Signed Nano Phillip

Technician



316 East 9<sup>th</sup> Street  
Owensboro, KY 42303  
USA

Part No. 650517 Rev.A

Pursuant to Section 28.1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.