



CERTIFICATION by CUSTODIAN of RECORDS  
 Tucson Police Department Crime Laboratory  
 1306 W. Miracle Mile, Tucson, AZ 85705

*MS*  
 11-25-15

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002265 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 11-5-2015 Cal 1 0.099 AC Cal 2 0.099 AC

Between the dates of 10-8-2015 and 11-5-2015 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.  
 The instrument was returned to service on -.

*Terry Gallegos*

Terry S. Gallegos


11-25-2015

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 25<sup>th</sup>, day of November, 2015.

*Elma Escarcega*

Notary Public

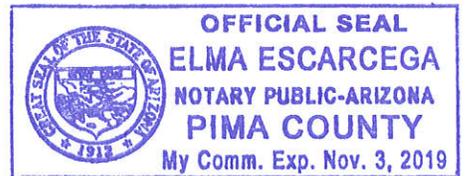

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1 of 7  
TS  
DB  
11-25-15

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department

DATE 11-5-2015 TIME 1345

INTOXILYZER SERIAL # 80-002265 LOCATION DAI-HQ

(  ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (  ) 1. Clock time check.
- (  ) 2. Date check.

OPERATIONAL TESTS

- (  ) 1. Alcohol-free subject test result 0. 000 AC.
- (  ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- (  ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (  ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 099 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE 

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002265  
Core Version: 8105.48  
11/05/2015

13:46:49

2077  
TS4  
JL

---

Standard Lot#: 11514100A2  
Last Changed By: T. GALLEGOS #28847

---

QAS: T. GALLEGOS #28847  
TUCSON PD

---

Operator: T. GALLEGOS #28847  
TUCSON PD

---

Subject: TEST, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	13:47:40
Diagnostic Test	Pass	13:48:15
Air Blank	0.000	13:48:42
0.100 Cal Check	0.099	13:49:03
Air Blank	0.000	13:49:33
Subject Test	0.000	13:50:02
Air Blank	0.000	13:50:29
Wait	ABT*	
Air Blank	0.000	13:51:04

\*Sequence Aborted

### Not a Successfully Completed Test Sequence

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002265  
Core Version: 8105.48  
11/05/2015

13:51:09

30f7  
TS4  
[Signature]

---

Standard Lot#: 11514100A2  
Last Changed By: T. GALLEGOS #28847

---

QAS: T. GALLEGOS #28847  
TUCSON PD

---

Operator: T. GALLEGOS #28847  
TUCSON PD

---

Subject: TEST, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	13:52:35
Diagnostic Test	Pass	13:53:10
Air Blank	0.000	13:53:38
0.100 Cal Check	RFI*	13:53:53
Air Blank	IPS**	13:54:00

\*RFI Detect  
\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002265  
Core Version: 8105.48  
11/05/2015

13:54:07

4 of 7  
TJh  
JL

---

Standard Lot#: 11514100A2  
Last Changed By: T. GALLEGOS #28847

---

QAS: T. GALLEGOS #28847  
TUCSON PD

---

Operator: T. GALLEGOS #28847  
TUCSON PD

---

Subject: TEST, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	13:54:55
Diagnostic Test	Pass	13:55:30
Air Blank	0.000	13:55:58
0.100 Cal Check	0.100	13:56:19
Air Blank	0.000	13:56:48
Subject Test	0.194*	13:59:56
Air Blank	0.000	14:00:30
Air Blank	0.000	14:00:57
Subject Test	XXX**	14:01:48
Air Blank	0.000	14:02:22

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

### Not a Successfully Completed Test Sequence

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002265  
Core Version: 8105.48  
11/05/2015 14:04:34

Soft  
TS4  
dd

-----  
Standard Lot#: 11514100A2  
Last Changed By: T. GALLEGOS #28847  
-----

QAS: T. GALLEGOS #28847  
TUCSON PD  
-----

Operator: T. GALLEGOS #28847  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	14:05:03
0.100 Cal Check	0.099	14:05:24
Air Blank	0.000	14:05:54

Preventative Maintenance  
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000  
STANDARD CALIBRATION CHECK PROCEDURE

60f7  
TS4  
AK  
11-25-15

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department  
DATE 11-5-2015 TIME 1345  
Intoxilyzer Serial # 80-002265 Location DUI-HQ

- (  ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C.
- (  ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (  ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (  ) 4. Air blank completed.
- (  ) 5. Calibration check completed. Test results 0. 099 AC.
- (  ) 6. Air blank completed.
- (  ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE 

DPS Form Exh G-2 (Rev. 05-01)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002265  
Core Version: 8105.48  
11/05/2015 13:45:02

7077  
TS4  
[Signature]

-----  
Standard Lot#: 11514100A2  
Last Changed By: T. GALLEGOS #28847  
-----

QAS: T. GALLEGOS #28847  
TUCSON PD  
-----

Operator: T. GALLEGOS #28847  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	13:45:46
0.100 Cal Check	0.099	13:46:07
Air Blank	0.000	13:46:37

Preventative Maintenance  
Performed -- Timer Reset



TUCSON POLICE DEPARTMENT  
CRIME LABORATORY

INTOXILYZER SERVICE RECORD  
FOR  
CALIBRATION GAS TANK CHANGE

10 of 2  
TS4  
11-25-15

Date: 11-5-2015  
QAS: Terry Gallegos 28847  
Intoxilyzer Serial Number: 80-002265  
New Lot Number: 15615100A1  
New Expiration Date: 07/05/2017

- (✓) 1. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 2. Remove the old gas tank and install the new one.
- (✓) 3. Hit the "ESC" key twice and swipe the QAS card.
- (✓) 4. Go to the "M-Maintenance" menu and select "S-Standard Change."  
Select "D" for dry gas and enter the following information:  
Enter Target Value: 0.100  
Enter Control Lot #: From the tank label  
Enter Exp Date: From the tank label (mm/dd/yyyy)  
Edit info Y/N? Y or N  
Update Std Info? Y ENTER  
"Please wait saving..."
- (✓) 5. Hit the "ESC" key twice.
- (✓) 6. Go to the "M-Maintenance" menu and select "T-Control Testing."  
Select "D" for dry control test. After selection is made press ENTER.
- (✓) 7. Air blank completed.
- (✓) 8. Calibration check completed. Test results 0. 098 AC.
- (✓) 9. Air blank completed.
- (✓) 10. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Terry Gallegos

2 of 2  
TS4  
*[Signature]*

INTOXILYZER 8000  
Location: TUCSON PD  
Serial Number: 80-002265  
Core Version: 8105.48  
11/05/2015 18:03:04

-----  
Standard Lot#: 15615100A1  
Last Changed By: T. GALLEGOS #28847  
-----

QAS: T. GALLEGOS #28847  
TUCSON PD  
-----

Operator: T. GALLEGOS #28847  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	18:03:40
0.100 Cal Check	0.098	18:04:01
Air Blank	0.000	18:04:30

Preventative Maintenance  
Performed -- Timer Reset