



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002265 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

WLL
3/2/2016

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 2-25-2016 Cal 1 0.099 AC Cal 2 0.098 AC

Between the dates of 1-28-2016 and 2-25-2016 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on
The instrument was returned to service on

[Handwritten signature of Terry S. Gallegos]

Terry S. Gallegos


3.1.2016

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 15th day of March, 2016.

[Handwritten signature of Elma Escarcega]

Notary Public

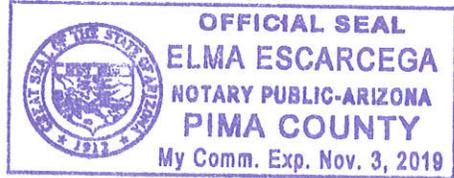

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

10 of 9
TS4

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

WLL
3/2/2016

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Terry Gallegos AGENCY Tucson Police Department

DATE 2-25-2016 TIME 1108

INTOXILYZER SERIAL # 80-002265 LOCATION DUI-HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 098 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE 

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002265
Core Version: 8105.48
02/25/2016

11:10:13

2 of 9
TS4
WLL

Standard Lot#: 15615100A1
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST,SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:11:06
Diagnostic Test	Pass	11:11:41
Air Blank	0.000	11:12:09
0.100 Cal Check	0.099	11:12:29
Air Blank	0.000	11:12:59
Subject Test	0.000	11:13:16
Air Blank	0.000	11:13:43
Wait	ABT*	
Air Blank	0.000	11:14:17

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002265
Core Version: 8105.48
02/25/2016

11:14:25

30f9
B4
WLL

Standard Lot#: 15615100A1
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST,SQAP

DOB: 11/11/1987

Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:16:12
Diagnostic Test	Pass	11:16:47
Air Blank	0.000	11:17:14
0.100 Cal Check	RFI*	11:17:31
Air Blank	IPS**	11:17:38

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002265
Core Version: 8105.48
02/25/2016

11:17:44

4 of 9
TS4
WLL

Standard Lot#: 15615100A1
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST,SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:18:35
Diagnostic Test	Pass	11:19:10
Air Blank	0.000	11:19:37
0.100 Cal Check	0.099	11:19:58
Air Blank	0.000	11:20:28
Subject Test	0.068*	11:23:35
Air Blank	0.000	11:24:09
Air Blank	0.000	11:24:36
Subject Test	XXX**	11:24:59
Air Blank	0.000	11:25:33

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002265
Core Version: 8105.48
02/25/2016 11:26:01

50f9
T54
WLL

Standard Lot#: 15615100A1
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	11:26:30
0.100 Cal Check	0.098	11:26:51
Air Blank	0.000	11:27:21

Preventative Maintenance
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

6079
TDH

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

WLL
3/2/2016

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department
DATE 2-25-2016 TIME 1108
Intoxilyzer Serial # 80-002265 Location DUI-HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C + 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 099 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE 

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002265
Core Version: 8105.48
02/25/2016 11:08:42

70F9
704

Standard Lot#: 15615100A1
Last Changed By: T. GALLEGOS #28847

WLL

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	11:09:10
0.100 Cal Check	0.099	11:09:31
Air Blank	0.000	11:10:01

Preventative Maintenance
Performed -- Timer Reset



TUCSON POLICE DEPARTMENT
CRIME LABORATORY

80f9
TS4

INTOXILYZER SERVICE RECORD
FOR
CALIBRATION GAS TANK CHANGE

WLL
3/2/2016

Date: 2-25-2016
QAS: Terry Gallegos 28847
Intoxilyzer Serial Number: 80-002265
New Lot Number: 12415100A4
New Expiration Date: 06-05-2017

- (✓) 1. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 2. Remove the old gas tank and install the new one.
- (✓) 3. Hit the "ESC" key twice and swipe the QAS card.
- (✓) 4. Go to the "M-Maintenance" menu and select "S-Standard Change."
Select "D" for dry gas and enter the following information:
Enter Target Value: 0.100
Enter Control Lot #: From the tank label
Enter Exp Date: From the tank label (mm/dd/yyyy)
Edit info Y/N? Y or N
Update Std Info? Y ENTER
"Please wait saving..."
- (✓) 5. Hit the "ESC" key twice.
- (✓) 6. Go to the "M-Maintenance" menu and select "T-Control Testing."
Select "D" for dry control test. After selection is made press ENTER.
- (✓) 7. Air blank completed.
- (✓) 8. Calibration check completed. Test results 0. 098 AC.
- (✓) 9. Air blank completed.
- (✓) 10. Remove printed record. Attach the record to the completed checklist.

SIGNATURE 

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002265
Core Version: 8105.48
02/25/2016 11:42:49

9 of 9
754

WLL

Standard Lot#: 12415100A4
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	11:43:18
0.100 Cal Check	0.098	11:43:39
Air Blank	0.000	11:44:08

Preventative Maintenance
Performed -- Timer Reset