



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

T54
12-13-07

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002345 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 11/19/07 Cal 1 0.100 AC Cal 2 0.100 AC

Between the dates of 10/19/07 and 11/19/07 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on —.
The instrument was returned to service on —.


Sarah L. Cione

Sarah L. Cione

11/29/07

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 29, day of November, 2007.

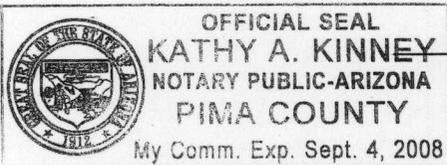

Kathy A Kinney

Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

1 of 7sc

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

BH
12-13-07

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 11/19/07 TIME 9:44am

INTOXILYZER SERIAL # 80-002345 LOCATION DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 100 AC

Instrument is operating properly and accurately. YES NO

COMMENTS none

SIGNATURE Sarah Cione

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002345
11/19/2007

09:48:01

2 of 7 sc
P4

Standard Lot#: 706602I
Last Changed By: J. MORGAN III #42130

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:49:40
Diagnostic Test	Pass	09:50:07
Air Blank	0.000	09:50:35
0.100 Cal Check	0.100	09:50:54
Air Blank	0.000	09:51:23
Subject Test	0.000	09:51:44
Air Blank	0.000	09:52:11
Five Minute Wait	ABT*	
Air Blank	0.000	09:53:03

*Sequence Aborted

Not a Successfully Completed Test Sequence

Pursuant to section 20-1527, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002345
11/19/2007 09:53:31

3 of 7 sc
174

Standard Lot#: 706602I
Last Changed By: J. MORGAN III #42130

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:54:33
Diagnostic Test	Pass	09:55:00
Air Blank	RFI*	09:55:12
Air Blank	IPS**	09:55:31

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

Pursuant to section 15-27, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002345
11/19/2007 09:56:01

4 of 7sc

AB

Standard Lot#: 706602I
Last Changed By: J. MORGAN III #42130

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP,TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:56:51
Diagnostic Test	Pass	09:57:19
Air Blank	0.000	09:57:46
0.100 Cal Check	0.100	09:58:05
Air Blank	0.000	09:58:34
Subject Test	0.168*	10:01:43
Air Blank	0.000	10:02:15
Air Blank	0.000	10:02:42
Subject Test	XXX**	10:03:12
Air Blank	0.000	10:03:45

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

Pursuant to section 1-217, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-002345
11/19/2007 10:04:42

5 of 7 sc

AK

Standard Lot#: 706602I
Last Changed By: J. MORGAN III #42130

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	10:05:10
Cal Check	0.100	10:05:29
Air Blank	0.000	10:05:58

Sarah Cione

Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6 of 7 SC
Ply

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 11/19/2007 TIME 9:44 am
Intoxilyzer Serial # 80-002345 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-002345
11/19/2007 09:44:24

7 of 7 sc
DH

Standard Lot#: 706602I
Last Changed By: J. MORGAN III #42130

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:44:51
Cal Check	0.100	09:45:10
Air Blank	0.000	09:45:39

Sarah Cione

Operator's Signature