



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

CU
2-20-08

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002345 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 2/14/08 Cal 1 0.101 AC Cal 2 0.101 AC

Between the dates of 1/16/08 and 2/14/08 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.
The instrument was returned to service on -.


Sarah Cione

Sarah L. Cione

2/15/08

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 15, day of February, 2008.

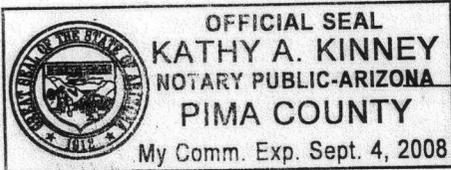

Kathy A Kinney
Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

1 of 8 sc
Ch
2-22-08

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 2/14/08 TIME 9:48

INTOXILYZER SERIAL # 80-002345 LOCATION DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 101 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Sarah Cione

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002345
02/14/2008 09:51:36

2 of 8 sc
CU
2-20-08

Standard Lot#: 725402I
Last Changed By: J. MORGAN III #42130

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:52:29
Diagnostic Test	Pass	09:52:57
Air Blank	0.000	09:53:24
0.100 Cal Check	0.101	09:53:43
Air Blank	0.000	09:54:12
Subject Test	0.000	09:54:31
Air Blank	0.000	09:54:59
Five Minute Wait	ABT*	
Air Blank	0.000	09:55:44

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002345
02/14/2008 09:57:17

3 of 85c
CU
2.22-08

Standard Lot#: 725402I
Last Changed By: J. MORGAN III #42130

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	09:58:14
Air Blank	IPS**	09:58:22

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location: TUCSON POLICE DEPT

Serial Number : 80-002345

02/14/2008

10:00:45

4 of 8 sc

lh

2-20-08

Standard Lot#: 725402I

Last Changed By: J. MORGAN III #42130

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, TEST

DOB: 11/11/1111

Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:01:36
Diagnostic Test	Pass	10:02:03
Air Blank	0.000	10:02:31
0.100 Cal Check	0.101	10:02:50
Air Blank	0.000	10:03:19
Subject Test	RGE*	10:03:53
Air Blank	0.000	10:04:29

*Range Exceeded

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number : 80-002345
02/14/2008 10:04:46

5 of 8 sc

lu
2-20-08

Standard Lot#: 725402I
Last Changed By: J. MORGAN III #42130

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:06:01
Diagnostic Test	Pass	10:06:28
Air Blank	0.000	10:06:56
0.100 Cal Check	0.101	10:07:15
Air Blank	0.000	10:07:44
Subject Test	0.028*	10:10:53
Air Blank	0.000	10:11:23
Air Blank	0.000	10:11:51
Subject Test	XXX**	10:12:11
Air Blank	0.000	10:12:44

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-002345
02/14/2008 10:14:18

6 of 85c

Ch

2-20-08

Standard Lot#: 725402I
Last Changed By: J. MORGAN III #42130

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	10:14:45
Cal Check	0.101	10:15:04
Air Blank	0.000	10:15:34

Sarah Cione
Operator's Signature

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000
STANDARD CALIBRATION CHECK PROCEDURE

7 of 8 sc
CM
2-20-08

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 2/14/08 TIME 9:48
Intoxilyzer Serial # 80-002345 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 101 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-002345
02/14/2008 09:48:58

8 of 8 SC

CU

2-20-08

Standard Lot#: 725402I
Last Changed By: J. MORGAN III #42130

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:49:25
Cal Check	0.101	09:49:44
Air Blank	0.000	09:50:13

Sarah Cione

Operator's Signature



TUCSON POLICE DEPARTMENT
CRIME LABORATORY

1 of 25c

cu

2-22-08

INTOXILYZER SERVICE RECORD
FOR
CALIBRATION GAS TANK CHANGE

Date: FEBRUARY 9, 2008

QAS: OFC. MORGAN #42130

Intoxilyzer Serial Number: 80-002345

New Lot Number: 727501I

New Expiration Date: OCTOBER 4, 2009

- (✓) 1. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 2. Remove the old gas tank and install the new one.
- (✓) 3. Hit the "ESC" key twice and swipe the QAS card.
- (✓) 4. Go to the "M-Maintenance" menu and select "S-Standard Change."
Select "D" for dry gas and enter the following information:
Enter Target Value: 0.100
Enter Control Lot #: From the tank label
Enter Exp Date: From the tank label (mm/dd/yyyy)
Edit info Y/N? Y or N
Update Std Info? Y ENTER
"Please wait saving..."
- (✓) 5. Hit the "ESC" key twice.
- (✓) 6. Go to the "M-Maintenance" menu and select "T-Control Testing."
Select "D" for dry control test. After selection is made press
ENTER.
- (✓) 7. Air blank completed.
- (✓) 8. Calibration check completed. Test results 0. 092 AC.
- (✓) 9. Air blank completed.
- (✓) 10. Remove and sign printed record. Attach the record to the completed checklist.

SIGNATURE *[Signature]* #42130

2 of 2 SC

Oh
2-20-08

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-002345
02/09/2008 13:47:52

Standard Lot#: 7275011
Last Changed By: J. MORGAN III #42130

QAS: J. MORGAN III #42130
TUCSON PD

Test	gr/210L	Time
Air Blank	0.000	13:48:19
Cal Check	0.092	13:48:38
Air Blank	0.000	13:49:08

J. Morgan III #42130

Operator's Signature



TUCSON POLICE DEPARTMENT
CRIME LABORATORY

1 of 2 sc
CU
2-20-08

INTOXILYZER SERVICE RECORD
FOR
CALIBRATION GAS TANK CHANGE

Date: FEBRUARY 10, 2008

QAS: OFC. MORGAN #42130

Intoxilyzer Serial Number: 80-002345

New Lot Number: 725402-I

New Expiration Date: SEPT. 12, 2009

- (✓) 1. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 2. Remove the old gas tank and install the new one.
- (✓) 3. Hit the "ESC" key twice and swipe the QAS card.
- (✓) 4. Go to the "M-Maintenance" menu and select "S-Standard Change."
Select "D" for dry gas and enter the following information:
 Enter Target Value: 0.100
 Enter Control Lot #: From the tank label
 Enter Exp Date: From the tank label (mm/dd/yyyy)
 Edit info Y/N? Y or N
 Update Std Info? Y ENTER
 "Please wait saving..."
- (✓) 5. Hit the "ESC" key twice.
- (✓) 6. Go to the "M-Maintenance" menu and select "T-Control Testing."
Select "D" for dry control test. After selection is made press
ENTER.
- (✓) 7. Air blank completed.
- (✓) 8. Calibration check completed. Test results 0. 098 AC.
- (✓) 9. Air blank completed.
- (✓) 10. Remove and sign printed record. Attach the record to the completed checklist.

SIGNATURE *J. Morgan* #42130


2 of 250
CM
2-20-08

INTOXILYZER 8000
Location: TUCSON POLICE DEPT
Serial Number: 80-002345
02/10/2006 21:52:37

Standard Lot#: 7254021
Last Changed By: J. MORGAN III #42130

OAS: J. MORGAN III #42130
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	21:53:04
Cal Check	1.098	21:53:23
Air Blank	0.000	21:53:52

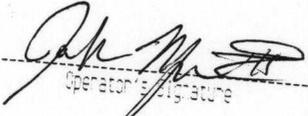
Operator's Signature