



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

124
6-3-08

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002345 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 5-20-08 Cal 1 0.104 AC Cal 2 0.103 AC

Between the dates of N/A and N/A the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.
The instrument was returned to service on 5-20-08.


Clare Conley

Clare L. Conley

5-20-08

ADPS QAS Permit # 33314

SUBSCRIBED AND SWORN to before me this 20, day of May, 2008.

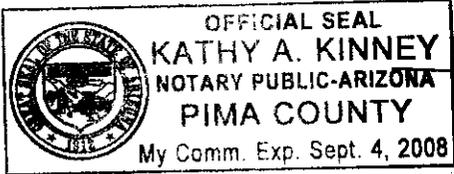

Kathy A Kinney
Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.



TUCSON POLICE DEPARTMENT
CRIME LABORATORY

1-6-12
cm

INTOXILYZER SERVICE RECORD

64
6-3-08

QAS: Clare Conley

Intoxilyzer Serial Number: 80-002345

Date Returned to Service: 5-20-08

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

2 of 12 Cu
TS4
6-3-08

QA SPECIALIST Clare Conley AGENCY Tucson Police Department

DATE 5-20-08 TIME 1154

INTOXILYZER SERIAL # 80-002345 LOCATION Crime Lab

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 103 AC

Instrument is operating properly and accurately. YES NO

COMMENTS None

SIGNATURE Clare Conley

3 of 12 Cu

PA

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002345
05/20/2008

11:54:56

Standard Lot#: 735401G
Last Changed By: C. CONLEY #34517

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:55:52
Diagnostic Test	Pass	11:56:19
Air Blank	0.000	11:56:47
0.100 Cal Check	0.102	11:57:05
Air Blank	0.000	11:57:35
Subject Test	0.000	11:57:51
Air Blank	ABT*	11:57:57
Air Blank	0.000	11:58:25

*Sequence Aborted

Not a Successfully Completed Test Sequence

Pursuant to section 28-1527, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

4 of 12 Cn

BA

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002345
05/20/2008

11:59:47

Standard Lot#: 735401G
Last Changed By: C. CONLEY #34517

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	12:00:45
Diagnostic Test	Pass	12:01:12
Air Blank	0.000	12:01:40
0.100 Cal Check	0.102	12:01:58
Air Blank	0.000	12:02:28
Subject Test	0.238*	12:05:36
Air Blank	0.000	12:06:09
Air Blank	0.000	12:06:36
Subject Test	XXX**	12:06:56
Air Blank	0.000	12:07:30

*Deficient Sample
 **Invalid Sample
 Begin new deprivation period

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

5 of 12 cu
Bl

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002345
05/20/2008

12:07:58

Standard Lot#: 735401G
Last Changed By: C. CONLEY #34517

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	12:08:35
Air Blank	IPS**	12:08:43

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

Pursuant to section 28-1527, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002345
05/20/2008


12:09:46

6 of 12
Cu
Ph

Standard Lot#: 735401G
Last Changed By: C. CONLEY #34517

QAS: C. CONLEY #34517
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	12:10:13
Cal Check	0.103	12:10:32
Air Blank	0.000	12:11:02

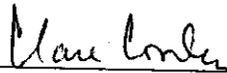

Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

7 of 12 Cu
754
6.3.08

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Clare Conley AGENCY Tucson Police Department
DATE 5-20-08 TIME 1154
Intoxilyzer Serial # 80-002345 Location Crime Lab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 150 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 104 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Clare Conley

DPS Form Exh G-2 (Rev. 05-01)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

8 of 12 cu
B4

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002345
05/20/2008 11:53:32

Standard Lot#: 735401G
Last Changed By: C. CONLEY #34517

QAS: C. CONLEY #34517
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	11:54:00
Cal Check	0.104	11:54:19
Air Blank	0.000	11:54:48

Care Conley
Operator's Signature



TUCSON POLICE DEPARTMENT
CRIME LABORATORY

9/12
CW
154
6-3-08

INTOXILYZER SERVICE RECORD
FOR
CALIBRATION GAS TANK CHANGE

Date: 5-20-08
 QAS: Clare Conley
 Intoxilyzer Serial Number: 80-002345
 New Lot Number: 7354016
 New Expiration Date: 12-28-09

- (✓) 1. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 2. Remove the old gas tank and install the new one.
- (✓) 3. Hit the "ESC" key twice and swipe the QAS card.
- (✓) 4. Go to the "M-Maintenance" menu and select "S-Standard Change." Select "D" for dry gas and enter the following information:
 Enter Target Value: 0.100
 Enter Control Lot #: From the tank label
 Enter Exp Date: From the tank label (mm/dd/yyyy)
 Edit info Y/N? Y or N
 Update Std Info? Y ENTER
 "Please wait saving..."
- (✓) 5. Hit the "ESC" key twice.
- (✓) 6. Go to the "M-Maintenance" menu and select "T-Control Testing." Select "D" for dry control test. After selection is made press ENTER.
- (✓) 7. Air blank completed.
- (✓) 8. Calibration check completed. Test results 0.102 AC.
- (✓) 9. Air blank completed.
- (✓) 10. Remove and sign printed record. Attach the record to the completed checklist.

SIGNATURE Clare Conley

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002345
05/20/2008

11:51:36

105612
lu
Tsh

Standard Lot#: 735401G
Last Changed By: C. CONLEY #34517

QAS: C. CONLEY #34517
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	11:52:04
Cal Check	0.102	11:52:23
Air Blank	0.000	11:52:52

Clare Conley

Operator's Signature

INTOX
PBT
TOXITEST

DATE: 4-17-08
RECEIVED VIA: LWS

CMI INC.
316 E. 9th St., Owensboro, KY 42303
270-685-6200 — Fax 270-685-6288
1-800-835-0690

SERVICE WORK ORDER
316457 CM
11/08/12

B
I
L
L
T
O

ATTN:

CITY STATE ZIP

CUSTOMER NO. 105075

PHONE

S
H
I
P
T
O

Tucson Police Dept.

ATTN: Terry Gallegos

270 S. Stone Ave.

CITY STATE ZIP
Tucson AZ 85701

PHONE 520-791-4499 ext 1443

MODEL # 8000 S/N 80-002343 DATE RCVD. P.O. # 20669

PROD. CODE EXTRA PARTS RCVD. REQ. BILL CODE SERVICE LOCATION CODE

DESCRIPTION OF PROBLEM: gas shell, BT, new battery

DSD Failure

QUOTE WARRANTY SERVICE DEPT. NEW

- PRELIMINARY TEST DISMANTLE, INSPECT, DIAGNOSTIC TEST
- POWER SUPPLY TEST, INSPECT REPLACE
- MAIN PWB's - REMOVE & CLEAN, REINSTALL, TEST
- PWB's DIAGNOSE, REPAIR REPLACE
- PWB's - SURFACE CLEAN IN PLACE
- FILTER WHEEL - REMOVE, CLEAN, REINSTALL, CHECK CELL INTERIOR
- CHOPPER WHEEL MOTOR: TEST & INSPECT REPLACE
- CELL - REMOVE, CLEAN, REINSTALL, REPLACE LENSES/GASKETS REQ'D
- IR LIGHT SOURCE REMOVE, REALIGN REPLACE
- PRINTER ASSY. DISASSEMBLE CLEAN & LUBRICATE
- REMOVE & REPLACE BELT REALIGNMENT
- REPLACE TEST PRINTER FUNCTIONS
- BATTERIES - REPLACE CPU RAM
- REPLACE PARTS REASSEMBLE CALIBRATION FINAL TEST
- MISC.
- PBT PRELIMINARY TESTING - PRESSURE, BATTERY, ZERO SET
- DISASSEMBLE REPLACE PARTS REASSEMBLE
- TEST CALIBRATE FINAL TEST
- TOXITEST INCOMING LEAK TEST VISUAL/MECH. INSPECTION
- TEST THERMOMETER & THERMOSTAT DISASSEMBLE
- REPLACE PARTS REASSEMBLE OUTGOING & LEAK TEST
- FINAL TEST
- MISC.

PARTS USED			
QTY	ID NO.	DESCRIPTION	CODE
1	122982	Source Assy	
1	441657	End block Gasket	
1	340208	3" check Valve	
1	330506	Keyboard cable	
1	441209	Breath hose insert	
1	441277	3" check Valve	
1	441144	RT Angle FIT	
1	151027	Keyboard, P.I. block	

COMMENTS: This unit needed a new source. There was an intermittent connection in the Keyboard, so I replaced the internal cable & keyboard. In addition, the cell needed two new check valves & a new end block Gasket. Finally, the unit needed a new breath hose insert to create a better seal.

HOURS: 5 hours

DATE: 5/12/08

DATE: 5-6-08

Pursuant to Section 28.1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department

SHIPPING DATE/CHAI

SPECIAL SHIPPING INSTRUCTIONS

12 of 12
04
B4

CERTIFICATE OF CALIBRATION

INTOXILYZER® 8000

This is to certify that the calibration of INTOXILYZER serial number 80-002345, manufactured by CMI, subsidiary of MPD, Inc. of Owensboro, Kentucky, was tested and found to be in compliance with the National Safety Administration Standard for Devices to Measure Breath Alcohol (F.R., No. 179 48705-48710 Sept. 17, 1993) and with International Recommendation OIML R126 Edition 1998 (E). Calibration solutions are traceable to NIST (NBS) standard material 1828.

Date 5-6-08

Signed



Technician

CMI INC.

316 EAST 9TH STREET
OWENSBORO, KY. 42303

Pursuant to Section 28.1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department