



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002345 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

TS4
8.1.08

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 7-25-08 Cal 1 0.101 AC Cal 2 0.100 AC

Between the dates of 6-27-08 and 7-25-08 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A
The instrument was returned to service on N/A


Clare Conley

Clare L. Conley

7-25-08

ADPS QAS Permit # 33314

SUBSCRIBED AND SWORN to before me this 25th, day of July, 20 08.

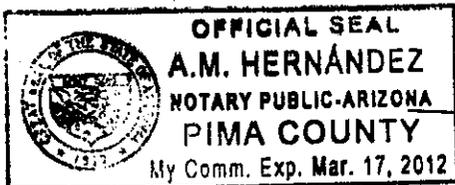

A.M. Hernandez

Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000
STANDARD QUALITY ASSURANCE PROCEDURE

1567
Cu
7/24/08

QA SPECIALIST Clare Conley AGENCY Tucson Police Department
DATE 7-25-08 TIME 0615
INTOXILYZER SERIAL # 80-002345 LOCATION Crime Lab

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 100 AC

Instrument is operating properly and accurately. YES NO

COMMENTS None.

SIGNATURE Clare Conley

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002345
07/25/2008

06:16:57

2 of 7
C
104

Standard Lot#: 735401G
Last Changed By: C. CONLEY #34517

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	06:17:48
Diagnostic Test	Pass	06:18:16
Air Blank	0.000	06:18:43
0.100 Cal Check	0.101	06:19:02
Air Blank	0.000	06:19:31
Subject Test	0.000	06:19:50
Air Blank	ABT*	06:19:56
Air Blank	0.000	06:20:24

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002345
07/25/2008

06:20:30

Standard Lot#: 735401G
Last Changed By: C. CONLEY #34517

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	06:21:05
Air Blank	IPS**	06:21:24

*RFI Detect
**Improper Sample

347
cu
154

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002345
07/25/2008

06:21:37

Standard Lot#: 735401G
Last Changed By: C. CONLEY #34517

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	06:22:29
Diagnostic Test	Pass	06:22:57
Air Blank	0.000	06:23:24
0.100 Cal Check	0.101	06:23:43
Air Blank	0.000	06:24:12
Subject Test	DEF*	06:27:21
Air Blank	0.000	06:27:48
Air Blank	0.000	06:28:16
Subject Test	XXX**	06:28:43
Air Blank	0.000	06:29:17

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

4867
C
TS4

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002345
07/25/2008

06:29:55

Standard Lot#: 735401G
Last Changed By: C. CONLEY #34517

QAS: C. CONLEY #34517
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	06:30:23
Cal Check	0.100	06:30:41
Air Blank	0.000	06:31:11

Clare Conley
Operator's Signature

5067
en
TS4

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6 of 7
An
134
8.1.08

QA Specialist Clare Conley AGENCY Tucson Police Department
DATE 7-25-08 TIME 0615
Intoxilyzer Serial # 80-202345 Location Crime Lab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 101 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Clare Conley

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002345
07/25/2008

06:15:34

Standard Lot#: 735401G
Last Changed By: C. CONLEY #34517

QAS: C. CONLEY #34517
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	06:16:01
Cal Check	0.101	06:16:20
Air Blank	0.000	06:16:49

Clare Conley

Operator's Signature

7567
CW
TS4