



CERTIFICATION by CUSTODIAN of RECORDS
 Tucson Police Department Crime Laboratory
 270 S. Stone Avenue Tucson, AZ 85701-1917

Ch
 11-21-08

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002345 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 11/14/08 Cal 1 0.103 AC Cal 2 0.103 AC

Between the dates of 10/17/08 and 11/14/08 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.
 The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione


11/20/08

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 20th day of November, 2008.

A. Hernandez

Notary Public

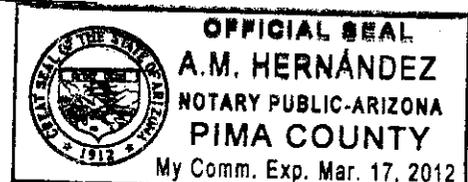

*Includes documents specified in A.A.C. R13-10-104.A.5

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

1 of 7sc

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000
STANDARD QUALITY ASSURANCE PROCEDURE

CU
11-21-08

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 11/14/08 TIME 9:29

INTOXILYZER SERIAL # 80-002345 LOCATION Team 3

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 103 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

2 of 7sc

CM

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002345
11/14/2008

09:30:41

Standard Lot#: 735401G
Last Changed By: C. CONLEY #34517

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP,TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:31:27
Diagnostic Test	Pass	09:31:54
Air Blank	0.000	09:32:21
0.100 Cal Check	0.103	09:32:40
Air Blank	0.000	09:33:10
Subject Test	0.000	09:33:27
Air Blank	0.000	09:33:55
Five Minute Wait	ABT*	
Air Blank	0.000	09:34:34

*Sequence Aborted

Not a Successfully Completed Test Sequence

3 of 7sc

DM

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002345
11/14/2008

09:34:40

Standard Lot#: 735401G
Last Changed By: C. CONLEY #34517

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	09:35:10
Air Blank	IPS**	09:35:18

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

4 of 7sc

cu

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002345
11/14/2008 09:35:23

Standard Lot#: 735401G
Last Changed By: C. CONLEY #34517

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SOAP, TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:36:07
Diagnostic Test	Pass	09:36:34
Air Blank	0.000	09:37:01
0.100 Cal Check	0.102	09:37:20
Air Blank	0.000	09:37:49
Subject Test	0.168*	09:40:58
Air Blank	0.000	09:41:30
Air Blank	0.000	09:41:57
Subject Test	XXX**	09:42:28
Air Blank	0.000	09:43:02

*Deficient Sample
**Invalid Sample
Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002345
11/14/2008

09:43:27

5 of 7 sc

cu

Standard Lot#: 735401G
Last Changed By: C. CONLEY #34517

QAS: S. CIONE #50616
TUCSON PD

Test g/210L Time

Air Blank 0.000 09:43:54
Cal Check 0.103 09:44:13
Air Blank 0.000 09:44:42

Sarah Cione
Operator's Signature

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

6 of 7 sc

Ch
11-21-08

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 11/14/08 TIME 9:29
Intoxilyzer Serial # 80-002345 Location Team 3

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 103 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002345
11/14/2008 09:29:17

7 of 7 sc

cu

Standard Lot#: 735401G
Last Changed By: C. CONLEY #34517

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:29:45
Cal Check	0.103	09:30:04
Air Blank	0.000	09:30:33

Sarah Cione

Operator's Signature