



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002345 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 03/11/10 Cal 1 0.102 AC Cal 2 0.101 AC

Between the dates of 02/11/10 and 03/11/10 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.
The instrument was returned to service on N/A.

Cheryl K. Straub

Cheryl K. Straub

03/11/10


ADPS QAS Permit # 33408

SUBSCRIBED AND SWORN to before me this 11th, day of March, 20 10.

A.M. Hernandez

Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

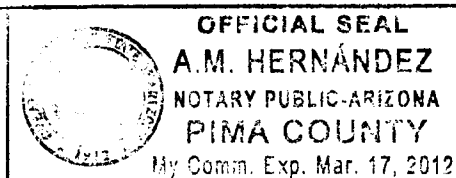

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

TS4
3.11.2010QA SPECIALIST Cheryl STRAUB AGENCY Tucson Police DepartmentDATE 03/11/10 TIME 07:05INTOXILYZER SERIAL # 80-002345 LOCATION TEAM 3

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
-
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0.
- 000
- AC.
-
- (✓) 2. Error recognition logic system functioning.
-
- Not a Successfully Completed Test Sequence printed
-
- (✓) 3. Proper sample recognition system.
-
- Not a Successfully Completed Test Sequence printed
-
- Deficient sample printed.
-
- (✓) 4. Standard Calibration Check standard 0.
- 100
- AC
-
- Result: 0.
- 101
- AC

Instrument is operating properly and accurately. YES ☒ NO ☐

COMMENTS

SIGNATURE

Cheryl Straub

287
CS

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002345
03/11/2010

07:07:42

PH

Standard Lot#: 642845
Last Changed By: C. STRAUB #50772

QAS: C. STRAUB #50772
TUCSON PD

Operator: C. STRAUB #50772
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	07:08:35
Diagnostic Test	Pass	07:09:02
Air Blank	0.000	07:09:30
0.100 Cal Check	0.101	07:09:49
Air Blank	0.000	07:10:18
Subject Test	0.000	07:10:38
Air Blank	ABT*	07:10:45
Air Blank	0.000	07:11:14

*Sequence Aborted

Not a Successfully Completed Test Sequence

377
CS

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002345
03/11/2010

07:11:20

54

Standard Lot#: 642845
Last Changed By: C. STRAUB #50772

QAS: C. STRAUB #50772
TUCSON PD

Operator: C. STRAUB #50772
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	07:11:54
Air Blank	IPS**	07:12:02

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

477
CS

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002345
03/11/2010

07:12:10

64

Standard Lot#: 642845
Last Changed By: C. STRAUB #50772

QAS: C. STRAUB #50772
TUCSON PD

Operator: C. STRAUB #50772
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	07:13:00
Diagnostic Test	Pass	07:13:28
Air Blank	0.000	07:13:55
0.100 Cal Check	0.102	07:14:14
Air Blank	0.000	07:14:44
Subject Test	0.096*	07:17:52
Air Blank	0.000	07:18:25
Air Blank	0.000	07:18:52
Subject Test	XXX**	07:19:23
Air Blank	0.000	07:19:57

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

597
CS

INTOXILYZER 8000

Location: ARIZONA

Serial Number: 80-002345

03/11/2010

07:20:14

TS4


Standard Lot#: 642845

Last Changed By: C. STRAUB #50772

QAS: C. STRAUB #50772

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	07:20:42
Cal Check	0.101	07:21:00
Air Blank	0.000	07:21:30

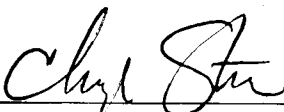

Operator's Signature

EXHIBIT G-2

677
CS

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

154
3.11.2010

QA Specialist Cheryl STRAUB AGENCY Tucson Police Department
DATE 03/11/10 TIME 07:05
Intoxilyzer Serial # 80-002345 Location TEAM 3

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution N/A AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 102 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Cheryl Straub

DPS Form Exh G-2 (Rev. 05-01)

727
CS

INTOXILYZER 8000

Location: ARIZONA

Serial Number: 80-002345

03/11/2010

07:05:34

154

Standard Lot#: 642845

Last Changed By: C. STRAUB #50772

QAS: C. STRAUB #50772

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	07:06:02
Cal Check	0.102	07:06:21
Air Blank	0.000	07:06:50



Operator's Signature