



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
270 S. Stone Avenue Tucson, AZ 85701-1917

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002345 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

TS4  
5.6.2010

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 05/06/10 Cal 1 0.103 AC Cal 2 0.101 AC

Between the dates of 04/08/10 and 05/06/10 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.  
The instrument was returned to service on N/A.

Cheryl K. Straub

Cheryl K. Straub


05/06/10

ADPS QAS Permit # 33408

SUBSCRIBED AND SWORN to before me this 6<sup>th</sup>, day of May, 2010.

A.M. Hernandez

Notary Public

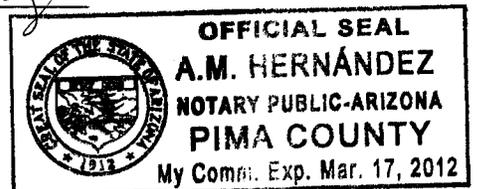

\*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

178  
CS

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000  
STANDARD QUALITY ASSURANCE PROCEDURE

TS4  
5.6.2010

QA SPECIALIST Cheryl STRAUB AGENCY Tucson Police Department  
DATE 05/06/10 TIME 06:50  
INTOXILYZER SERIAL # 80-002345 LOCATION TEAM 3

( ✓ ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- ( ✓ ) 1. Clock time check.
- ( ✓ ) 2. Date check.

OPERATIONAL TESTS

- ( ✓ ) 1. Alcohol-free subject test result 0. 000 AC.
- ( ✓ ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- ( ✓ ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- ( ✓ ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 101 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS When performing mouth alcohol TEST (INVALID  
SAMPLE TEST) RANGE WAS EXCEEDED; TEST REPEATED.  
CS.

SIGNATURE Cheryl Straub

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

27805

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002345  
05/06/2010

06:51:39

SM

---

Standard Lot#: 642845  
Last Changed By: C. STRAUB #50772

---

QAS: C. STRAUB #50772  
TUCSON PD

---

Operator: C. STRAUB #50772  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	06:52:32
Diagnostic Test	Pass	06:52:59
Air Blank	0.000	06:53:26
0.100 Cal Check	0.102	06:53:45
Air Blank	0.000	06:54:15
Subject Test	0.000	06:54:54
Air Blank	ABT*	06:55:01
Air Blank	0.000	06:55:29

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

398  
ES

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002345  
05/06/2010

06:55:36

1/24

Standard Lot#: 642845  
Last Changed By: C. STRAUB #50772

QAS: C. STRAUB #50772  
TUCSON PD

Operator: C. STRAUB #50772  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	06:56:07
Air Blank	IPS**	06:56:15

\*RFI Detect  
\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

488  
CS

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002345  
05/06/2010

06:56:20

154

Standard Lot#: 642845  
Last Changed By: C. STRAUB #50772

QAS: C. STRAUB #50772  
TUCSON PD

Operator: C. STRAUB #50772  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	06:57:09
Diagnostic Test	Pass	06:57:37
Air Blank	0.000	06:58:04
0.100 Cal Check	0.102	06:58:23
Air Blank	0.000	06:58:52
Subject Test	0.255*	07:02:01
Air Blank	0.000	07:02:33
Air Blank	0.000	07:03:01
Subject Test	RGE**	07:03:19
Air Blank	0.000	07:03:56

\*Deficient Sample  
\*\*Range Exceeded

**Not a Successfully Completed Test Sequence**

578  
CS

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002345  
05/06/2010

07:04:05

BY

---

Standard Lot#: 642845  
Last Changed By: C. STRAUB #50772

---

QAS: C. STRAUB #50772  
TUCSON PD

---

Operator: C. STRAUB #50772  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	07:04:58
Diagnostic Test	Pass	07:05:25
Air Blank	0.000	07:05:53
0.100 Cal Check	0.101	07:06:11
Air Blank	0.000	07:06:41
Subject Test	XXX*	07:06:59
Air Blank	0.000	07:07:33

\*Invalid Sample  
Begin new deprivation period

**Not a Successfully Completed Test Sequence**

678  
CS

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-002345  
05/06/2010 07:08:08

154

-----  
Standard Lot#: 642845  
Last Changed By: C. STRAUB #50772  
-----

QAS: C. STRAUB #50772  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	07:08:35
Cal Check	0.101	07:08:54
Air Blank	0.000	07:09:24



Operator's Signature

778 CS

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

TJH  
5.6.2010

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Cheryl STRAUB AGENCY Tucson Police Department  
DATE 05/06/10 TIME 06:50  
Intoxilyzer Serial # 80-002345 Location TEAM 3

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution N/A AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 103 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Cheryl Straub

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

888  
CS

TS4

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-002345  
05/06/2010

06:50:13

-----  
Standard Lot#: 642845  
Last Changed By: C. STRAUB #50772  
-----

QAS: C. STRAUB #50772  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	06:50:41
Cal Check	0.103	06:50:59
Air Blank	0.000	06:51:29



Operator's Signature