



CERTIFICATION by CUSTODIAN of RECORDS
 Tucson Police Department Crime Laboratory
 270 S. Stone Avenue Tucson, AZ 85701-1917

TS4
 6-7-2010

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002345 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 6/3/10 Cal 1 0.102 AC Cal 2 0.101 AC

Between the dates of 5/6/10 and 6/3/10 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.
 The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione


6/3/10

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 3rd, day of June, 2010.

A.M. Hernandez

Notary Public

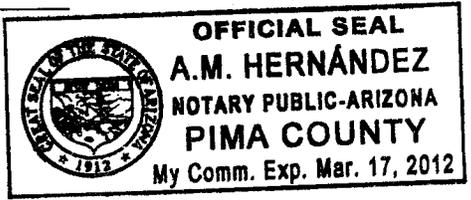

*Includes documents specified in A.A.C. R13-10-104.A.5

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

1 of 7 sc
By
6-7-2010

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000
STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department
DATE 6/3/10 TIME 900
INTOXILYZER SERIAL # 80-002345 LOCATION Team 3

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 101 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Sarah Cione

2 of 7 sc
134

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002345
06/03/2010 09:02:56

Standard Lot#: 642845
Last Changed By: C. STRAUB #50772

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SOAP, SOAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:03:40
Diagnostic Test	Pass	09:04:09
Air Blank	0.000	09:04:36
0.100 Cal Check	0.102	09:04:55
Air Blank	0.000	09:05:24
Subject Test	0.000	09:05:43
Air Blank	0.000	09:06:10
Five Minute Wait	ABT*	
Air Blank	0.000	09:06:45

*Sequence Aborted

Not a Successfully
Completed
Test Sequence

3 of 7 sc
B4

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002345
06/03/2010 09:06:53

Standard Lot#: 642845
Last Changed By: C. STRAUB #50772

GAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SOAP, SOAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:07:41
Diagnostic Test	Pass	09:08:09
Air Blank	RFI*	09:08:22
Air Blank	IPS**	09:08:31

*RFI Detect
**Improper Sample

Not a Successfully
Completed
Test Sequence

4 of 7 sc
BY

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002345
06/03/2010 09:08:39

Standard Lot#: 642845
Last Changed By: C. STRAUB #50772

OAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SOAP, SOAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:09:23
Diagnostic Test	Pass	09:09:52
Air Blank	0.000	09:10:19
0.100 Cal Check	0.102	09:10:38
Air Blank	0.000	09:11:08
Subject Test	0.074*	09:14:16
Air Blank	0.000	09:14:47
Air Blank	0.000	09:15:15
Subject Test	XXX**	09:15:35
Air Blank	0.000	09:16:07

*Deficient Sample
**Invalid Sample
Begin new deprivation period

Not a Successfully
Completed
Test Sequence

5 of 7sc

154

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002345
06/03/2010 09:16:27

Standard Lot#: 642845
Last Changed By: C. STRAUB #50772

OAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:16:54
Cal Check	0.101	09:17:13
Air Blank	0.000	09:17:43

Sarah Cione
Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

6 of 7 SC

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000
STANDARD CALIBRATION CHECK PROCEDURE

TS4
6-7-2010

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 6/3/10 TIME 900
Intoxilyzer Serial # 80-002345 Location Team 3

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 102 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

7 of 7sc

B4

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002345
06/03/2010 09:00:57

Standard Lot#: 642845
Last Changed By: C. STRAUB #50772

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:01:25
Cal Check	0.102	09:01:43
Air Blank	0.000	09:02:13

Sarah Cione
Operator's Signature
