



CERTIFICATION by CUSTODIAN of RECORDS  
 Tucson Police Department Crime Laboratory  
 270 S. Stone Avenue Tucson, AZ 85701-1917

TS6  
 9-2-2010

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002345 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 08/26/10 Cal 1 0.100 AC Cal 2 0.099 AC

Between the dates of 07/29/10 and 08/26/10 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.  
 The instrument was returned to service on N/A.

Cheryl K. Straub

Cheryl K. Straub


08/31/10

ADPS QAS Permit # 33408

SUBSCRIBED AND SWORN to before me this 31, day of August, 2010.

A.M. Hernandez

Notary Public

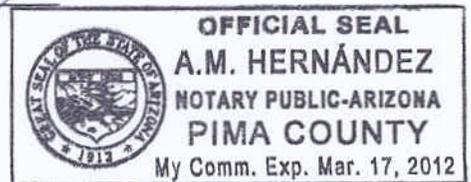

\*Includes documents specified in A.A.C. R13-10-104.A.5

187  
CS

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000  
STANDARD QUALITY ASSURANCE PROCEDURE

TS4  
9-2-2010

QA SPECIALIST Cheryl STRAUB AGENCY Tucson Police Department  
DATE 08/26/10 TIME 07:00  
INTOXILYZER SERIAL # 80-002345 LOCATION TEAM 3

(  ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (  ) 1. Clock time check.
- (  ) 2. Date check.

OPERATIONAL TESTS

- (  ) 1. Alcohol-free subject test result 0. 000 AC.
- (  ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- (  ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (  ) 4. Standard Calibration Check standard 0. 100 AC.  
Result: 0. 099 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE Cheryl Straub

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

287  
CS

TS4

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002345  
08/26/2010 07:02:31

Standard Lot#: 642845  
Last Changed By: C. STRAUB #50772

QAS: C. STRAUB #50772  
TUCSON PD

Operator: C. STRAUB #50772  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	07:03:19
Diagnostic Test	Pass	07:03:46
Air Blank	0.000	07:04:13
0.100 Cal Check	0.100	07:04:32
Air Blank	0.000	07:05:02
Subject Test	0.000	07:05:34
Air Blank	ABT*	07:05:41
Air Blank	0.000	07:06:09

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

317  
CS

154

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002345  
08/26/2010 07:06:20

---

Standard Lot#: 642845  
Last Changed By: C. STRAUB #50772

---

QAS: C. STRAUB #50772  
TUCSON PD

---

Operator: C. STRAUB #50772  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	RFI*	07:06:51
Air Blank	IPS**	07:07:01

\*RFI Detect  
\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

497  
CS

TH

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002345  
08/26/2010

07:07:08

Standard Lot#: 642845  
Last Changed By: C. STRAUB #50772

QAS: C. STRAUB #50772  
TUCSON PD

Operator: C. STRAUB #50772  
TPD

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	07:08:17
Diagnostic Test	Pass	07:08:45
Air Blank	0.000	07:09:12
0.100 Cal Check	0.101	07:09:31
Air Blank	0.000	07:10:00
Subject Test	0.126*	07:13:09
Air Blank	0.000	07:13:41
Air Blank	0.000	07:14:08
Subject Test	XXX**	07:14:30
Air Blank	0.000	07:15:04

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

577  
CS

TS4

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-002345  
08/26/2010 07:15:30

-----  
Standard Lot#: 642845  
Last Changed By: C. STRAUB #50772

-----  
QAS: C. STRAUB #50772  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	07:15:57
Cal Check	0.099	07:16:16
Air Blank	0.000	07:16:45

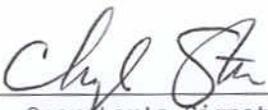
  
-----  
Operator's Signature

EXHIBIT G-2

697  
CS

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

TS4  
9-2-2010

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Cheryl STRAUB AGENCY Tucson Police Department  
DATE 08/26/10 TIME 07:00  
Intoxilyzer Serial # 80-002345 Location TEAM 3

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution N/A AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Cheryl Straub

DPS Form Exh G-2 (Rev. 05-01)

707  
CS


INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-002345  
08/26/2010 07:00:55

TSG

-----  
Standard Lot#: 642845  
Last Changed By: C. STRAUB #50772  
-----

QAS: C. STRAUB #50772  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	07:01:22
Cal Check	0.100	07:01:41
Air Blank	0.000	07:02:11

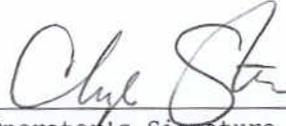

Operator's Signature