



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002345 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

By
9.27.2010

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 09/23/10 Cal 1 0.101 AC Cal 2 0.099 AC

Between the dates of 08/26/10 and 09/23/10 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.
The instrument was returned to service on N/A.

Cheryl K. Straub

Cheryl K. Straub


09/27/10

ADPS QAS Permit # 33408

SUBSCRIBED AND SWORN to before me this 27, day of September, 2010.

A.M. Hernandez

Notary Public

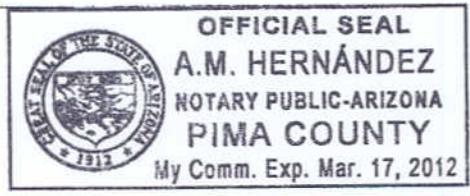

*Includes documents specified in A.A.C. R13-10-104.A.5

178 CS

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000
STANDARD QUALITY ASSURANCE PROCEDURE

TSG
9.27.2010

QA SPECIALIST Cheryl STRAUB AGENCY Tucson Police Department
DATE 09/23/10 TIME 06:40
INTOXILYZER SERIAL # 80-002345 LOCATION TEAM 3

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 099 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

RFE during deficient Sample Test, TEST
REPEATED!

SIGNATURE

Cheryl Straub

278
CS

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002345
09/23/2010

06:41:32

134

Standard Lot#: 642845
Last Changed By: C. STRAUB #50772

QAS: C. STRAUB #50772
TUCSON PD

Operator: C. STRAUB #50772
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	06:42:40
Diagnostic Test	Pass	06:43:07
Air Blank	0.000	06:43:34
0.100 Cal Check	0.101	06:43:53
Air Blank	0.000	06:44:23
Subject Test	0.000	06:44:42
Air Blank	ABT*	06:44:50
Air Blank	0.000	06:45:18

*Sequence Aborted

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

318
CS

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002345
09/23/2010

06:45:24

154

Standard Lot#: 642845
Last Changed By: C. STRAUB #50772

QAS: C. STRAUB #50772
TUCSON PD

Operator: C. STRAUB #50772
TUCSON PD

Subject: SOAP, SOAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	06:45:59
Air Blank	IPS**	06:46:09

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

498
ES

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002345
09/23/2010

06:46:15

TS4

Standard Lot#: 642845
Last Changed By: C. STRAUB #50772

QAS: C. STRAUB #50772
TUCSON PD

Operator: C. STRAUB #50772
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	06:47:04
Diagnostic Test	Pass	06:47:32
Air Blank	RFI*	06:47:42
Air Blank	0.000	06:48:11

*RFI Detect

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

528
CS

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002345
09/23/2010

06:48:17

TS4

Standard Lot#: 642845
Last Changed By: C. STRAUB #50772

QAS: C. STRAUB #50772
TUCSON PD

Operator: C. STRAUB #50772
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	06:49:10
Diagnostic Test	Pass	06:49:38
Air Blank	0.000	06:50:05
0.100 Cal Check	0.100	06:50:24
Air Blank	0.000	06:50:53
Subject Test	0.100*	06:54:02
Air Blank	0.000	06:54:34
Air Blank	0.000	06:55:02
Subject Test	XXX**	06:55:29
Air Blank	0.000	06:56:02

*Deficient Sample
**Invalid Sample
Begin new deprivation period

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

698
CS


INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002345
09/23/2010 06:56:20

TS4

Standard Lot#: 642845
Last Changed By: C. STRAUB #50772

QAS: C. STRAUB #50772
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	06:56:48
Cal Check	0.099	06:57:06
Air Blank	0.000	06:57:36

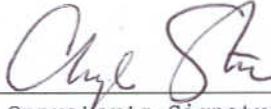

Operator's Signature

EXHIBIT G-2

788
esTHIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000
STANDARD CALIBRATION CHECK PROCEDURETS4
9.27.2010QA Specialist Cheryl STRAUB AGENCY Tucson Police Department
DATE 09/23/10 TIME 06:40
Intoxilyzer Serial # 80-002345 Location TEAM 3

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution N/A AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 101 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Cheryl Straub

DPS Form Exh G-2 (Rev. 05-01)

8/8
CS


INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002345
09/23/2010 06:40:06

TSy

Standard Lot#: 642845
Last Changed By: C. STRAUB #50772

QAS: C. STRAUB #50772
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	06:40:33
Cal Check	0.101	06:40:52
Air Blank	0.000	06:41:22

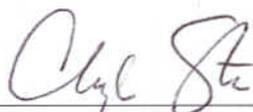

Operator's Signature