



CERTIFICATION by CUSTODIAN of RECORDS
 Tucson Police Department Crime Laboratory
 270 S. Stone Avenue Tucson, AZ 85701-1917

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002345 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

154
 10.25.2010

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 10/21/10 Cal 1 0.101 AC Cal 2 0.100 AC

Between the dates of 9/23/10 and 10/21/10 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A
 The instrument was returned to service on N/A

Cheryl K. Straub

Cheryl K. Straub


10/21/10

ADPS QAS Permit # 33408

SUBSCRIBED AND SWORN to before me this 21, day of October, 20 10.

Amy Hernandez

Notary Public

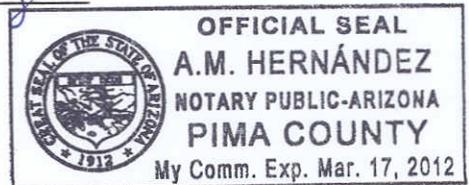

*Includes documents specified in A.A.C. R13-10-104.A.5

187 es

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

TS4
10-25-2010

QA SPECIALIST Cheryl STRAUB AGENCY Tucson Police Department

DATE 10/21/10 TIME 06:58

INTOXILYZER SERIAL # 80-002345 LOCATION TEAM 3

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 100 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Cheryl Straub

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

237
CS

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002345
10/21/2010 07:00:08

134

Standard Lot#: 642845
Last Changed By: C. STRAUB #50772

QAS: C. STRAUB #50772
TUCSON PD

Operator: C. STRAUB #50772
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	07:00:59
Diagnostic Test	Pass	07:01:26
Air Blank	0.000	07:01:54
0.100 Cal Check	0.101	07:02:12
Air Blank	0.000	07:02:42
Subject Test	0.000	07:03:05
Air Blank	ABT*	07:03:12
Air Blank	0.000	07:03:41

*Sequence Aborted

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

367
CS

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002345
10/21/2010

07:03:46

B4

Standard Lot#: 642845
Last Changed By: C. STRAUB #50772

QAS: C. STRAUB #50772
TUCSON PD

Operator: C. STRAUB #50772
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	07:04:25
Air Blank	IPS**	07:04:35

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

477
CS

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002345
10/21/2010

07:04:45

TS4

Standard Lot#: 642845
Last Changed By: C. STRAUB# #50772

QAS: C. STRAUB #50772
TUCSON PD

Operator: C. STRAUB #50772
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	07:05:38
Diagnostic Test	Pass	07:06:05
Air Blank	0.000	07:06:32
0.100 Cal Check	0.101	07:06:51
Air Blank	0.000	07:07:21
Subject Test	0.021*	07:10:29
Air Blank	0.000	07:11:02
Air Blank	0.000	07:11:29
Subject Test	XXX**	07:11:51
Air Blank	0.000	07:12:25

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

577
CS


INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002345
10/21/2010 07:12:42

SH

Standard Lot#: 642845
Last Changed By: C. STRAUB #50772

QAS: C. STRAUB #50772
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	07:13:10
Cal Check	0.100	07:13:28
Air Blank	0.000	07:13:58

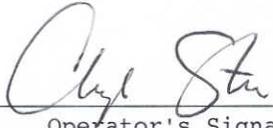

Operator's Signature

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

TSY
10-25-2010

QA Specialist Cheryl STRAUB AGENCY Tucson Police Department
DATE 10/21/10 TIME 06:58
Intoxilyzer Serial # 80-002345 Location TEAM 3

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution N/A AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 101 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Cheryl Straub

DPS Form Exh G-2 (Rev. 05-01)

787
CS


INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002345
10/21/2010 06:58:40

BY

Standard Lot#: 642845
Last Changed By: C. STRAUB #50772

QAS: C. STRAUB #50772
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	06:59:08
Cal Check	0.101	06:59:27
Air Blank	0.000	06:59:56

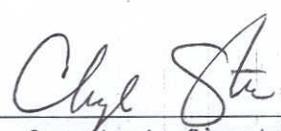

Operator's Signature