



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
270 S. Stone Avenue Tucson, AZ 85701-1917

B4  
2-10-2011

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002345 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 2/9/11 Cal 1 0.102 AC Cal 2 0.102 AC

Between the dates of - and - the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.  
The instrument was returned to service on 2/9/11.

Sarah Cione

Sarah L. Cione


2/10/11

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 10, day of February, 20 11.

A.M. Hernandez

Notary Public

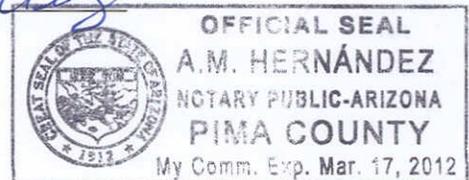

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

1 of 8 JC

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

BH  
2-10-2011

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 2/9/11 TIME 1009

INTOXILYZER SERIAL # 80-002345 LOCATION Crime Lab

( ✓ ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- ( ✓ ) 1. Clock time check.
- ( ✓ ) 2. Date check.

OPERATIONAL TESTS

- ( ✓ ) 1. Alcohol-free subject test result 0. 000 AC.
- ( ✓ ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- ( ✓ ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- ( ✓ ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 102 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

2 of 8 sc

104

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002345  
02/09/2011 10:11:15

Standard Lot#: 20910100A1  
Last Changed By: C. STRAUB #50772

QAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:12:00
Diagnostic Test	Pass	10:12:27
Air Blank	0.000	10:12:54
0.100 Cal Check	0.101	10:13:13
Air Blank	0.000	10:13:43
Subject Test	0.000	10:14:16
Air Blank	0.000	10:14:43
Five Minute Wait	ABT*	
Air Blank	0.000	10:15:35

\*Sequence Aborted

### Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

3 of 8 SC

Bh

INTOXILYZER 8000  
 Location: ARIZONA  
 Serial Number : 80-002345  
 02/09/2011 10:15:40

---

Standard Lot#: 20910100A1  
 Last Changed By: C. STRAUB #50772

---

QAS: S. CIONE #50616  
 TUCSON PD

---

Operator: S. CIONE #50616  
 TUCSON PD

---

Subject: SQAP, SQAP  
 DOB: 11/11/1111  
 Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:16:28
Diagnostic Test	Pass	10:16:55
Air Blank	0.000	10:17:22
0.100 Cal Check	0.102	10:17:41
Air Blank	RFI*	10:17:50
Air Blank	IPS**	10:17:59

\*RFI Detect  
 \*\*Improper Sample

**Not a Successfully Completed Test Sequence**

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

4 of 8 sc

124

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002345  
02/09/2011 10:18:07

Standard Lot#: 20910100A1  
Last Changed By: C. STRAUB #50772

QAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:18:50
Diagnostic Test	Pass	10:19:17
Air Blank	0.000	10:19:45
0.100 Cal Check	0.102	10:20:03
Air Blank	0.000	10:20:33
Subject Test	0.134*	10:23:41
Air Blank	0.000	10:24:13
Air Blank	0.000	10:24:40
Subject Test	XXX**	10:24:56
Air Blank	0.000	10:25:29

\*Deficient Sample  
\*\*Invalid Sample  
Begin new deprivation period

**Not a Successfully Completed Test Sequence**

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-002345  
02/09/2011 10:25:46

5 of 8 SC

BH

-----  
Standard Lot#: 20910100A1  
Last Changed By: C. STRAUB.#50772  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	10:26:14
Cal Check	0.102	10:26:32
Air Blank	0.000	10:27:02

Sarah Cione  
Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

6 of 8 sc

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

TS4  
2.10.2011

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department  
DATE 2/9/11 TIME 1009  
Intoxilyzer Serial # 80-002345 Location Crime Lab

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 102 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-002345  
02/09/2011 10:09:53

7 of 8 sc

BS4

-----  
Standard Lot#: 20910100A1  
Last Changed By: C. STRAUB #50772  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	10:10:21
Cal Check	0.102	10:10:39
Air Blank	0.000	10:11:09

*Sarah Cione*

Operator's Signature



TUCSON POLICE DEPARTMENT  
CRIME LABORATORY

8 of 8 JC

INTOXILYZER SERVICE RECORD

Tbl  
2-10-2011

QAS: Sarah Cione

Intoxilyzer Serial Number: 80-002345

Date Returned to Service: 2/9/2011

- Breath hose was not heating up to acceptable temperature. Replaced breath hose, which then heated to proper temperature. Did not have to send instrument to CMI for repair.