



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
270 S. Stone Avenue Tucson, AZ 85701-1917

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002345 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

TS4  
5.31.2011

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 5/26/11 Cal 1 0.101 AC Cal 2 0.101 AC

Between the dates of 5/5/11 and 5/26/11  
the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -  
The instrument was returned to service on -


Sarah Cione

Sarah L. Cione

5/30/11

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 30 day of MAY, 2011.

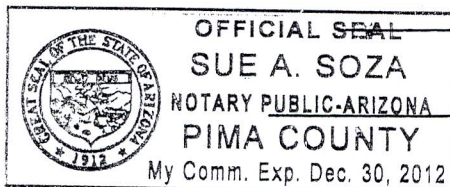

Sue A Soza

Notary Public

\*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department  
DATE 5/26/11 TIME 938  
INTOXILYZER SERIAL # 80-002345 LOCATION DUI HQ

(✓) 1. Display reads "PUSH BUTTON TO START"

## DIAGNOSTIC TESTS

(✓) 1. Clock time check.  
(✓) 2. Date check.

## OPERATIONAL TESTS

(✓) 1. Alcohol-free subject test result 0. 000 AC.  
(✓) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed  
(✓) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.  
(✓) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 101 AC

Instrument is operating properly and accurately. YES ✓ NO       

## COMMENTS

\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE Sarah Cione

2 of 7sc

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002345  
05/26/2011

09:40:50

TS4

Standard Lot#: 20910100A1  
Last Changed By: C. STRAUB #50772

QAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:41:32
Diagnostic Test	Pass	09:41:59
Air Blank	0.000	09:42:27
0.100 Cal Check	0.101	09:42:45
Air Blank	0.000	09:43:15
Subject Test	0.000	09:43:40
Air Blank	0.000	09:44:08
Five Minute Wait	ABT*	
Air Blank	0.000	09:44:47

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002345  
05/26/2011

09:44:52

Standard Lot#: 20910100A1  
Last Changed By: C. STRAUB #50772

QAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	09:45:15
Air Blank	IPS**	09:45:23

\*RFI Detect

\*\*Improper Sample

**Not a Successfully Completed Test Sequence**



4 of 7sc

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002345  
05/26/2011

09:46:10

Standard Lot#: 20910100A1  
Last Changed By: C. STRAUB #50772

QAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:46:53
Diagnostic Test	Pass	09:47:20
Air Blank	0.000	09:47:48
0.100 Cal Check	0.102	09:48:06
Air Blank	0.000	09:48:36
Subject Test	0.010*	09:51:44
Air Blank	0.000	09:52:16
Air Blank	0.000	09:52:43
Subject Test	XXX**	09:53:06
Air Blank	0.000	09:53:39

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000

Location: ARIZONA

Serial Number: 80-002345

05/26/2011

09:53:55

-----  
Standard Lot#: 20910100A1

Last Changed By: C. STRAUB #50772  
-----

QAS: S. CIONE #50616

TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	09:54:23
Cal Check	0.101	09:54:42
Air Blank	0.000	09:55:11

*Sarah Cione*

Operator's Signature

5 of 7 sc

TS4

EXHIBIT G-2

6 of 7 sc

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

T34  
5.31.2011

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department  
DATE 5/26/11 TIME 938  
Intoxilyzer Serial # 80-002345 Location DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach  $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 101 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000

Location: ARIZONA

Serial Number: 80-002345

05/26/2011

09:38:26

-----  
Standard Lot#: 20910100A1

Last Changed By: C. STRAUB #50772  
-----

QAS: S. CIONE #50616

TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	09:38:53
Cal Check	0.101	09:39:12
Air Blank	0.000	09:39:42

*Sarah Cione*

\_\_\_\_\_  
Operator's Signature