



# CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
270 S. Stone Avenue Tucson, AZ 85701-1917

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-00 2345 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

BH  
7.25.2011

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 7-21-11 Cal 1 0.099 AC Cal 2 0.099 AC

Between the dates of 6-23-11 and 7-21-11 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.  
The instrument was returned to service on N/A.

Clare Conley

Clare L. Conley


7-21-11

ADPS QAS Permit # 33314

SUBSCRIBED AND SWORN to before me this 21<sup>st</sup>, day of July, 20 11.

Susan A. Hitchye

Notary Public

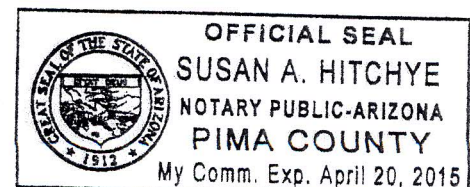

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1567  
Cn  
By  
7.25.2011

QA SPECIALIST Clare Conley AGENCY Tucson Police Department

DATE 7-21-11 TIME 0955

INTOXILYZER SERIAL # 80-002345 LOCATION ODW DUI HQ

( ✓ ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

( ✓ ) 1. Clock time check.

( ✓ ) 2. Date check.

OPERATIONAL TESTS

( ✓ ) 1. Alcohol-free subject test result 0. 000 AC.

( ✓ ) 2. Error recognition logic system functioning.

Not a Successfully Completed Test Sequence printed

( ✓ ) 3. Proper sample recognition system.

Not a Successfully Completed Test Sequence printed

Deficient sample printed.

( ✓ ) 4. Standard Calibration Check standard 0. 100 AC

Result: 0. 099 AC

Instrument is operating properly and accurately. YES ✓ NO       

COMMENTS None.  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE Clare Conley

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-002345  
07/21/2011 10:14:05

2067  
m

-----  
Standard Lot#: 20910100A1  
Last Changed By: C. STRAUB #50772  
-----

B4

-----  
QAS: C. CONLEY #34517  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	10:14:33
Cal Check	0.099	10:14:52
Air Blank	0.000	10:15:21

Clare Conley  
Operator's Signature

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002345  
07/21/2011

09:57:28

Standard Lot#: 20910100A1  
Last Changed By: C. STRAUB #50772

QAS: C. CONLEY #34517  
TUCSON PD

Operator: C. CONLEY #34517  
TUCSON PD

Subject: TEST, SQAP  
DOB: 11/11/1111  
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:58:51
Diagnostic Test	Pass	09:59:19
Air Blank	0.000	09:59:46
0.100 Cal Check	0.099	10:00:05
Air Blank	0.000	10:00:34
Subject Test	0.000	10:01:06
Air Blank	ABT*	10:01:13
Air Blank	0.000	10:01:41

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**



INTOXILYZER 8000

Location: ARIZONA

Serial Number : 80-002345

07/21/2011

10:02:35

Standard Lot#: 20910100A1

Last Changed By: C. STRAUB

#50772

QAS: C. CONLEY

#34517

TUCSON PD

Operator: C. CONLEY

#34517

TUCSON PD

Subject: TEST, SQAP

DOB: 11/11/1111

Sex: M

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:03:31
Diagnostic Test	Pass	10:03:58
Air Blank	0.000	10:04:26
0.100 Cal Check	0.100	10:04:44
Air Blank	RFI*	10:04:56
Air Blank	IPS**	10:05:04

\*RFI Detect

\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002345  
07/21/2011

10:05:12

Standard Lot#: 20910100A1  
Last Changed By: C. STRAUB #50772

QAS: C. CONLEY #34517  
TUCSON PD

Operator: C. CONLEY #34517  
TUCSON PD

Subject: TEST, SQAP  
DOB: 11/11/1111  
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:06:03
Diagnostic Test	Pass	10:06:30
Air Blank	0.000	10:06:57
0.100 Cal Check	0.099	10:07:16
Air Blank	0.000	10:07:45
Subject Test	0.259*	10:10:54
Air Blank	0.000	10:11:26
Air Blank	0.000	10:11:54
Subject Test	XXX**	10:13:14
Air Blank	0.000	10:13:47

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

TS4  
7-25-2011

QA Specialist Clare Conley AGENCY Tucson Police Department  
DATE 7-21-11 TIME 0955  
Intoxilyzer Serial # 80-002345 Location ODW DUI HQ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach  $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 099 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Clare Conley

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-002345  
07/21/2011

09:55:52

-----  
Standard Lot#: 20910100A1  
Last Changed By: C. STRAUB #50772  
-----

QAS: C. CONLEY #34517  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	09:56:19
Cal Check	0.099	09:56:38
Air Blank	0.000	09:57:07

Clare Conley  
Operator's Signature

7567  
on

154