



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002345 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

TS4
8.2.2011

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 07/28/11 Cal 1 0.099 AC Cal 2 0.099 AC

Between the dates of 07/21/11 and 07/28/11
the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A
The instrument was returned to service on N/A

Cheryl K. Straub

Cheryl K. Straub


08/01/11

ADPS QAS Permit # 33408

SUBSCRIBED AND SWORN to before me this 1, day of August, 2011.

Sue A. Soza

Notary Public

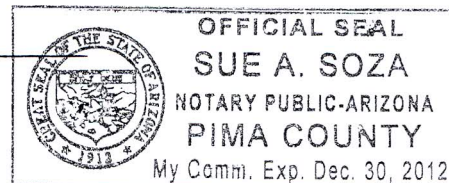

*Includes documents specified in A.A.C. R13-10-104.A.5

187 CS

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

B4
8-2-2011

QA SPECIALIST Cheryl STRAUB AGENCY Tucson Police Department
DATE 07/28/2011 TIME 09:03
INTOXILYZER SERIAL # 80-002345 LOCATION TEAMS

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

(✓) 1. Clock time check.
(✓) 2. Date check.

OPERATIONAL TESTS

(✓) 1. Alcohol-free subject test result 0. 000 AC.
(✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
(✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
(✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 099 AC

Instrument is operating properly and accurately. YES ✓ NO

COMMENTS

SIGNATURE

Cheryl Straub

287
CS

INTOXILYZER 8000

Location: ARIZONA

Serial Number : 80-002345

07/28/2011

09:05:29

154

Standard Lot#: 20910100A1

Last Changed By: C. STRAUB

#50772

QAS: C. STRAUB

#50772

TUCSON PD

Operator: C. STRAUB

#50772

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:06:19
Diagnostic Test	Pass	09:06:47
Air Blank	0.000	09:07:14
0.100 Cal Check	0.099	09:07:33
Air Blank	0.000	09:08:02
Subject Test	0.000	09:08:30
Air Blank	ABT*	09:08:36
Air Blank	0.000	09:09:05

*Sequence Aborted

Not a Successfully Completed Test Sequence

327
CS

INTOXILYZER 8000

Location: ARIZONA

Serial Number : 80-002345

07/28/2011

09:09:14

154

Standard Lot#: 20910100A1

Last Changed By: C. STRAUB

#50772

QAS: C. STRAUB

#50772

TUCSON PD

Operator: C. STRAUB

#50772

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	09:10:03
Air Blank	IPS**	09:10:12

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

487
CS

INTOXILYZER 8000

Location: ARIZONA

Serial Number : 80-002345

07/28/2011

09:10:46

BM

Standard Lot#: 20910100A1

Last Changed By: C. STRAUB

#50772

QAS: C. STRAUB

#50772

TUCSON PD

Operator: C. STRAUB

#50772

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1111

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	09:11:37
Diagnostic Test	Pass	09:12:05
Air Blank	0.000	09:12:32
0.100 Cal Check	0.099	09:12:51
Air Blank	0.000	09:13:20
Subject Test	0.181*	09:16:29
Air Blank	0.000	09:17:01
Air Blank	0.000	09:17:28
Subject Test	XXX**	09:17:52
Air Blank	0.000	09:18:26

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

507 CS

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002345
07/28/2011

09:18:44

BH

Standard Lot#: 20910100A1
Last Changed By: C. STRAUB #50772

QAS: C. STRAUB #50772
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:19:12
Cal Check	0.099	09:19:30
Air Blank	0.000	09:20:00



Operator's Signature

677
CS

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

By
8.2.2011

QA Specialist Cheryl STRAUB AGENCY Tucson Police Department
DATE 07/28/2011 TIME 09:03
Intoxilyzer Serial # 80-002345 Location TEAM 5

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution N/A AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 0.99 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Cheryl Straub

DPS Form Exh G-2 (Rev. 05-01)

72 +
CS

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002345
07/28/2011

09:03:55

BH

Standard Lot#: 20910100A1
Last Changed By: C. STRAUB #50772

QAS: C. STRAUB #50772
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	09:04:22
Cal Check	0.099	09:04:41
Air Blank	0.000	09:05:11


Operator's Signature