



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

TS4
12-21-2011

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002345 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 12-16-11 Cal 1 0.100 AC Cal 2 0.100 AC

Between the dates of 11-17-11 and 12-16-11 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.
The instrument was returned to service on N/A.


Clare Conley

Clare L. Conley

12-16-11

ADPS QAS Permit # 33314

SUBSCRIBED AND SWORN to before me this 16, day of December, 2011.

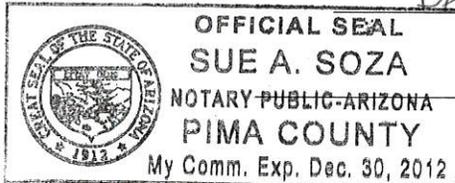

Sue A Soza

Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1 of 7

cm

154
12-21-2011

QA SPECIALIST Clare Conley AGENCY Tucson Police Department

DATE 12-16-11 TIME 0725

INTOXILYZER SERIAL # 80-002345 LOCATION DD

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

(✓) 1. Clock time check.

(✓) 2. Date check.

OPERATIONAL TESTS

(✓) 1. Alcohol-free subject test result 0. 000 AC.

(✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed

(✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.

(✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 100 AC

Instrument is operating properly and accurately. YES NO

COMMENTS Added three minutes to the internal clock.

SIGNATURE Clare Conley

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002345
12/16/2011

07:27:02

2 of 7
CM
BM

Standard Lot#: 20910100A1
Last Changed By: C. STRAUB #50772

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	07:27:47
Diagnostic Test	Pass	07:28:15
Air Blank	0.000	07:28:42
0.100 Cal Check	0.101	07:29:01
Air Blank	0.000	07:29:30
Subject Test	0.000	07:29:47
Air Blank	ABT*	07:29:54
Air Blank	RFI**	07:30:01

*Sequence Aborted
**RFI Detect

Not a Successfully Completed Test Sequence

3567
CN

104

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002345
12/16/2011

07:30:05

Standard Lot#: 20910100A1
Last Changed By: C. STRAUB #50772

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	IPS*	07:30:31
Air Blank	0.000	07:30:59

*Improper Sample

Not a Successfully Completed Test Sequence

4087
cu

P4

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002345
12/16/2011 07:31:11

Standard Lot#: 20910100A1
Last Changed By: C. STRAUB #50772

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	07:31:56
Diagnostic Test	Pass	07:32:23
Air Blank	0.000	07:32:50
0.100 Cal Check	0.100	07:33:09
Air Blank	0.000	07:33:39
Subject Test	0.107*	07:36:47
Air Blank	0.000	07:37:19
Air Blank	0.000	07:37:47
Subject Test	XXX**	07:38:11
Air Blank	0.000	07:38:44

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

5-4-07
cm
cm

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002345
12/16/2011 07:39:00

104

Standard Lot#: 20910100A1
Last Changed By: C. STRAUB #50772

QAS: C. CONLEY #34517
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	07:39:27
Cal Check	0.100	07:39:46
Air Blank	0.000	07:40:15

Clare Conley

Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6 ~~5017~~ ^{cu}
cu

TSY
12.21.2011

QA Specialist Clare Conley AGENCY Tucson Police Department
DATE 12-16-11 TIME 0725
Intoxilyzer Serial # 80-002345 Location ODJ

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Clare Conley

DPS Form Exh G-2 (Rev. 05-01)

7
C#34517
m
124

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002345
12/16/2011 07:25:40

Standard Lot#: 20910100A1
Last Changed By: C. STRAUB #50772

QAS: C. CONLEY #34517
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	07:26:07
Cal Check	0.100	07:26:26
Air Blank	0.000	07:26:55



Operator's Signature