



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
1306 W. Miracle Mile, Tucson, AZ 85705

TS4  
4-11-2012

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002345 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 4/5/12 Cal 1 0.102 AC Cal 2 0.101 AC

Between the dates of 3/9/12 and 4/5/12 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.  
The instrument was returned to service on -.

Sarah Cione

Sarah L. Cione


4/9/12

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 9th, day of April, 2012.

Elma Escarcega

Notary Public

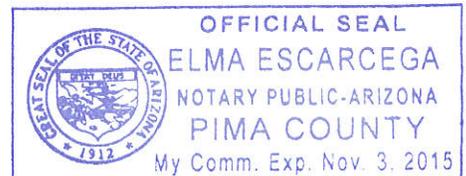

\*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

1 of 7sc

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

TS4

4.11.2012

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 4/5/12 TIME 1524

INTOXILYZER SERIAL # 80-002345 LOCATION ODD-Tm 5

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 101 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

---



---



---



---

SIGNATURE Sarah Cione

2 of 7<sub>sc</sub>

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002345  
04/05/2012 15:25:54

TS4

Standard Lot#: 20910100A1  
Last Changed By: C. STRAUB #50772

QAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:26:44
Diagnostic Test	Pass	15:27:11
Air Blank	0.000	15:27:38
0.100 Cal Check	0.102	15:27:57
Air Blank	0.000	15:28:27
Subject Test	0.000	15:28:46
Air Blank	0.000	15:29:13
Five Minute Wait	ABT*	
Air Blank	0.000	15:29:53

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

3 of 7 sc

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002345  
04/05/2012 15:30:06

B4

Standard Lot#: 20910100A1  
Last Changed By: C. STRAUB #50772

QAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	15:30:38
Air Blank	IPS**	15:30:46

\*RFI Detect  
\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

4 of 7 sc

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number : 80-002345  
04/05/2012 15:30:51

134

Standard Lot#: 20910100A1  
Last Changed By: C. STRAUB #50772

QAS: S. CIONE #50616  
TUCSON PD

Operator: S. CIONE #50616  
TUCSON PD

Subject: SQAP, SQAP  
DOB: 11/11/1111  
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	15:31:36
Diagnostic Test	Pass	15:32:04
Air Blank	0.000	15:32:31
0.100 Cal Check	0.101	15:32:50
Air Blank	0.000	15:33:20
Subject Test	0.034*	15:36:28
Air Blank	0.000	15:36:59
Air Blank	0.000	15:37:27
Subject Test	XXX**	15:37:45
Air Blank	0.000	15:38:18

\*Deficient Sample  
\*\*Invalid Sample  
Begin new deprivation period

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-002345  
04/05/2012 15:38:38

5 of 7 sc  
B4

-----  
Standard Lot#: 20910100A1  
Last Changed By: C. STRAUB #50772  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	15:39:05
Cal Check	0.101	15:39:24
Air Blank	0.000	15:39:53

*Sarah Cione*

Operator's Signature

EXHIBIT G-2

6 of 7 SC

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

B4  
4-11-2012

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department  
DATE 4/5/12 TIME 1524  
Intoxilyzer Serial # 80-002345 Location ODD-Tm 5

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 102 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

7 of 7sc

INTOXILYZER 8000  
Location: ARIZONA  
Serial Number: 80-002345  
04/05/2012 15:24:28

BL

-----  
Standard Lot#: 20910100A1  
Last Changed By: C. STRAUB #50772  
-----

QAS: S. CIONE #50616  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	15:24:56
Cal Check	0.102	15:25:15
Air Blank	0.000	15:25:44

Sarah Cione  
Operator's Signature