



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

SC
7/9/12

I hereby certify that these documents are a true and correct copy of the
periodic maintenance* for the Intoxilyzer 8000 Serial Number
80-002345 maintained by the Tucson Police Department
Crime Laboratory pursuant to the requirements of the Arizona Department
of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following
dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2,
for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 6-28-2012 Cal 1 0.100 AC Cal 2 0.099 AC

Between the dates of 5-31-2012 and 6-28-2012
the instrument was operating properly and accurately.

The instrument was taken out of service for repair on
The instrument was returned to service on

[Handwritten signature of Terry S. Gallegos]

Terry S. Gallegos

7-3-2012

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 3rd, day of
July, 2012.

[Handwritten signature of Notary Public]

Notary Public



*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1 of 7
TSg
SC
7/9/12

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department
DATE 6-28-2012 TIME 1119
INTOXILYZER SERIAL # 80-002345 LOCATION ODD

() 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- () 1. Clock time check.
- () 2. Date check.

OPERATIONAL TESTS

- () 1. Alcohol-free subject test result 0. 000 AC.
- () 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- () 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- () 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 099 AC

Instrument is operating properly and accurately. YES NO

COMMENTS _____

SIGNATURE Terry Gallegos

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002345
06/28/2012

11:21:09

2077
TS4

SC

Standard Lot#: 23511100A2
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST,SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:22:02
Diagnostic Test	Pass	11:22:29
Air Blank	0.000	11:22:56
0.100 Cal Check	0.099	11:23:15
Air Blank	0.000	11:23:45
Subject Test	0.000	11:24:12
Air Blank	0.000	11:24:40
Five Minute Wait	ABT*	
Air Blank	0.000	11:25:13

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002345
06/28/2012

11:25:43

3087
TS4
SC

Standard Lot#: 23511100A2
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST,SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:26:38
Diagnostic Test	Pass	11:27:05
Air Blank	0.000	11:27:33
0.100 Cal Check	RFI*	11:27:46
Air Blank	IPS**	11:27:56

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002345
06/28/2012 11:28:05

4 of 7
TS4

Standard Lot#: 23511100A2
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

SC

Subject: TEST,SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:28:56
Diagnostic Test	Pass	11:29:23
Air Blank	0.000	11:29:50
0.100 Cal Check	0.099	11:30:09
Air Blank	0.000	11:30:39
Subject Test	0.140*	11:33:47
Air Blank	0.000	11:34:19
Air Blank	0.000	11:34:46
Subject Test	XXX**	11:35:27
Air Blank	0.000	11:36:00

*Deficient Sample
**Invalid Sample
Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002345
06/28/2012 11:36:21


5057
B4

SC

Standard Lot#: 23511100A2
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	11:36:48
Cal Check	0.099	11:37:07
Air Blank	0.000	11:37:37

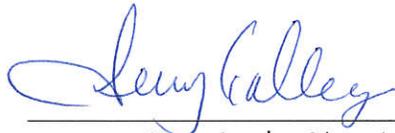

Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6 of 7
TSK
SC
7/9/12

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department
DATE 6.28.2012 TIME 1119
Intoxilyzer Serial # 80-002345 Location ODD

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE 

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002345
06/28/2012 11:19:00


7 of 7
B4

Standard Lot#: 23511100A2
Last Changed By: T. GALLEGOS #28847

50

QAS: T. GALLEGOS #28847
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	11:19:27
Cal Check	0.100	11:19:46
Air Blank	0.000	11:20:16

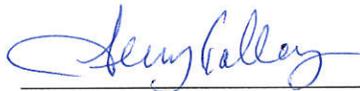

Operator's Signature