



CERTIFICATION by CUSTODIAN of RECORDS
 Tucson Police Department Crime Laboratory
 1306 W. Miracle Mile, Tucson, AZ 85705

SC
7/31/12

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002345 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 7-26-2012 Cal 1 0.100 AC Cal 2 0.099 AC

Between the dates of 6-28-2012 and 7-26-2012 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on _____.
 The instrument was returned to service on _____.

Terry S. Gallegos

Terry S. Gallegos


7-31-2012

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 31, day of July, 2012.

Ann L. Beecroft

Notary Public

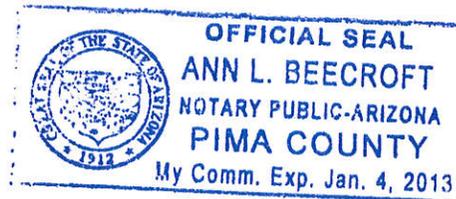

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1 of 7
TS4
SC
7/31/12

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department
DATE 7-26-2012 TIME 1104
INTOXILYZER SERIAL # 80-002345 LOCATION DDD

() 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

() 1. Clock time check.
() 2. Date check.

OPERATIONAL TESTS

() 1. Alcohol-free subject test result 0. 000 AC.
() 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
() 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
() 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 099 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE 

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002345
07/26/2012

11:05:37

2077
TS4

5c

Standard Lot#: 23511100A2
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|------------------|--------|----------|
| Air Blank | 0.000 | 11:06:39 |
| Diagnostic Test | Pass | 11:07:06 |
| Air Blank | 0.000 | 11:07:33 |
| 0.100 Cal Check | 0.099 | 11:07:52 |
| Air Blank | 0.000 | 11:08:22 |
| Subject Test | 0.000 | 11:08:58 |
| Air Blank | 0.000 | 11:09:25 |
| Five Minute Wait | ABT* | |
| Air Blank | 0.000 | 11:10:00 |

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002345
07/26/2012

11:10:18

30f7
TS4

5c

Standard Lot#: 23511100A2
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 11:11:14 |
| Diagnostic Test | Pass | 11:11:42 |
| Air Blank | 0.000 | 11:12:09 |
| 0.100 Cal Check | RFI* | 11:12:23 |
| Air Blank | IPS** | 11:12:33 |

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002345
07/26/2012

11:12:43

40f7
TS4

Standard Lot#: 23511100A2
Last Changed By: T. GALLEGOS #28847

SC

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 11:13:34 |
| Diagnostic Test | Pass | 11:14:01 |
| Air Blank | 0.000 | 11:14:29 |
| 0.100 Cal Check | 0.100 | 11:14:47 |
| Air Blank | 0.000 | 11:15:17 |
| Subject Test | 0.105* | 11:18:25 |
| Air Blank | 0.000 | 11:18:57 |
| Air Blank | 0.000 | 11:19:24 |
| Subject Test | XXX** | 11:19:46 |
| Air Blank | 0.000 | 11:20:19 |

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence


INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002345
07/26/2012 11:20:35

5077
TS4
SC

Standard Lot#: 23511100A2
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

| Test | g/210L | Time |
|-----------|--------|----------|
| Air Blank | 0.000 | 11:21:03 |
| Cal Check | 0.099 | 11:21:21 |
| Air Blank | 0.000 | 11:21:51 |

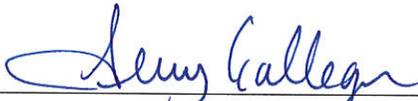

Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6057
TS4
SC
7/31/12

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department
DATE 7-26-2012 TIME 1104
Intoxilyzer Serial # 80-002345 Location ODD

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Terry Gallegos

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002345
07/26/2012 11:04:12

7 of 7
D4

Standard Lot#: 23511100A2
Last Changed By: T. GALLEGOS #28847

5c

QAS: T. GALLEGOS #28847
TUCSON PD

| Test | g/210L | Time |
|-----------|--------|----------|
| Air Blank | 0.000 | 11:04:39 |
| Cal Check | 0.100 | 11:04:58 |
| Air Blank | 0.000 | 11:05:28 |



Operator's Signature