



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

SC
1/14/14

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002345 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 1-9-2014 Cal 1 0.100 AC Cal 2 0.098 AC

Between the dates of 12.12.2013 and 1.9.2014 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on .
The instrument was returned to service on .

Terry S. Gallegos

Terry S. Gallegos


1.13.2014

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 13th, day of January, 2014.

Edma Escarcega

Notary Public

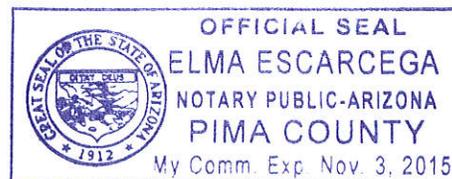

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1 of 7
TS4 SC
1/14/14

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department
DATE 1-9-2014 TIME 1215
INTOXILYZER SERIAL # 80-002345 LOCATION ODD

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS


- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 098 AC

Instrument is operating properly and accurately. YES NO

COMMENTS _____

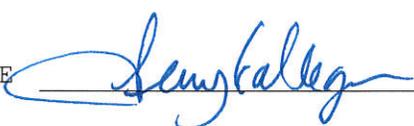SIGNATURE 

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location:
Serial Number: 80-002345
Core Version: 8105.48
01/09/2014 12:17:11

2 of 7
TS4
SC

Standard Lot#: 25913100A2
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	12:18:04
Diagnostic Test	Pass	12:18:39
Air Blank	0.000	12:19:06
0.100 Cal Check	0.100	12:19:27
Air Blank	0.000	12:19:56
Subject Test	0.000	12:20:30
Air Blank	0.000	12:20:57
Wait	ABT*	
Air Blank	0.000	12:21:32

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location:
Serial Number: 80-002345
Core Version: 8105.48
01/09/2014 12:21:55

3057
B4

Standard Lot#: 25913100A2
Last Changed By: T. GALLEGOS #28847

5C

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	12:22:43
Diagnostic Test	Pass	12:23:18
Air Blank	0.000	12:23:46
0.100 Cal Check	RFI*	12:23:58
Air Blank	IPS**	12:24:06

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location:

Serial Number: 80-002345

Core Version: 8105.48

01/09/2014

12:24:10

4 of 7
TS4

Standard Lot#: 25913100A2

Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847

TUCSON PD

Operator: T. GALLEGOS #28847

TUCSON PD

Subject: TEST,SQAP

DOB: 11/11/1987

Sex: F

Weight: 111

5c

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	12:25:07
Diagnostic Test	Pass	12:25:42
Air Blank	0.000	12:26:10
0.100 Cal Check	0.099	12:26:31
Air Blank	0.000	12:27:00
Subject Test	0.229*	12:30:08
Air Blank	0.000	12:30:42
Air Blank	0.000	12:31:10
Subject Test	XXX**	12:31:33
Air Blank	0.000	12:32:06

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location:
Serial Number: 80-002345
Core Version: 8105.48
01/09/2014 12:32:28

50f 7
TB4
5C

Standard Lot#: 25913100A2
Last Changed By: T. GALLEGOS #28847

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	12:32:58
0.100 Cal Check	0.098	12:33:19
Air Blank	0.000	12:33:48

Preventative Maintenance
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6077
TS4
SC
1/14/14

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department
DATE 1-9-2014 TIME 1215
Intoxilyzer Serial # 80-002345 Location ODD

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Terry Gallegos

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

