



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

TS4
3.11.2015

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002345 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 3-6-15 Cal 1 0.102 AC Cal 2 0.101 AC

Between the dates of 2-6-15 and 3-6-15 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.
The instrument was returned to service on N/A.

Amy L. Autret
Amy L. Autret


3-9-15

ADPS QAS Permit # 33718

SUBSCRIBED AND SWORN to before me this Open, day of March, 2015.

Elma Escarcega

Notary Public

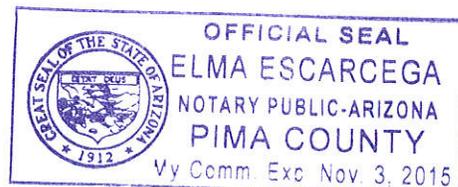

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

10F7
JMO

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

TS4
3.11.2015

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Amy Autret 53546 AGENCY Tucson Police Department

DATE 3-6-15 TIME 1134

INTOXILYZER SERIAL # 80-002345 LOCATION 000

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.


OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 101 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE

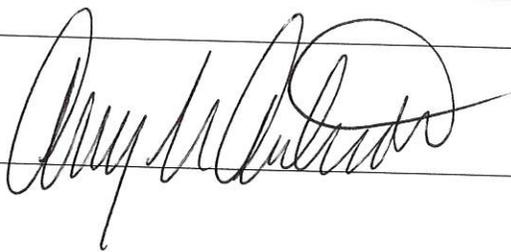

2 of 7
AA
TS4

INTOXILYZER 8000
Location:
Serial Number: 80-002345
Core Version: 8105.48
03/06/2015

11:36:47

Standard Lot#: 25414100A4
Last Changed By: A. AUTRET #53546

QAS: A. AUTRET #53546
TUCSON PD

Operator: A. AUTRET #53546
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:37:45
Diagnostic Test	Pass	11:38:20
Air Blank	0.000	11:38:47
0.100 Cal Check	0.101	11:39:08
Air Blank	0.000	11:39:38
Subject Test	0.000	11:39:59
Air Blank	0.000	11:40:26
Wait	ABT*	
Air Blank	0.000	11:41:01

*Sequence Aborted

Not a Successfully Completed Test Sequence

Pursuant to section 23-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

3 of 7
HA
TSh

INTOXILYZER 8000
Location:
Serial Number: 80-002345
Core Version: 8105.48
03/06/2015 11:41:21

Standard Lot#: 25414100A4
Last Changed By: A. AUTRET #53546

QAS: A. AUTRET #53546
TUCSON PD

Operator: A. AUTRET #53546
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:42:19
Diagnostic Test	Pass	11:42:54
Air Blank	0.000	11:43:22
0.100 Cal Check	RFI*	11:43:34
Air Blank	IPS**	11:43:42

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

Pursuant to section 23-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

4 of 7
ASB
TSh

INTOXILYZER 8000
Location:
Serial Number: 80-002345
Core Version: 8105.48
03/06/2015 11:44:24

Standard Lot#: 25414100A4
Last Changed By: A. AUTRET #53546

QAS: A. AUTRET #53546
TUCSON PD

Operator: A. AUTRET #53546
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:45:27
Diagnostic Test	Pass	11:46:02
Air Blank	0.000	11:46:29
0.100 Cal Check	0.102	11:46:50
Air Blank	0.000	11:47:20
Subject Test	0.178*	11:50:27
Air Blank	0.000	11:50:59
Air Blank	0.000	11:51:27
Subject Test	XXX**	11:51:48
Air Blank	0.000	11:52:22

*Deficient Sample
**Invalid Sample
Begin new deprivation period

Not a Successfully Completed Test Sequence

6 of 7
MA

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

TS4
3.11.2015

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Amy Antret 53546 AGENCY Tucson Police Department
DATE 3-6-15 TIME 1134
Intoxilyzer Serial # 80-002345 Location 000

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 102 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE 

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

7 of 7
AS
Bh

INTOXILYZER 8000
Location:
Serial Number: 80-002345
Core Version: 8105.48
03/06/2015 11:34:34

Standard Lot#: 25414100A4
Last Changed By: A. AUTRET #53546

QAS: A. AUTRET #53546
TUCSON PD

Operator: A. AUTRET #53546
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	11:35:05
0.100 Cal Check	0.102	11:35:26
Air Blank	0.000	11:35:55

Preventative Maintenance
Performed -- Timer Reset