



CERTIFICATION by CUSTODIAN of RECORDS
 Tucson Police Department Crime Laboratory
 1306 W. Miracle Mile, Tucson, AZ 85705

Ab
4/22/15

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002345 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 4.1.2015 Cal 1 0.100 AC Cal 2 0.100 AC

Between the dates of 3.6.2015 and 4.1.2015 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.
 The instrument was returned to service on -.

Terry Gallegos

Terry S. Gallegos


4.16.2015

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 16th, day of April, 2015.

Elma Escarcega

Notary Public

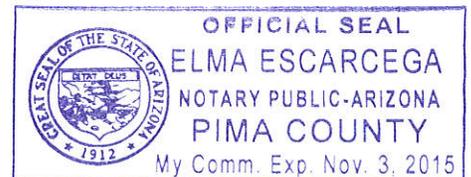

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

10 of 7
TS4
Ab
4/22/15

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department
DATE 4.1.2015 TIME 1041
INTOXILYZER SERIAL # 80-002345 LOCATION ODD

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 100 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE

Terry Gallegos

INTOXILYZER 8000
Location:
Serial Number: 80-002345
Core Version: 8105.48
04/01/2015

10:43:56

2 of 7
754
AD

Standard Lot#: 25414100A4
Last Changed By: A. AUTRET #53546

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST,SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:44:56
Diagnostic Test	Pass	10:45:31
Air Blank	0.000	10:45:58
0.100 Cal Check	0.100	10:46:19
Air Blank	0.000	10:46:48
Subject Test	0.000	10:47:12
Air Blank	0.000	10:47:40
Wait	ABT*	
Air Blank	0.000	10:48:26

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location:
Serial Number: 80-002345
Core Version: 8105.48
04/01/2015

10:48:45

3 of 7
TS4
JW

Standard Lot#: 25414100A4
Last Changed By: A. AUTRET #53546

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Subject: TEST,SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:49:35
Diagnostic Test	Pass	10:50:10
Air Blank	0.000	10:50:37
0.100 Cal Check	RFI*	10:50:52
Air Blank	IPS**	10:51:00

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000

Location:

Serial Number: 80-002345

Core Version: 8105.48

04/01/2015

10:51:04

4 of 7
TS4
AB

Standard Lot#: 25414100A4

Last Changed By: A. AUTRET

#53546

QAS: T. GALLEGOS

#28847

TUCSON PD

Operator: T. GALLEGOS

#28847

TUCSON PD

Subject: TEST,SQAP

DOB: 11/11/1987

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:52:02
Diagnostic Test	Pass	10:52:37
Air Blank	0.000	10:53:05
0.100 Cal Check	0.101	10:53:25
Air Blank	0.000	10:53:55
Subject Test	0.161*	10:57:02
Air Blank	0.000	10:57:34
Air Blank	0.000	10:58:02
Subject Test	XXX**	10:58:36
Air Blank	0.000	10:59:09

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

5087
TS4
JD

INTOXILYZER 8000
Location:
Serial Number: 80-002345
Core Version: 8105.48
04/01/2015 10:59:51

Standard Lot#: 25414100A4
Last Changed By: A. AUTRET #53546

QAS: T. GALLEGOS #28847
TUCSON PD

Operator: T. GALLEGOS #28847
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	11:00:22
0.100 Cal Check	0.100	11:00:43
Air Blank	0.000	11:01:13

Preventative Maintenance
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6077
TS4
AB
4/22/15

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department

DATE 4.1.2015 TIME 1041

Intoxilyzer Serial # 80-002345 Location ODD

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE 

DPS Form Exh G-2 (Rev. 05-01)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000

Location:

Serial Number: 80-002345

Core Version: 8105.48

04/01/2015

10:41:41

Standard Lot#: 25414100A4

Last Changed By: A. AUTRET #53546

QAS: T. GALLEGOS #28847

TUCSON PD

Operator: T. GALLEGOS #28847

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	10:42:20
0.100 Cal Check	0.100	10:42:40
Air Blank	0.000	10:43:10

Preventative Maintenance
Performed -- Timer Reset

70F7
TS4
JAO