



CERTIFICATION by CUSTODIAN of RECORDS
 Tucson Police Department Crime Laboratory
 1306 W. Miracle Mile, Tucson, AZ 85705

TS4
 7-2-2015


I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002345 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 6-26-2015 Cal 1 0.100 AC Cal 2 0.098 AC

Between the dates of 5-29-2015 and 6-26-2015 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on _____.
 The instrument was returned to service on _____.

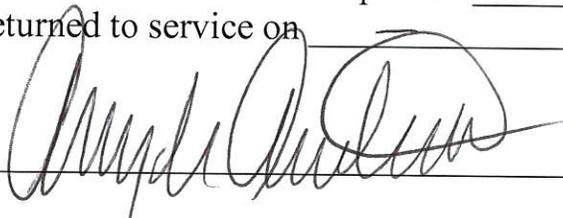

 Amy L. Autret


7/2/15

ADPS QAS Permit # 33718

SUBSCRIBED AND SWORN to before me this 2nd, day of July, 2015.



Notary Public

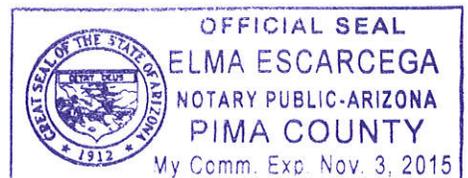

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

1 of 7
Bh
7-2-2015

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Amy Autret 53546 AGENCY Tucson Police Department

DATE 6-26-2015 TIME 1204

INTOXILYZER SERIAL # 80-002345 LOCATION 000

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.


OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 098 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE

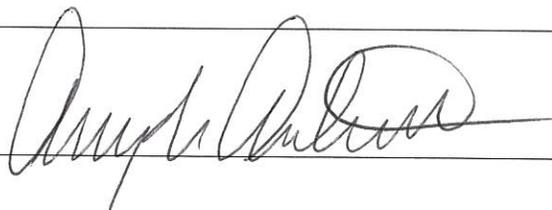

INTOXILYZER 8000

Location:

Serial Number: 80-002345

Core Version: 8105.48

06/26/2015

12:05:44

20F7
JLD

B4

Standard Lot#: 25414100A4

Last Changed By: A. AUTRET #53546

QAS: A. AUTRET #53546

TUCSON PD

Operator: A. AUTRET #53546

TUCSON PD

Subject: SQAP, SQAP

DOB: 11/11/1987

Sex: F

Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	12:06:35
Diagnostic Test	Pass	12:07:10
Air Blank	0.000	12:07:37
0.100 Cal Check	0.100	12:07:58
Air Blank	0.000	12:08:27
Subject Test	0.000	12:08:42
Air Blank	0.000	12:09:10
Wait	ABT*	
Air Blank	0.000	12:09:44

*Sequence Aborted

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location:
Serial Number: 80-002345
Core Version: 8105.48
06/26/2015

12:09:49

30F7
JL
TSH

Standard Lot#: 25414100A4
Last Changed By: A. AUTRET #53546

QAS: A. AUTRET #53546
TUCSON PD

Operator: A. AUTRET #53546
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	12:10:37
Diagnostic Test	Pass	12:11:12
Air Blank	0.000	12:11:39
0.100 Cal Check	RFI*	12:11:52
Air Blank	IPS**	12:12:00

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location:
Serial Number: 80-002345
Core Version: 8105.48
06/26/2015

12:12:06

4 of 7
JL
TG

Standard Lot#: 25414100A4
Last Changed By: A. AUTRET #53546

QAS: A. AUTRET #53546
TUCSON PD

Operator: A. AUTRET #53546
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	12:12:57
Diagnostic Test	Pass	12:13:32
Air Blank	0.000	12:13:59
0.100 Cal Check	0.100	12:14:20
Air Blank	0.000	12:14:49
Subject Test	0.147*	12:17:57
Air Blank	0.000	12:18:29
Air Blank	0.000	12:18:56
Subject Test	XXX**	12:19:20
Air Blank	0.000	12:19:54

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location:
Serial Number: 80-002345
Core Version: 8105.48
06/26/2015 12:20:08

5 of 7
JL
Tsh

Standard Lot#: 25414100A4
Last Changed By: A. AUTRET #53546

QAS: A. AUTRET #53546
TUCSON PD

Operator: A. AUTRET #53546
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	12:20:38
0.100 Cal Check	0.098	12:20:58
Air Blank	0.000	12:21:28

Preventative Maintenance
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

6 of 7
TSJ
7.2.2015

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Amy Autret 53546 AGENCY Tucson Police Department
DATE 6-26-2015 TIME 1204
Intoxilyzer Serial # 80-002345 Location ODA

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE 

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location:
Serial Number: 80-002345
Core Version: 8105.48
06/26/2015 12:04:14

Standard Lot#: 25414100A4
Last Changed By: A. AUTRET #53546

QAS: A. AUTRET #53546
TUCSON PD

Operator: A. AUTRET #53546
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	12:04:43
0.100 Cal Check	0.100	12:05:04
Air Blank	0.000	12:05:34

Preventative Maintenance
Performed -- Timer Reset

7087
JLB

TSh

7-2-2015