



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

TS4
8.25.2015

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002345 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 8/21/2015 Cal 1 0.101 AC Cal 2 0.100 AC

Between the dates of 7/24/2015 and 8/21/2015 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.
The instrument was returned to service on -.

Amy L. Autret

Amy L. Autret


8/24/2015

ADPS QAS Permit # 33718

SUBSCRIBED AND SWORN to before me this 24th, day of August, 2015.

Elma Escarcega

Notary Public

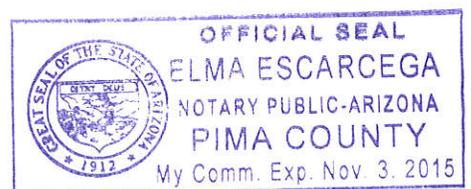

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

1 of 8
JLS
TSG
8.25.2015

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Amy Autret 53546 AGENCY Tucson Police Department

DATE 8/21/2015 TIME 1333

INTOXILYZER SERIAL # 80-002345 LOCATION 000

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS


- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 100 AC

Instrument is operating properly and accurately. YES NO

COMMENTS _____

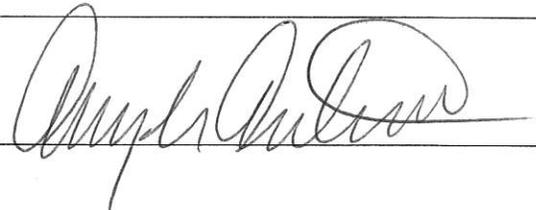SIGNATURE 

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location:
Serial Number: 80-002345
Core Version: 8105.48
08/21/2015 13:35:57

2018
JL
TS4

Standard Lot#: 25414100A4
Last Changed By: A. AUTRET #53546

QAS: A. AUTRET #53546
TUCSON PD

Operator: A. AUTRET #53546
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	13:36:58
Diagnostic Test	Pass	13:37:33
Air Blank	0.000	13:38:00
0.100 Cal Check	0.100	13:38:21
Air Blank	0.000	13:38:50
Subject Test	0.000	13:39:12
Air Blank	0.000	13:39:40
Wait	ABT*	
Air Blank	0.000	13:40:14

*Sequence Aborted

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location:
Serial Number: 80-002345
Core Version: 8105.48
08/21/2015 13:40:46

3 of 8
JL
TSG

Standard Lot#: 25414100A4
Last Changed By: A. AUTRET #53546

QAS: A. AUTRET #53546
TUCSON PD

Operator: A. AUTRET #53546
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	13:41:48
Diagnostic Test	Pass	13:42:23
Air Blank	0.000	13:42:50
0.100 Cal Check	RFI*	13:43:06
Air Blank	IPS**	13:43:15

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location:
Serial Number: 80-002345
Core Version: 8105.48
08/21/2015

13:43:38

4 of 8
Add
739

Standard Lot#: 25414100A4
Last Changed By: A. AUTRET #53546

QAS: A. AUTRET #53546
TUCSON PD

Operator: A. AUTRET #53546
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	13:44:31
Diagnostic Test	Pass	13:45:06
Air Blank	0.000	13:45:34
0.100 Cal Check	0.100	13:45:54
Air Blank	0.000	13:46:24
Subject Test	0.195*	13:49:31
Air Blank	0.000	13:50:04
Air Blank	0.000	13:50:31
Subject Test	XXX**	13:50:56
Air Blank	0.000	13:51:29

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location:
Serial Number: 80-002345
Core Version: 8105.48
08/21/2015 13:51:45

5 of 8
JL
TBH

Standard Lot#: 25414100A4
Last Changed By: A. AUTRET #53546

QAS: A. AUTRET #53546
TUCSON PD

Operator: A. AUTRET #53546
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	13:52:15
0.100 Cal Check	0.100	13:52:35
Air Blank	0.000	13:53:05

Preventative Maintenance
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

6 of 8
JLB

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

TSH
8.25.2015

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Amy Autret 53546 AGENCY Tucson Police Department

DATE 8/21/2015 TIME 1333

Intoxilyzer Serial # 80-002345 Location 000

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 101 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE 

DPS Form Exh G-2 (Rev. 05-01)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location:
Serial Number: 80-002345
Core Version: 8105.48
08/21/2015 13:32:17

70F8
AD

Standard Lot#: 25414100A4
Last Changed By: A. AUTRET #53546

TS4

QAS: A. AUTRET #53546
TUCSON PD

Operator: A. AUTRET #53546
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	13:32:46
0.100 Cal Check	0.100	13:33:07
Air Blank	RFI*	13:33:15
Air Blank	0.000	13:33:44

*RFI Detect

Repeat cal check -
RFI Detect -
officer walked past
instrument while
transmitting on
radio. AD

INTOXILYZER 8000
Location:
Serial Number: 80-002345
Core Version: 8105.48
08/21/2015 13:33:47

8018
dld
B4

Standard Lot#: 25414100A4
Last Changed By: A. AUTRET #53546

QAS: A. AUTRET #53546
TUCSON PD

Operator: A. AUTRET #53546
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	13:34:17
0.100 Cal Check	0.101	13:34:38
Air Blank	0.000	13:35:07

Preventative Maintenance
Performed -- Timer Reset