



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
1306 W. Miracle Mile, Tucson, AZ 85705

*Handwritten initials and date: 10/5/15*

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002345 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 9-17-2015 Cal 1 0.099 AC Cal 2 0.099 AC

Between the dates of 8-21-2015 and 9-17-2015 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on                     .  
The instrument was returned to service on                     .

*Handwritten signature of Terry S. Gallegos*

Terry S. Gallegos


9-29-2015

ADPS QAS Permit # 33169

SUBSCRIBED AND SWORN to before me this 29<sup>th</sup>, day of September, 20 15.

*Handwritten signature of Elma Escarcega*

Notary Public

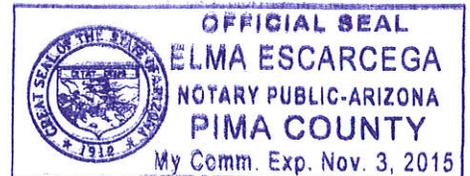

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

1 of 7  
TS4  
add

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Terry Gallegos 28847 AGENCY Tucson Police Department  
DATE 9-17-2015 TIME 1246  
INTOXILYZER SERIAL # 80-002345 LOCATION ODD

(  ) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

(  ) 1. Clock time check.  
(  ) 2. Date check.

OPERATIONAL TESTS

(  ) 1. Alcohol-free subject test result 0. 000 AC.  
(  ) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed  
(  ) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.  
(  ) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 099 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE Terry Gallegos

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000  
Location:  
Serial Number: 80-002345  
Core Version: 8105.48  
09/17/2015 12:48:16

2087  
TS4  
ole

---

Standard Lot#: 25414100A4  
Last Changed By: A. AUTRET #53546

---

QAS: T. GALLEGOS #28847  
TUCSON PD

---

Operator: T. GALLEGOS #28847  
TUCSON PD

---

Subject: TEST, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	12:49:18
Diagnostic Test	Pass	12:49:53
Air Blank	0.000	12:50:21
0.100 Cal Check	0.100	12:50:42
Air Blank	0.000	12:51:11
Subject Test	0.000	12:51:34
Air Blank	0.000	12:52:01
Wait	ABT*	
Air Blank	0.000	12:52:41

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location:  
Serial Number: 80-002345  
Core Version: 8105.48  
09/17/2015

12:52:55

30f7  
By  
[Signature]

---

Standard Lot#: 25414100A4  
Last Changed By: A. AUTRET #53546

---

QAS: T. GALLEGOS #28847  
TUCSON PD

---

Operator: T. GALLEGOS #28847  
TUCSON PD

---

Subject: TEST, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	12:53:47
Diagnostic Test	Pass	12:54:22
Air Blank	0.000	12:54:49
0.100 Cal Check	RFI*	12:55:03
Air Blank	IPS**	12:55:10

\*RFI Detect  
\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location:  
Serial Number: 80-002345  
Core Version: 8105.48  
09/17/2015 12:55:24

4 of 7  
TS4  
[Signature]

---

Standard Lot#: 25414100A4  
Last Changed By: A. AUTRET #53546

---

QAS: T. GALLEGOS #28847  
TUCSON PD

---

Operator: T. GALLEGOS #28847  
TUCSON PD

---

Subject: TEST, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	12:56:11
Diagnostic Test	Pass	12:56:46
Air Blank	0.000	12:57:13
0.100 Cal Check	0.099	12:57:34
Air Blank	0.000	12:58:03
Subject Test	0.156*	13:01:11
Air Blank	0.000	13:01:43
Air Blank	0.000	13:02:11
Subject Test	XXX**	13:02:35
Air Blank	0.000	13:03:08

\*Deficient Sample  
\*\*Invalid Sample  
Begin new deprivation period

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000

Location:

Serial Number: 80-002345

Core Version: 8105.48

09/17/2015

13:03:30

50f 7  
TJ  
D

-----  
Standard Lot#: 25414100A4

Last Changed By: A. AUTRET #53546  
-----

QAS: T. GALLEGOS #28847

TUCSON PD  
-----

Operator: T. GALLEGOS #28847

TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	13:03:59
0.100 Cal Check	0.099	13:04:20
Air Blank	0.000	13:04:50

Preventative Maintenance  
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

6077  
Dy  
Mo  
10/5/15

ARIZONA DEPARTMENT OF PUBLIC SAFETY  
STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000  
STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Terry Gallegos 28847 AGENCY Tucson Police Department  
DATE 9-17-2015 TIME 1246  
Intoxilyzer Serial # 80-002345 Location ODD

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 099 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE 

DPS Form Exh G-2 (Rev. 05-01)

