



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

TS4
1.13.2016

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002345 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 1/8/16 Cal 1 0.101 AC Cal 2 0.099 AC

Between the dates of 12/11/15 and 1/8/16 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on —.
The instrument was returned to service on —.

Amy L. Autret


1/11/16

ADPS QAS Permit # 33718

SUBSCRIBED AND SWORN to before me this 11, day of

January, 2016.

Notary Public

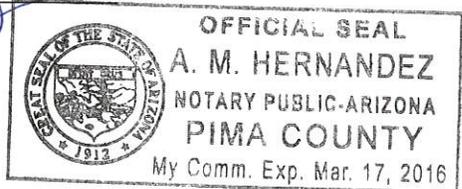

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1047
T84
1.13.2016

QA SPECIALIST Amy Autret 53546 AGENCY Tucson Police Department

DATE 1/8/16 TIME 1035

INTOXILYZER SERIAL # 80-002345 LOCATION 000

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 099 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE



Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location:
Serial Number: 80-002345
Core Version: 8105.48
01/08/2016

10:37:38

2 of 7
JLA
TS4

Standard Lot#: 25414100A4
Last Changed By: A. AUTRET #53546

QAS: A. AUTRET #53546
TUCSON PD

Operator: A. AUTRET #53546
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:38:35
Diagnostic Test	Pass	10:39:10
Air Blank	0.000	10:39:37
0.100 Cal Check	0.101	10:39:58
Air Blank	0.000	10:40:28
Subject Test	0.000	10:40:49
Air Blank	0.000	10:41:16
Wait	ABT*	
Air Blank	0.000	10:41:52

*Sequence Aborted

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location:
Serial Number: 80-002345
Core Version: 8105.48
01/08/2016 10:42:02

3 of 7
JLD
154

Standard Lot#: 25414100A4
Last Changed By: A. AUTRET #53546

QAS: A. AUTRET #53546
TUCSON PD

Operator: A. AUTRET #53546
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:42:59
Diagnostic Test	Pass	10:43:34
Air Blank	0.000	10:44:01
0.100 Cal Check	RFI*	10:44:16
Air Blank	IPS**	10:44:25

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location:
Serial Number: 80-002345
Core Version: 8105.48
01/08/2016 10:44:29

4 of 7
JLS
154

Standard Lot#: 25414100A4
Last Changed By: A. AUTRET #53546

QAS: A. AUTRET #53546
TUCSON PD

Operator: A. AUTRET #53546
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:45:26
Diagnostic Test	Pass	10:46:01
Air Blank	0.000	10:46:28
0.100 Cal Check	0.101	10:46:49
Air Blank	0.000	10:47:18
Subject Test	0.148*	10:50:26
Air Blank	0.000	10:50:58
Air Blank	0.000	10:51:25
Subject Test	XXX**	10:51:55
Air Blank	0.000	10:52:28

*Deficient Sample
**Invalid Sample
Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location:
Serial Number: 80-002345
Core Version: 8105.48
01/08/2016 10:52:41

5 of 7
JLB
T84

Standard Lot#: 25414100A4
Last Changed By: A. AUTRET #53546

QAS: A. AUTRET #53546
TUCSON PD

Operator: A. AUTRET #53546
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	10:53:09
0.100 Cal Check	0.099	10:53:30
Air Blank	0.000	10:54:00

Preventative Maintenance
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

6 of 7
T34
1.13.2016

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Amy Autret 53546 AGENCY Tucson Police Department

DATE 1/8/16 TIME 1035

Intoxilyzer Serial # 80-002345 Location 000

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 101 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE 

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location:
Serial Number: 80-002345
Core Version: 8105.48
01/08/2016 10:35:25

2 of 7
TS4

Standard Lot#: 25414100A4
Last Changed By: A. AUTRET #53546

QAS: A. AUTRET #53546
TUCSON PD

Operator: A. AUTRET #53546
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	10:35:54
0.100 Cal Check	0.101	10:36:15
Air Blank	0.000	10:36:44

Preventative Maintenance
Performed -- Timer Reset