



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
1306 W. Miracle Mile, Tucson, AZ 85705

TS4
7-1-2016


I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002345 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 6-24-16 Cal 1 0.100 AC Cal 2 0.099 AC

Between the dates of 5-27-16 and 6-24-16 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on _____.
The instrument was returned to service on _____.

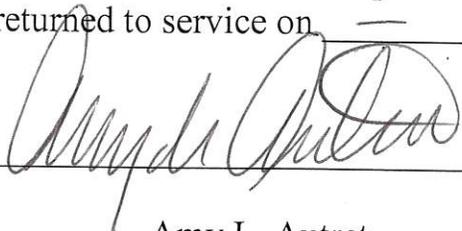


Amy L. Autret


6-29-16

ADPS QAS Permit # 33718

SUBSCRIBED AND SWORN to before me this 29th, day of June, 2016.

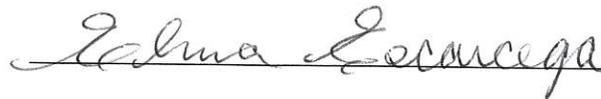

Notary Public

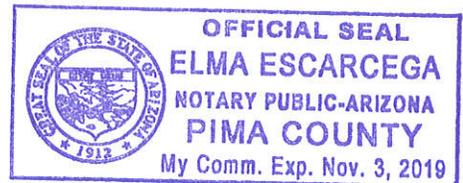

*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000
STANDARD QUALITY ASSURANCE PROCEDURE

1 of 7
MS
TSY
7-1-2016

QA SPECIALIST Amy Autret 53546 AGENCY Tucson Police Department
DATE 6-24-2016 TIME 1021
INTOXILYZER SERIAL # 80-002345 LOCATION ODD

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.


OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 099 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE

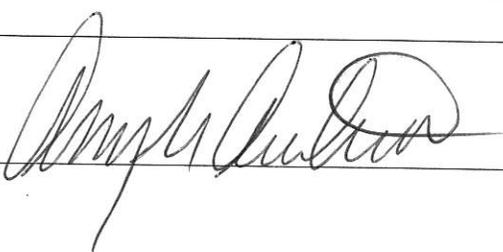

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location:
Serial Number: 80-002345
Core Version: 8105.48
06/24/2016 10:23:19

2 of 7
JJA
TSh

Standard Lot#: 25414100A4
Last Changed By: A. AUTRET #53546

QAS: A. AUTRET #53546
TUCSON PD

Operator: A. AUTRET #53546
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:24:14
Diagnostic Test	Pass	10:24:49
Air Blank	0.000	10:25:16
0.100 Cal Check	0.100	10:25:37
Air Blank	0.000	10:26:07
Subject Test	0.000	10:26:23
Air Blank	0.000	10:26:51
Wait	ABT*	
Air Blank	0.000	10:27:26

*Sequence Aborted

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location:
Serial Number: 80-002345
Core Version: 8105.48
06/24/2016

10:27:30

327
JLD
T34

Standard Lot#: 25414100A4
Last Changed By: A. AUTRET #53546

QAS: A. AUTRET #53546
TUCSON PD

Operator: A. AUTRET #53546
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:28:27
Diagnostic Test	Pass	10:29:02
Air Blank	0.000	10:29:29
0.100 Cal Check	RFI*	10:29:42
Air Blank	IPS**	10:29:51

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location:
Serial Number: 80-002345
Core Version: 8105.48
06/24/2016

10:29:55

4 of 7
JLB
TS4

Standard Lot#: 25414100A4
Last Changed By: A. AUTRET #53546

QAS: A. AUTRET #53546
TUCSON PD

Operator: A. AUTRET #53546
TUCSON PD

Subject: SQAP, SQAP
DOB: 11/11/1987
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:30:44
Diagnostic Test	Pass	10:31:19
Air Blank	0.000	10:31:47
0.100 Cal Check	0.100	10:32:08
Air Blank	0.000	10:32:37
Subject Test	0.333*	10:35:44
Air Blank	0.000	10:36:17
Air Blank	0.000	10:36:44
Subject Test	XXX**	10:37:11
Air Blank	0.000	10:37:44

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000

Location:

Serial Number: 80-002345

Core Version: 8105.48

06/24/2016

10:37:57

5 of 7
JW
TSH

Standard Lot#: 25414100A4

Last Changed By: A. AUTRET #53546

QAS: A. AUTRET #53546

TUCSON PD

Operator: A. AUTRET #53546

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	10:38:26
0.100 Cal Check	0.099	10:38:47
Air Blank	0.000	10:39:16

Preventative Maintenance
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6 of 7
JLB
TS4
7-1-2016

QA Specialist Amy Auter 53546 AGENCY Tucson Police Department
DATE 6-24-2016 TIME 10211
Intoxilyzer Serial # 80-002345 Location 000

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

Amy Auter

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000

Location:

Serial Number: 80-002345

Core Version: 8105.48

06/24/2016

10:21:50

7 of 7
JA
TS4

Standard Lot#: 25414100A4

Last Changed By: A. AUTRET #53546

QAS: A. AUTRET #53546

TUCSON PD

Operator: A. AUTRET #53546

TUCSON PD

Test	g/210L	Time
Air Blank	0.000	10:22:19
0.100 Cal Check	0.100	10:22:40
Air Blank	0.000	10:23:09

Preventative Maintenance
Performed -- Timer Reset