



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
1306 W. Miracle Mile, Tucson, AZ 85705

TS4  
7-27-2016

I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002345 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 7-22-16 Cal 1 0.101 AC Cal 2 0.099 AC

Between the dates of 6-24-16 and 7-22-16 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on —.  
The instrument was returned to service on —.


Amy L. Autret

7-26-16

ADPS QAS Permit # 33718

SUBSCRIBED AND SWORN to before me this 26<sup>th</sup>, day of July, 2016.

Notary Public

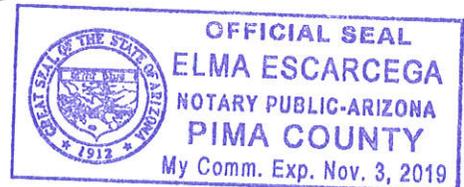

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

10F7  
Old  
By 7-27-2016

QA SPECIALIST Amy Autret 53546 AGENCY Tucson Police Department

DATE 7-22-16 TIME 1028

INTOXILYZER SERIAL # 80-002345 LOCATION 000

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 099 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS

---



---



---



---

SIGNATURE

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000

Location:

Serial Number: 80-002345

Core Version: 8105.48

07/22/2016

10:30:04

*2 of 7  
Add  
TS4*

---

Standard Lot#: 25414100A4

Last Changed By: A. AUTRET

#53546

---

QAS: A. AUTRET

#53546

TUCSON PD

---

Operator: A. AUTRET

#53546

TUCSON PD

---

Subject: SQAP, SQAP

DOB: 11/11/1987

Sex: F

Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	10:31:02
Diagnostic Test	Pass	10:31:37
Air Blank	0.000	10:32:04
0.100 Cal Check	0.101	10:32:25
Air Blank	0.000	10:32:54
Subject Test	0.000	10:33:11
Air Blank	0.000	10:33:38
Wait	ABT*	
Air Blank	0.000	10:34:15

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000

Location:

Serial Number: 80-002345

Core Version: 8105.48

07/22/2016

10:34:18

3.0 F7  
dlb  
784

---

Standard Lot#: 25414100A4

Last Changed By: A. AUTRET #53546

---

QAS: A. AUTRET #53546

TUCSON PD

---

Operator: A. AUTRET #53546

TUCSON PD

---

Subject: SQAP, SQAP

DOB: 11/11/1987

Sex: F

Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	10:35:13
Diagnostic Test	Pass	10:35:49
Air Blank	0.000	10:36:16
0.100 Cal Check	RFI*	10:36:29
Air Blank	IPS**	10:36:38

\*RFI Detect

\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location:  
Serial Number: 80-002345  
Core Version: 8105.48  
07/22/2016

10:36:45

4 of 7  
JJA  
T34

---

Standard Lot#: 25414100A4  
Last Changed By: A. AUTRET #53546

---

QAS: A. AUTRET #53546  
TUCSON PD

---

Operator: A. AUTRET #53546  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	10:37:36
Diagnostic Test	Pass	10:38:12
Air Blank	0.000	10:38:39
0.100 Cal Check	0.101	10:39:00
Air Blank	0.000	10:39:29
Subject Test	0.208*	10:42:37
Air Blank	0.000	10:43:10
Air Blank	0.000	10:43:37
Subject Test	XXX**	10:44:02
Air Blank	0.000	10:44:36

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000

Location:

Serial Number: 80-002345

Core Version: 8105.48

07/22/2016

10:45:00

-----  
Standard Lot#: 25414100A4

Last Changed By: A. AUTRET #53546  
-----

QAS: A. AUTRET #53546

TUCSON PD  
-----

Operator: A. AUTRET #53546

TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	10:45:29
0.100 Cal Check	0.099	10:45:50
Air Blank	0.000	10:46:20

Preventative Maintenance  
Performed -- Timer Reset

5 of 7  
TSG  
[Signature]

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6 of 7  
add  
TS4  
7-27-2016


QA Specialist Amy Autrict 53546 AGENCY Tucson Police Department

DATE 7-22-16 TIME 1028

Intoxilyzer Serial # 80-002345 Location 000

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 101 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE

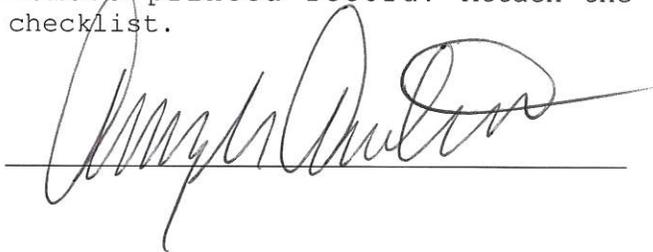

DPS Form Exh G-2 (Rev. 05-01)

