



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory  
1306 W. Miracle Mile, Tucson, AZ 85705

TS4  
8.22.2016


I hereby certify that these documents are a true and correct copy of the periodic maintenance\* for the Intoxilyzer 8000 Serial Number 80-002345 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 8-19-16 Cal 1 0.100 AC Cal 2 0.100 AC

Between the dates of 7-22-16 and 8-19-16 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on \_\_\_\_\_.  
The instrument was returned to service on \_\_\_\_\_.

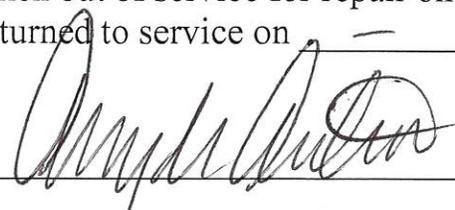  
\_\_\_\_\_

Amy L. Autret


8-22-16

ADPS QAS Permit # 33718

SUBSCRIBED AND SWORN to before me this 22nd, day of August, 2016.



Notary Public

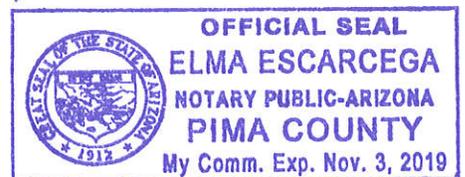

\*Includes documents specified in A.A.C. R13-10-104.A.5

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1 of 7  
dls  
B9  
8.22.2016

QA SPECIALIST Amy Autret 53546 AGENCY Tucson Police Department

DATE 8-19-16 TIME 1002

INTOXILYZER SERIAL # 80-002345 LOCATION OPD

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.  
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.  
Not a Successfully Completed Test Sequence printed  
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC  
Result: 0. 100 AC

Instrument is operating properly and accurately. YES  NO

COMMENTS  
\_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_

SIGNATURE 

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000

Location:

Serial Number: 80-002345

Core Version: 8105.48

08/19/2016

10:03:39

2 of 7  
JW  
BH

---

Standard Lot#: 25414100A4

Last Changed By: A. AUTRET #53546

---

QAS: A. AUTRET #53546

TUCSON PD

---

Operator: A. AUTRET #53546

TUCSON PD

---

Subject: SQAP, SQAP

DOB: 11/11/1987

Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	10:04:28
Diagnostic Test	Pass	10:05:03
Air Blank	0.000	10:05:30
0.100 Cal Check	0.101	10:05:51
Air Blank	0.000	10:06:21
Subject Test	0.000	10:06:36
Air Blank	0.000	10:07:03
Wait	ABT*	
Air Blank	0.000	10:07:38

\*Sequence Aborted

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location:  
Serial Number: 80-002345  
Core Version: 8105.48  
08/19/2016

10:07:42

3 of 7  
JLB  
B4

---

Standard Lot#: 25414100A4  
Last Changed By: A. AUTRET #53546

---

QAS: A. AUTRET #53546  
TUCSON PD

---

Operator: A. AUTRET #53546  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	10:08:33
Diagnostic Test	Pass	10:09:08
Air Blank	0.000	10:09:36
0.100 Cal Check	RFI*	10:09:49
Air Blank	IPS**	10:09:58

\*RFI Detect  
\*\*Improper Sample

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000  
Location:  
Serial Number: 80-002345  
Core Version: 8105.48  
08/19/2016

10:10:02

4 of 7  
dls  
B4

---

Standard Lot#: 25414100A4  
Last Changed By: A. AUTRET #53546

---

QAS: A. AUTRET #53546  
TUCSON PD

---

Operator: A. AUTRET #53546  
TUCSON PD

---

Subject: SQAP, SQAP  
DOB: 11/11/1987  
Sex: F Weight: 111

---

15 Minute Deprivation Period? Yes

---

Test	g/210L	Time
Air Blank	0.000	10:10:52
Diagnostic Test	Pass	10:11:27
Air Blank	0.000	10:11:54
0.100 Cal Check	0.101	10:12:15
Air Blank	0.000	10:12:44
Subject Test	0.179*	10:15:52
Air Blank	0.000	10:16:24
Air Blank	0.000	10:16:52
Subject Test	XXX**	10:17:19
Air Blank	0.000	10:17:53

\*Deficient Sample

\*\*Invalid Sample

Begin new deprivation period

**Not a Successfully Completed Test Sequence**

INTOXILYZER 8000

Location:

Serial Number: 80-002345

Core Version: 8105.48

08/19/2016

10:18:11

5 of 7  
JLH  
BH

-----  
Standard Lot#: 25414100A4

Last Changed By: A. AUTRET #53546

-----  
QAS: A. AUTRET #53546

TUCSON PD

-----  
Operator: A. AUTRET #53546

TUCSON PD

-----  
Test                      g/210L                      Time  
-----  
Air Blank                      0.000                      10:18:40  
0.100 Cal Check              0.100                      10:19:01  
Air Blank                      0.000                      10:19:31

Preventative Maintenance  
Performed -- Timer Reset

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY  
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES  
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

6 of 7  
dls  
Bq  
8.22.2016

QA Specialist Amy Autret 53546 AGENCY Tucson Police Department  
DATE 8-19-16 TIME 1002  
Intoxilyzer Serial # 80-002345 Location 000

- ( ✓ ) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.  
OR  
Pour a standard alcohol concentration solution \_\_\_\_\_ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach  $34^{\circ}\text{C} \pm 0.2^{\circ}\text{C}$
- ( ✓ ) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- ( ✓ ) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- ( ✓ ) 4. Air blank completed.
- ( ✓ ) 5. Calibration check completed. Test results 0. 100 AC.
- ( ✓ ) 6. Air blank completed.
- ( ✓ ) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE



DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000  
Location:  
Serial Number: 80-002345  
Core Version: 8105.48  
08/19/2016 10:02:11

7 of 7  
als  
BSG

-----  
Standard Lot#: 25414100A4  
Last Changed By: A. AUTRET #53546  
-----

QAS: A. AUTRET #53546  
TUCSON PD  
-----

Operator: A. AUTRET #53546  
TUCSON PD  
-----

Test	g/210L	Time
Air Blank	0.000	10:02:39
0.100 Cal Check	0.100	10:03:00
Air Blank	0.000	10:03:30

Preventative Maintenance  
Performed -- Timer Reset