



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

TJY
12-21-07

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002759 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 12/13/07 Cal 1 0.102 AC Cal 2 0.100 AC

Between the dates of 11/21/07 and 12/13/07 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -.
The instrument was returned to service on -.


Sarah L. Cione

Sarah L. Cione

12/19/07

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 19, day of December, 2007.

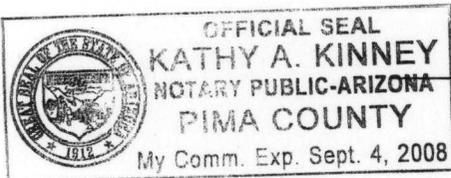

Kathy A Kinney

Notary Public

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

1 of 7 sc

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

134

STANDARD QUALITY ASSURANCE PROCEDURE

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department

DATE 12/13/07 TIME 14:52

INTOXILYZER SERIAL # 80-002759 LOCATION Crime Lab

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 100 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Sarah Cione

2 of 7 sc

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002759
12/13/2007

14:54:13

BY

Standard Lot#: 725401I
Last Changed By: C. CONLEY #34517

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, TEST
DOB: 11/11/1111
Sex: F Weight: 1.11

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|------------------|--------|----------|
| Air Blank | 0.000 | 14:55:01 |
| Diagnostic Test | Pass | 14:55:29 |
| Air Blank | 0.000 | 14:55:56 |
| 0.100 Cal Check | 0.101 | 14:56:15 |
| Air Blank | 0.000 | 14:56:44 |
| Subject Test | 0.000 | 14:57:01 |
| Air Blank | 0.000 | 14:57:28 |
| Five Minute Wait | ABT* | |
| Air Blank | 0.000 | 14:58:07 |

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002759
12/13/2007

14:58:29

3 of 7 sc

TSY

Standard Lot#: 725401I
Last Changed By: C. CONLEY #34517

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------|--------|----------|
| Air Blank | RFI* | 14:59:09 |
| Air Blank | IPS** | 14:59:32 |

*RFI Detect
**Improper Sample

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002759
12/13/2007

14:59:56

4 of 7 sc

TSY

Standard Lot#: 725401I
Last Changed By: C. CONLEY #34517

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 15:00:56 |
| Diagnostic Test | Pass | 15:01:23 |
| Air Blank | 0.000 | 15:01:50 |
| 0.100 Cal Check | 0.101 | 15:02:09 |
| Air Blank | 0.000 | 15:02:39 |
| Subject Test | 0.072* | 15:05:47 |
| Air Blank | 0.000 | 15:06:21 |
| Air Blank | 0.000 | 15:06:48 |
| Subject Test | XXX** | 15:07:11 |
| Air Blank | 0.000 | 15:07:46 |

*Deficient Sample

**Invalid Sample

Begin new deprivation period

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002759
12/13/2007 15:08:15

5 of 7 sc

BH

Standard Lot#: 725401I
Last Changed By: C. CONLEY #34517

QAS: S. CIONE #50616
TUCSON PD

| Test | g/210L | Time |
|-----------|--------|----------|
| Air Blank | 0.000 | 15:08:42 |
| Cal Check | 0.100 | 15:09:01 |
| Air Blank | 0.000 | 15:09:31 |

Sarah Cione

Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

6 of 7sc

ARIZONA DEPARTMENT OF PUBLIC SAFETY

174

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD CALIBRATION CHECK PROCEDURE

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 12/13/07 TIME 14:52
Intoxilyzer Serial # 80-002759 Location Crime Lab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 102 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002759
12/13/2007. 14:52:42

7 of 7

154

Standard Lot#: 725401I
Last Changed By: C. CONLEY #34517

QAS: S. CIONE #50616
TUCSON PD

| Test | g/210L | Time |
|-----------|--------|----------|
| Air Blank | 0.000 | 14:53:09 |
| Cal Check | 0.102 | 14:53:28 |
| Air Blank | 0.000 | 14:53:57 |

Sarah Cione
Operator's Signature