



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

754
5-2-08

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-202759 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 5-2-08 Cal 1 0.102 AC Cal 2 0.101 AC

Between the dates of 4-4-08 and 5-2-08 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.
The instrument was returned to service on N/A.

Clare Conley


Clare L. Conley

5-2-08

ADPS QAS Permit # 33314

SUBSCRIBED AND SWORN to before me this 2, day of May, 2008.

Kathy A Kinney
Notary Public

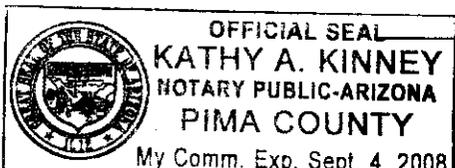

My Comm. Exp. Sept. 4, 2008

*Includes documents specified in A.A.C. R13-10-104.A.5

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

1 of 7 ar
TS4
5.2.08

QA SPECIALIST Clare Conley AGENCY Tucson Police Department

DATE 5-2-08 TIME 0819

INTOXILYZER SERIAL # 80-002759 LOCATION ODM

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 101 AC

Instrument is operating properly and accurately. YES NO

COMMENTS None.

SIGNATURE Clare Conley

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002759
05/02/2008

08:20:23

2 of 7 cu

TS4

Standard Lot#: 725401I
Last Changed By: C. CONLEY #34517

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP

DOB: 11/11/1111

Sex: M Weight: 111

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 08:21:10 |
| Diagnostic Test | Pass | 08:21:38 |
| Air Blank | 0.000 | 08:22:05 |
| 0.100 Cal Check | 0.102 | 08:22:24 |
| Air Blank | 0.000 | 08:22:53 |
| Subject Test | 0.000 | 08:23:08 |
| Air Blank | ABT* | 08:23:17 |
| Air Blank | 0.000 | 08:23:46 |

*Sequence Aborted

Not a Successfully Completed Test Sequence

Pursuant to section 28-1527, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002759
05/02/2008

08:23:53

387 cu
134

Standard Lot#: 725401I
Last Changed By: C. CONLEY #34517

OAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP

DOB: 11/11/1111

Sex: M Weight: 111

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------|--------|----------|
| Air Blank | RFI* | 08:24:23 |
| Air Blank | IPS** | 08:24:31 |

*RFI Detect

**Improper Sample

Not a Successfully Completed Test Sequence

Pursuant to section 28-1527, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002759
05/02/2008

08:24:48

4 of 7 cu
TS4

Standard Lot#: 725401I
Last Changed By: C. CONLEY #34517

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP

DOB: 11/11/1111

Sex: M Weight: 111

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 08:25:39 |
| Diagnostic Test | Pass | 08:26:06 |
| Air Blank | 0.000 | 08:26:34 |
| 0.100 Cal Check | 0.102 | 08:26:52 |
| Air Blank | 0.000 | 08:27:22 |
| Subject Test | 0.233* | 08:30:30 |
| Air Blank | 0.000 | 08:31:03 |
| Air Blank | 0.000 | 08:31:30 |
| Subject Test | XXX** | 08:32:07 |
| Air Blank | 0.000 | 08:32:41 |

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

Pursuant to section 28-1527, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002759
05/02/2008

08:32:59

507 cu
TS4

Standard Lot#: 725401I
Last Changed By: C. CONLEY #34517

QAS: C. CONLEY #34517
TUCSON PD

| Test | g/210L | Time |
|-----------|--------|----------|
| Air Blank | 0.000 | 08:33:27 |
| Cal Check | 0.101 | 08:33:46 |
| Air Blank | 0.000 | 08:34:15 |

Clare Conley

Operator's Signature

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000
STANDARD CALIBRATION CHECK PROCEDURE

6867 cu
TS4
5-2-08

QA Specialist Clare Conley AGENCY Tucson Police Department
DATE 5-2-08 TIME 0819
Intoxilyzer Serial # 80-002759 Location ODM

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 102 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Clare Conley

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002759
05/02/2008

08:19:00

7 of 7 cu

154

Standard Lot#: 725401I
Last Changed By: C. CONLEY #34517

QAS: C. CONLEY #34517
TUCSON PD

| Test | g/210L | Time |
|-----------|--------|----------|
| Air Blank | 0.000 | 08:19:28 |
| Cal Check | 0.102 | 08:19:46 |
| Air Blank | 0.000 | 08:20:16 |

Clare Conley
Operator's Signature