



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

OK
7-10-08

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002759 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 6/27/08 Cal 1 0.098 AC Cal 2 0.095 AC

Between the dates of 5/30/08 and 6/27/08 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on -
The instrument was returned to service on -

Sarah Cione


Sarah L. Cione

7/3/08

ADPS QAS Permit # 33355

SUBSCRIBED AND SWORN to before me this 3 day of July, 2008.

Kathy A Kinney
Notary Public

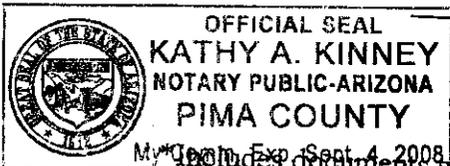

Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

1 of 7 sc

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000
STANDARD QUALITY ASSURANCE PROCEDURE

Ch
7-10-08

QA SPECIALIST Sarah Cione AGENCY Tucson Police Department
DATE 6/27/08 TIME 10:54
INTOXILYZER SERIAL # 80-002759 LOCATION Crime Lab

- (✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 098 AC

Instrument is operating properly and accurately. YES NO

COMMENTS

SIGNATURE Sarah Cione

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

2 of 75C

OK

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002759
06/27/2008

10:56:01

Standard Lot#: 725401I
Last Changed By: C. CONLEY #34517

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	10:57:01
Diagnostic Test	Pass	10:57:28
Air Blank	0.000	10:57:55
0.100 Cal Check	0.099	10:58:14
Air Blank	0.000	10:58:44
Subject Test	0.000	10:59:24
Air Blank	0.000	10:59:51
Five Minute Wait	ABT*	
Air Blank	0.000	11:00:25

*Sequence Aborted

Not a Successfully Completed Test Sequence

3 of 7sc

CM

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002759
06/27/2008 11:00:44

Standard Lot#: 725401I
Last Changed By: C. CONLEY #34517

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	11:01:08
Air Blank	IPS**	11:01:16

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002759
06/27/2008

11:01:40

4 of 7sc

cm

Standard Lot#: 725401I
Last Changed By: C. CONLEY #34517

QAS: S. CIONE #50616
TUCSON PD

Operator: S. CIONE #50616
TUCSON PD

Subject: SQAP, TEST
DOB: 11/11/1111
Sex: F Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	11:02:24
Diagnostic Test	Pass	11:02:51
Air Blank	0.000	11:03:18
0.100 Cal Check	0.099	11:03:37
Air Blank	0.000	11:04:07
Subject Test	0.066*	11:07:15
Air Blank	0.000	11:07:47
Air Blank	0.000	11:08:14
Subject Test	XXX**	11:08:47
Air Blank	0.000	11:09:20

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002759
06/27/2008

11:09:55

Standard Lot#: 725401I
Last Changed By: C. CONLEY #34517

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	11:10:22
Cal Check	0.098	11:10:41
Air Blank	0.000	11:11:11

Sarah Cione

Operator's Signature

5 of 7 sc

Ch

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

6 of 7 SC

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000
STANDARD CALIBRATION CHECK PROCEDURE

OK
7-10-08

QA Specialist Sarah Cione AGENCY Tucson Police Department
DATE 6/27/08 TIME 10:54
Intoxilyzer Serial # 80-002759 Location Crime Lab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 098 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Sarah Cione

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002759
06/27/2008

10:54:32

Standard Lot#: 725401I
Last Changed By: C. CONLEY #34517

QAS: S. CIONE #50616
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	10:54:59
Cal Check	0.098	10:55:18
Air Blank	0.000	10:55:47

Sarah Cione

Operator's Signature

7 of 75c

OK