



CERTIFICATION by CUSTODIAN of RECORDS
Tucson Police Department Crime Laboratory

TSY
9.25.08

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002759 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 9-19-08 Cal 1 0.097 AC Cal 2 0.097 AC

Between the dates of 8-22-08 and 9-19-08 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.
The instrument was returned to service on N/A.

Clare Conley

Clare L. Conley


9-25-08

ADPS QAS Permit # 33314

SUBSCRIBED AND SWORN to before me this 25th, day of September, 2008.

A.M. Hernandez

Notary Public

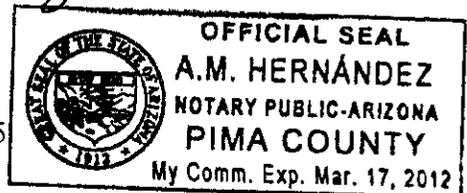

*Includes documents specified in A.A.C. R13-10-104.A.5

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104 (A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000
STANDARD QUALITY ASSURANCE PROCEDURE

1587
Qu
756
9.25.08

QA SPECIALIST Clare Conley AGENCY Tucson Police Department

DATE 9-19-08 TIME 0756

INTOXILYZER SERIAL # 80-002759 LOCATION Crime Lab

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 097 AC

Instrument is operating properly and accurately. YES NO

COMMENTS None.

SIGNATURE Clare Conley

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002759
09/19/2008

07:58:34

267
em
TSy

Standard Lot#: 725401I
Last Changed By: C. CONLEY #34517

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	07:59:26
Diagnostic Test	Pass	07:59:53
Air Blank	0.000	08:00:21
0.100 Cal Check	0.097	08:00:39
Air Blank	0.000	08:01:09
Subject Test	0.000	08:01:28
Air Blank	ABT*	08:01:35
Air Blank	0.000	08:02:03

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002759
09/19/2008

08:02:12

317
cu
dy

Standard Lot#: 725401I
Last Changed By: C. CONLEY #34517

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	RFI*	08:02:43
Air Blank	IPS**	08:02:51

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number : 80-002759
09/19/2008

08:03:21

467
On
104

Standard Lot#: 725401I
Last Changed By: C. CONLEY #34517

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

Test	g/210L	Time
Air Blank	0.000	08:04:07
Diagnostic Test	Pass	08:04:35
Air Blank	0.000	08:05:02
0.100 Cal Check	0.098	08:05:21
Air Blank	0.000	08:05:50
Subject Test	0.167*	08:08:59
Air Blank	0.000	08:09:31
Air Blank	0.000	08:09:58
Subject Test	XXX**	08:10:33
Air Blank	0.000	08:11:06

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002759
09/19/2008 08:11:28

Standard Lot#: 725401I
Last Changed By: C. CONLEY #34517

QAS: C. CONLEY #34517
TUCSON PD

5567
Cu
154

Test	g/210L	Time
Air Blank	0.000	08:11:56
Cal Check	0.097	08:12:14
Air Blank	0.000	08:12:44

Clare Conley

Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000
STANDARD CALIBRATION CHECK PROCEDURE

b567 Cu
TSH 9.25.08

QA Specialist Clare Conley AGENCY Tucson Police Department
DATE 9-19-08 TIME 0756
Intoxilyzer Serial # 80-002759 Location Crime Lab

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 097 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Clare Conley

DPS Form Exh G-2 (Rev. 05-01)

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: ARIZONA
Serial Number: 80-002759
09/19/2008 07:56:09

Standard Lot#: 725401I
Last Changed By: C. CONLEY #34517

QAS: C. CONLEY #34517
TUCSON PD

Test	g/210L	Time
Air Blank	0.000	07:56:36
Cal Check	0.097	07:56:55
Air Blank	0.000	07:57:24

7 of 7 Ch

TS4
9.25.08

Clare Conley

Operator's Signature