



CERTIFICATION by CUSTODIAN of RECORDS

Tucson Police Department Crime Laboratory
270 S. Stone Avenue Tucson, AZ 85701-1917

TS4
4-10-09

I hereby certify that these documents are a true and correct copy of the periodic maintenance* for the Intoxilyzer 8000 Serial Number 80-002759 maintained by the Tucson Police Department Crime Laboratory pursuant to the requirements of the Arizona Department of Public Safety for a quality assurance program.

A quality assurance specialist (QAS) tested the instrument on the following dates. The acceptable accuracy limits, as defined in A.A.C. R13-10-104.A.2, for a 0.100AC solution are 0.090AC to 0.110AC.

Date: 4-2-09 Cal 1 0.100 AC Cal 2 0.100 AC

Between the dates of 3-6-09 and 4-2-09 the instrument was operating properly and accurately.

The instrument was taken out of service for repair on N/A.
The instrument was returned to service on N/A.

Clare Conley

Clare L. Conley


4-9-09

ADPS QAS Permit # 33314

SUBSCRIBED AND SWORN to before me this 9th, day of April, 2009.

A.M. Hernandez

Notary Public

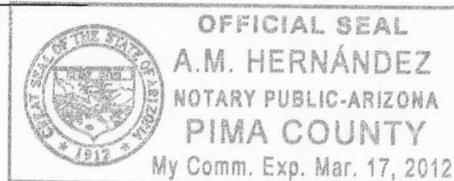

*Includes documents specified in A.A.C. R13-10-104.A.5

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

EXHIBIT G-4

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

ARIZONA DEPARTMENT OF PUBLIC SAFETY

STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000

STANDARD QUALITY ASSURANCE PROCEDURE

187
cu
TS4
4.10.09

QA SPECIALIST Clare Conley AGENCY Tucson Police Department

DATE 4-2-09 TIME 0630

INTOXILYZER SERIAL # 80-00 2759 LOCATION ODG

(✓) 1. Display reads "PUSH BUTTON TO START"

DIAGNOSTIC TESTS

- (✓) 1. Clock time check.
- (✓) 2. Date check.

OPERATIONAL TESTS

- (✓) 1. Alcohol-free subject test result 0. 000 AC.
- (✓) 2. Error recognition logic system functioning.
Not a Successfully Completed Test Sequence printed
- (✓) 3. Proper sample recognition system.
Not a Successfully Completed Test Sequence printed
Deficient sample printed.
- (✓) 4. Standard Calibration Check standard 0. 100 AC
Result: 0. 100 AC

Instrument is operating properly and accurately. YES NO

COMMENTS None.

SIGNATURE Clare Conley

Pursuant to section 28-1327, Arizona Revised Statutes, this document is a certified duplicate of the information contained in computer storage devices of the Tucson Police Department.

INTOXILYZER 8000
Location: TUCSON PD
Serial Number : 80-002759
04/02/2009

06:32:13

287
cu
104

Standard Lot#: 808101G
Last Changed By: S. CIONE #50616

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 06:33:01 |
| Diagnostic Test | Pass | 06:33:28 |
| Air Blank | 0.000 | 06:33:55 |
| 0.100 Cal Check | 0.101 | 06:34:14 |
| Air Blank | 0.000 | 06:34:44 |
| Subject Test | 0.000 | 06:35:00 |
| Air Blank | ABT* | 06:35:07 |
| Air Blank | 0.000 | 06:35:35 |

*Sequence Aborted

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number : 80-002759
04/02/2009

06:35:39

387
Cu
154

Standard Lot#: 808101G
Last Changed By: S. CIONE #50616

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------|--------|----------|
| Air Blank | RFI* | 06:36:09 |
| Air Blank | IPS** | 06:36:18 |

*RFI Detect
**Improper Sample

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number : 80-002759
04/02/2009

06:36:32

4.67
Cu
134

Standard Lot#: 808101G
Last Changed By: S. CIONE #50616

QAS: C. CONLEY #34517
TUCSON PD

Operator: C. CONLEY #34517
TUCSON PD

Subject: TEST, SQAP
DOB: 11/11/1111
Sex: M Weight: 111

15 Minute Deprivation Period? Yes

| Test | g/210L | Time |
|-----------------|--------|----------|
| Air Blank | 0.000 | 06:37:20 |
| Diagnostic Test | Pass | 06:37:47 |
| Air Blank | 0.000 | 06:38:15 |
| 0.100 Cal Check | 0.101 | 06:38:33 |
| Air Blank | 0.000 | 06:39:03 |
| Subject Test | 0.274* | 06:42:11 |
| Air Blank | 0.000 | 06:42:43 |
| Air Blank | 0.000 | 06:43:11 |
| Subject Test | XXX** | 06:43:45 |
| Air Blank | 0.000 | 06:44:19 |

*Deficient Sample

**Invalid Sample

Begin new deprivation period

Not a Successfully Completed Test Sequence

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002759
04/02/2009 06:44:35

5067
cu
124

Standard Lot#: 808101G
Last Changed By: S. CIONE #50616

QAS: C. CONLEY #34517
TUCSON PD

| Test | g/210L | Time |
|-----------|--------|----------|
| Air Blank | 0.000 | 06:45:03 |
| Cal Check | 0.100 | 06:45:21 |
| Air Blank | 0.000 | 06:45:51 |

Clare Conley

Operator's Signature

EXHIBIT G-2

THIS REPORT PREPARED PURSUANT TO DUTY IMPOSED BY
A.A.C. R13-10-104(A)

6567
cu

ARIZONA DEPARTMENT OF PUBLIC SAFETY
STANDARD QUALITY ASSURANCE PROCEDURES
INTOXILYZER MODEL 8000
STANDARD CALIBRATION CHECK PROCEDURE

TS4
4.10.09

QA Specialist Clare Conley AGENCY Tucson Police Department
DATE 4-2-09 TIME 0630
Intoxilyzer Serial # 80-002759 Location ODE

- (✓) 1. Ensure that gas tank is attached to instrument and contains a standard alcohol concentration solution 0. 100 AC.
OR
Pour a standard alcohol concentration solution _____ AC, into a clean dry simulator and assemble the simulator. Ensure that a tight seal has been made. Turn on the simulator and allow temperature to reach 34°C ± 0.2°C
- (✓) 2. Intoxilyzer 8000 display reads "PUSH BUTTON TO START"
- (✓) 3. Go to the "Control Testing Menu". Select "D" for dry control test or "W" for wet control test. After selection is made press ENTER.
- (✓) 4. Air blank completed.
- (✓) 5. Calibration check completed. Test results 0. 100 AC.
- (✓) 6. Air blank completed.
- (✓) 7. Remove printed record. Attach the record to the completed checklist.

SIGNATURE Clare Conley

DPS Form Exh G-2 (Rev. 05-01)

INTOXILYZER 8000
Location: TUCSON PD
Serial Number: 80-002759
04/02/2009 06:30:52

7587
cu

Standard Lot#: 808101G
Last Changed By: S. CIONE #50616

BSH

QAS: C. CONLEY #34517
TUCSON PD

| Test | g/210L | Time |
|-----------|--------|----------|
| Air Blank | 0.000 | 06:31:20 |
| Cal Check | 0.100 | 06:31:38 |
| Air Blank | 0.000 | 06:32:08 |

Clare Conley

Operator's Signature